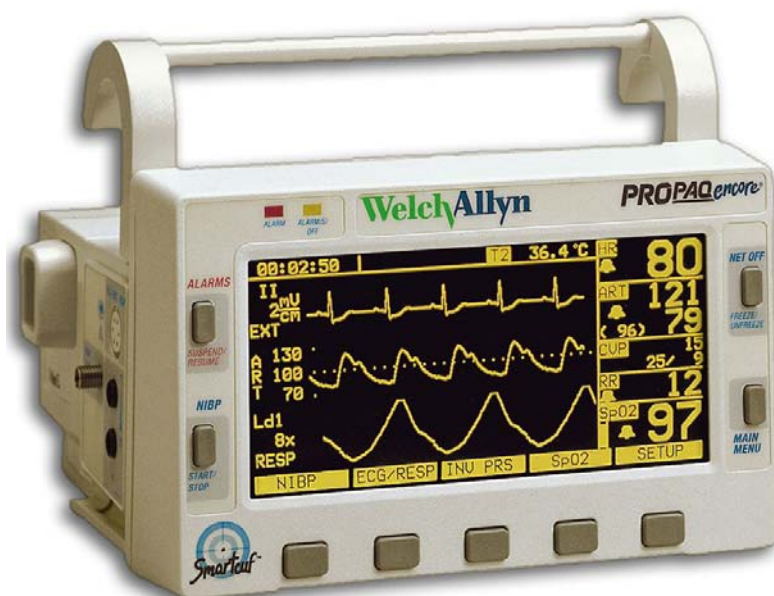


# Propaq® Encore

## Vital Signs Monitor

## Reference Guide



Model 202EL, 204EL, 206EL  
Software version 2.4X

**WelchAllyn®**

Advancing Frontline Care™

**Copyright © 2004 by Welch Allyn.** Welch Allyn® is a registered trademark of Welch Allyn. Propaq®, Propaq Encore®, and Acuity® are registered trademarks of Welch Allyn. Welch Allyn is protected under various patents and patents pending. Masimo® and SET® are registered trademarks and Signal Extraction Technology™ is a trademark of Masimo Corporation. NELLCOR®, DURASENSOR®, C-LOCK® and OXISENSOR® are registered trademarks of Nellcor Puritan Bennett, Inc. LIFEPAK 5® and LIFEPAK 6s® are registered trademarks of Medtronic Physio-Control Corporation.

**Copyright Notice:** Software in this Product is Copyright© 2004 by Welch Allyn or its vendors. All rights are reserved. The software is protected by United States of America copyright laws and international treaty provisions applicable all over the world. Under such laws, you are licensed to use the copy of the software incorporated with this instrument as intended in the operation of the product in which it is embedded, but the software may not be copied, decompiled, reverse-engineered, disassembled or otherwise reduced to human-perceivable form. This is not a sale of the software or any copy of the software; all right, title and ownership of the software remains with Welch Allyn or its vendors. Welch Allyn will make available specifications necessary for interoperability of this software on request; however, users should be aware that use of Welch Allyn hardware and software with devices or software not sold by Welch Allyn or its authorized dealers and affiliates may lead to erroneous results and consequent danger in patient care, and may also void Welch Allyn's warranty.

**Disclaimers:** Welch Allyn cautions the reader of this manual:

- This manual may be wholly or partially subject to change without notice.
- All rights are reserved. No one is permitted to reproduce or duplicate, in any form, the whole or part of this manual without permission from Welch Allyn.
- Welch Allyn will not be responsible for any injury to the user or other person(s) that may result from accidents during operation of the Propaq monitor.
- Welch Allyn assumes no responsibility for usage not in accordance with this manual that results in illegal or improper use of the Propaq monitor.
- No implied license: Possession or purchase of this device does not convey any express or implied license to use the device with unauthorized sensors or cables which would, alone, or in combination with this device, fall within the scope of one or more of the patents relating to this device.

**Welch Allyn Technical Support:**

<i>USA</i>	1-800-535-6663	<i>France</i>	(+33) 1-60-09-33-66	<i>Australia</i>	(+61) 2-9638-3000
<i>Latin America</i>	(+1) 305-669-9591	<i>Germany</i>	(+49) 7477-927-173	<i>Singapore</i>	(+65) 6291-0882
<i>European Call Center</i>	(+353) 469-067-790	<i>Canada</i>	1-800-561-8797	<i>Japan</i>	(+81) 3-5212-7391
<i>United Kingdom</i>	0-207-365-6780	<i>South Africa</i>	(+27) 11-777-7509	<i>China</i>	(+86) 21-6327-9631

**For information concerning this document or any Welch Allyn product, contact:**

**Welch Allyn**  
Customer Service  
8500 SW Creekside Place  
Beaverton, Oregon 97008-7107 USA

**Within USA, toll free:**

Phone: (800) 289-2500  
Phone Technical Services: (800) 535-6663

**WorldWide:**

Phone: (503) 530-7500  
Fax: (503) 526-4200  
Fax Technical Services: (503) 526-4970

**Internet:** <http://www.monitoring.welchallyn.com>  
E-mail Technical Services:  
[solutions@monitoring.welchallyn.com](mailto:solutions@monitoring.welchallyn.com)  
E-mail Marketing Dept.:  
[marketing@monitoring.welchallyn.com](mailto:marketing@monitoring.welchallyn.com)

**Welch Allyn European Customer Service**

IPA Business Park, Dublin Road  
Navan, County Meath, Ireland  
Phone: 353-46-67700  
Fax: 353-46-27128

**Welch Allyn U.K. Ltd.**

Cublington Road, Aston Abbotts  
Buckinghamshire HP22 4ND, England  
Phone: 44-1296-682140  
Fax: 44-1296-682104

**Welch Allyn Italia**

Via Napo Torriani, 29, 20124 Milan, Italy  
Phone: 39-02-6699-291  
Fax: 39-02-6671-3599

**Welch Allyn GmbH: Germany**

Postfach 31, Zollerstrasse 2-4  
72417 Jungingen, Germany  
Phone: 49-7477-92-710  
Fax: 49-7477-92-7190

**Welch Allyn: France**

814 Rue Charles de Gaulle  
77100 Mareuil les Meaux, France  
Phone: 01-6009-3366  
Fax: 01-6009-6797

**Welch Allyn: Pacific**

P.O. Box 39-293 Howick, Auckland, New Zealand  
Phone: 64-9-532-9524  
Fax: 64-9-532-9526

**Welch Allyn: Asia**

Room1002, 10/F Tung Sun Comm. Centre  
194-200 Lockhart Road, Wanchai, H.K.  
Phone: (852) 9016-7812  
Fax: (852) 2535-5650

**Welch Allyn, Inc.: Latin America**

MD International, 11300 NW 41st Street  
Miami, FL 33172 USA  
Phone: (305) 669-9003  
Fax: (305) 669-8951

Reorder Part No: 810-1719-XX

Manual Part No: 810-0640-03 Rev. A 8/04

Printed in USA

# Contents

---

<b>General Information . . . . .</b>	<b>9</b>
Safety Summary . . . . .	9
Symbols . . . . .	11
Propaq Encore Documentation . . . . .	13
The Documentation Set . . . . .	13
About This Reference Guide . . . . .	13
Statement of Expectations of the Reader . . . . .	13
 <b>Getting Started . . . . .</b>	 <b>15</b>
Introducing the Propaq Encore. . . . .	15
Intended Use . . . . .	15
Propaq Encore Models and Options . . . . .	15
Expansion Module . . . . .	16
Propaq Encore Pulse Oximetry Option (SpO <sub>2</sub> ) . . . . .	16
Capnography (CO <sub>2</sub> ) Options . . . . .	16
Impedance Pneumography (RESP) Option . . . . .	17
Printer Option . . . . .	17
Propaq-to-Acuity Option . . . . .	17
Modem-Propaq Option . . . . .	17
HP-compatible Side Panel Option . . . . .	17
Using the Propaq Encore . . . . .	18
System Controls (Right Side Panel) . . . . .	18
Alarm Lights . . . . .	19
Power-up Screen . . . . .	20
Patient Connections . . . . .	22
Propaq Encore Display . . . . .	24
Propaq Encore Buttons . . . . .	25
Propaq Encore Menus . . . . .	26
Monitor Setup . . . . .	29
Selecting Waveforms for Display . . . . .	31
Setting the Time and Date . . . . .	31
Changing the Date Format, Filter, and Units . . . . .	32
Setting the Current, Custom, and Powerup Modes . . . . .	33
Printer Functions . . . . .	35
Learning the Propaq Encore . . . . .	37
Using In-Service Mode . . . . .	37

What You Can Do With In-Service Mode .....	37
What You Cannot Do With In-Service Mode .....	38

## **Patient Monitoring ..... 39**

ECG/RESP .....	39
Intended Use—Impedance Pneumography (RESP) .....	39
Intended Use—ECG .....	40
ECG Connector and Applicable Accessories .....	41
Preparation .....	42
How ECG/RESP is Displayed .....	45
Using the Propaq Encore With Pacemaker Patients .....	47
Using the Filter to Better Display a Waveform .....	48
ECG Messages .....	48
RESP Messages .....	49
Invasive Pressure .....	50
Intended Use .....	50
Invasive Pressure Connectors and Transducers .....	50
Preparation .....	51
Rezeroing a Transducer .....	52
How Invasive Pressure is Displayed .....	53
Invasive Pressure Messages .....	54
NIBP .....	56
Intended Use .....	56
Improve NIBP Accuracy with Smartcuf™ .....	57
NIBP Connector and Cuffs .....	58
Preparation .....	59
Important Information About Automatic Measurements .....	62
NIBP Messages .....	62
Temperature .....	66
Intended Use .....	66
Temperature Connectors and Probes .....	66
Preparation .....	66
How Temperature is Displayed .....	66
Temperature Messages .....	67
Pulse Oximetry (SpO <sub>2</sub> ) .....	68
Perform SpO <sub>2</sub> Monitoring with Masimo Option .....	69
Perform SpO <sub>2</sub> Monitoring with Nellcor Option .....	70
Perform SpO <sub>2</sub> “Spot-Check” Monitoring .....	73
Capnography (CO <sub>2</sub> ) .....	75
Intended Use .....	75

CO2 Measurements and Display .....	76
CO2 Display Menus and Status Window .....	78
Mainstream CO2 Monitoring .....	80
How to Set Up the CO2 Channel and Alarm Limits .....	81
Sidestream CO2 Monitoring .....	83
CO2 Messages .....	86
 <b>Alarms and Limits .....</b>	 <b>91</b>
Alarms Window and Menu .....	91
Alarm Status Window .....	91
Alarms Menu .....	91
Adjusting the Alarm Tone Volume .....	91
Types of Propaq Alarms .....	93
Patient Alarms: Definitions and Indications .....	93
Life-Threatening Alarms .....	94
Temporary Alarms .....	94
Responding to Patient Alarms .....	95
Responding to Life-Threatening Alarms .....	96
Alarm Prints .....	96
Alarm Holdoffs .....	97
Propaq Audio and Visual Alarm Holdoffs .....	97
Propaq Audio Alarm Holdoff with Acuity .....	97
Setting Alarm Limits .....	98
Adjust Limits with STAT SET .....	98
Turning On and Off All Limits .....	99
Changing Individual Limits .....	99
Customizing Propaq Encore Alarm Limits .....	100
Equipment Alerts .....	101
Definitions and Indications .....	101
Responding to Equipment Alerts .....	101
Power-up Equipment Alert: Program Fault, Settings Lost .....	102
 <b>Trends .....</b>	 <b>103</b>
The Trend Status Window and Menu .....	103
How Trends are Accumulated .....	104
NIBP Trends .....	104
Displaying Trends .....	104
Selecting a Trend (NXT TRND) .....	105

<b>Printing</b>	<b>107</b>
Printing Patient Data	107
Printing Waveforms	107
Printing NIBP Measurements	109
Printing the Apnea Ticket	110
Printing When a Patient Alarm Occurs	111
OxyCRG	112
Printing Trends	113
Printing a Single Trend	113
Printing Several Trends	113
Automatic Trend Prints	114
 <b>Acuity Central Monitoring System</b>	 <b>115</b>
Intended Use	115
Connecting to the Acuity System	115
Key-press Route to Acuity Menu	117
Press NET OFF to Disconnect from Acuity	117
Printing at Acuity	118
Network Alert Message	118
 <b>Power Sources</b>	 <b>119</b>
Power Adapter Intended Use	119
Verifying Proper Power Adapter Configuration	120
Replacing the Power Adapter Fuses	120
Battery Care	123
Recharging Time	123
Monitor Functions Resumed	123
Operating Times Using Battery Power	123
Monitor Functions Based on Battery Voltage	124
Checking Battery Voltage	124
Replacing the Fuse	124
 <b>Care and Maintenance</b>	 <b>125</b>
Inspect and Clean the Monitor and Accessories	125
Maintenance	126
Service Interval Recommendations	126
Product Recycling	127
Monitor Care	127
Environmental Operating and Storage Limits	127

Extended Storage Precautions .....	127
Printer Maintenance .....	128
Loading Paper .....	128
Customer Services .....	130
Ordering and Customer Service .....	130
Technical Service .....	130
Repacking .....	130
<b>Defibrillator Synchronization .....</b>	<b>131</b>
Defib Sync Connector and Cables .....	131
Installing the Defib Sync Interface to the LIFEPAK 5 .....	131
Removing the Defib Sync Interface from the LIFEPAK 5 .....	132
Synchronized Cardioversion Using the LIFEPAK 5 .....	133
Synchronized Cardioversion Using the LIFEPAK 6s .....	134
Defib Sync Message .....	136
<b>Specifications .....</b>	<b>137</b>
ECG .....	137
Real-Time ECG Analog/Defib Sync .....	139
Impedance Pneumography (RESP) .....	140
Invasive Pressure .....	141
NIBP .....	142
Temperature .....	143
Pulse Oximetry (SpO <sub>2</sub> ) .....	144
Capnography (CO <sub>2</sub> ) .....	146
Alarms .....	150
Trends .....	151
Display .....	151
Monitor (Environmental) .....	152
Monitor (Physical) .....	153
Printer .....	154
Power .....	155
Power Adapters .....	156
Factory Default Settings .....	158
In-Service Simulated Values .....	160

**Glossary . . . . . 161**

**Index . . . . . 171**



# General Information

---

## ***Safety Summary***

This Safety Summary should be read by all Propaq Encore users. Specific warnings and cautions will be found throughout the Propaq Encore documentation where they apply.



**Caution** Caution: Federal (U.S.A.) law restricts this device to sale, distribution, or use by or on the order of a licensed medical practitioner.



**Warning** This monitor is to be operated by qualified personnel only. The operator of this monitor should read this entire manual, the monitor *Reference Guide* or *Directions For Use*, and all accessory *Directions For Use* before operating the monitor.

Place the Propaq monitor and accessories in locations where they cannot harm the patient if they fall from their shelf or mount. Lift the monitor only by its handle; do not lift it by any attached cables.

Do not connect more than one patient to a monitor. Do not connect more than one monitor to a patient.

Do not use the Propaq Encore in an MRI suite or a hyperbaric chamber.

Do not autoclave the Propaq. Autoclave accessories only if the manufacturer's instructions clearly approve it. Many accessories can be severely damaged by autoclaving.

Inspect the power adapter cord periodically for fraying or other damage, and replace the adapter as needed. Do not operate the apparatus from mains power with a damaged power adapter cord or plug.

When using a power adapter with this monitor, be sure to connect the power adapter to a three-wire, grounded, hospital-grade receptacle. Do not under any circumstances attempt to remove the grounding conductor from the power plug of the power adapter. Do not plug the power adapter into an extension cord. If there is any doubt about the integrity of the protective earth ground of the receptacle for the power adapter, do not plug in the power adapter; operate the monitor only on battery power. Contact your biomedical engineering department for assistance in identifying the proper power receptacle and making appropriate power connections.



**Warning** Make frequent electrical and visual checks on cables and electrode wires.

Avoid electrosurgery burns at monitoring sites by ensuring proper connection of the electrosurgery return circuit so that the return paths cannot be made through monitoring electrodes and probes.

During defibrillation, keep the discharge paddles away from ECG and other electrodes, as well as other conductive parts in contact with the patient. Avoid contact with any accessories connected to the Propaq's left side panel.

To ensure patient safety, the conductive parts of the ECG electrodes (including associated connectors) and other patient-applied parts should not contact other conductive parts, including earth ground, at any time.

Do not operate this product in the presence of flammable anesthetics or other flammable substance in combination with air, oxygen-enriched environments, or nitrous oxide; explosion can result.

Within certain governmental jurisdictions, all interconnected accessory equipment must be labeled by an approved testing laboratory. After interconnection with accessory equipment, risk (leakage) current and grounding requirements must be maintained.

To ensure conformance to risk (leakage) current requirements when operating from an ac mains power source, use only a Welch Allyn 503-0054 series power adapter.

This monitor should only be repaired by qualified service personnel. The operator should not attempt to open the monitor case or perform any maintenance on the monitor except for procedures explicitly described in this manual that can be performed by operators such as inspection and cleaning.

To ensure patient safety, use only accessories recommended or supplied by Welch Allyn. For a list of those accessories, see the Welch Allyn *Products and Accessories* book that accompanied this manual (P/N 810-0409-XX). Accessories must be used according to your hospital's standards and the manufacturer's recommendations. Always refer to the manufacturer's Directions for Use.



**Warning** Safe interconnection between the Propaq monitor and other devices must comply with applicable medical systems safety standards such as IEC 60601-1-1. Within certain governmental jurisdictions, all interconnected accessory equipment must be labeled by an approved testing laboratory. After interconnection with accessory equipment, risk (leakage) current and grounding requirements must be maintained.

As with all medical equipment, carefully route the patient cabling to reduce the possibility of patient entanglement or strangulation.

A product that has been dropped or severely abused should be checked by qualified service personnel to verify proper operation and acceptable risk (leakage) current values.

The pulse oximetry channel should NOT be used as an apnea monitor.

If the monitor detects an unrecoverable problem, an error message window appears containing an error number and a short message. Report such errors to Welch Allyn.

The Propaq Encore should be serviced only by a Welch Allyn service technician while under warranty. The *Propaq Encore Service Manual* (P/N 810-0696-XX) is available from Welch Allyn to assist the biomedical engineer during post-warranty period service.

## Symbols

The following symbols may appear on the Propaq Encore monitor or accessories or documentation. These internationally recognized symbols are defined by the International Electrotechnical Commission, IEC 878 and IEC 417A.














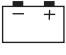



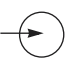


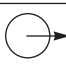

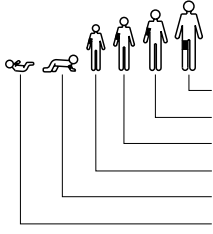

Notes identify useful tips or points about operation.



**Caution** Caution statements identify conditions or practices that could result in damage to the equipment or other property.



**Warning** Warning statements identify conditions or practices that could result in personal injury.

	Off (Standby)		Patient connections are Type CF, isolated for direct cardiac application, and protected against defibrillation.		Transformer meets requirements of a short-circuit-proof safety-isolating power transformer.
	On				Alternating current
	For continued fire protection, use only the specified fuse.		Patient connections are Type BF, and protected against defibrillation.		For indoor use only (on power adapter only)
	Direct current		Patient connections are Type B.		Caution: Refer to Reference Guide and accompanying documentation.
	Battery charging when green indicator illuminated	IPX1	Enclosure Protection Drip proof: Classification IPX1 per IEC Publication 529		The CE Mark signifies the device has met all essential requirements of European Medical Device Directive 93/42/EEC for a Class 1 product. (This symbol is on the Universal Power Adapter.)
	The CE Mark and Notified Body Registration Number signify the device has met all essential requirements of European Medical Device Directive 93/42/EEC.	 NRTL/C Evaluated to CSA 601-1 and UL2601-1	The Canadian Standards Association has evaluated this device according to CSA 601-1 and Underwriters Laboratory Standard UL 2601-1. (This symbol is on the Universal Power Adapter.)		
			Input port		Temperature sensor input
	Two way communication port		Output port		Single-use only (not reusable).
 <p>NIBP cuff sizes: Thigh Large adult Adult Small adult Child Infant</p>					Apply the NIBP cuff as shown.

## ***Propaq Encore Documentation***

### **The Documentation Set**

The Propaq Encore documentation set consists of documents for the clinician, the biomedical technician, and the department head or purchaser of accessories for the Propaq Encore monitors.

This *Propaq Encore Reference Guide* contains important safety and operating information for the clinician.

The *Propaq Encore Service Manual* (P/N 810-0696-XX) contains information on how to properly maintain the Propaq Encore through routine calibration, inspection, and maintenance.

The Welch Allyn *Products & Accessories* booklet (P/N 810-0409-XX) provides a comprehensive list of accessories recommended for Propaq Encore monitors and options.

### **About This Reference Guide**

This Reference Guide provides descriptions and operating information for the Propaq Encore models 202EL, 204EL, and 206EL, including all available options at the time of this manual's printing.

### **Statement of Expectations of the Reader**

This Reference Guide was written for the clinician. Although this guide may describe some monitoring techniques, Welch Allyn expects that you are a trained clinician who knows how to take and interpret a patient's vital signs. The Propaq Encore has been designed as a quality monitor; however, inherent limitations require that good clinical judgment always prevails.



# Getting Started

---

## *Introducing the Propaq Encore*

### Intended Use

Before using the Propaq Encore on a patient, be sure you understand the Safety Summary at the front of this book. It provides important information about safely using the Propaq Encore. The Propaq Encore monitor is intended to be used by skilled clinicians for multiparameter vital signs monitoring of neonatal, pediatric, and adult patients in health care facility bedside applications; as well as for intra- and interfacility transport.

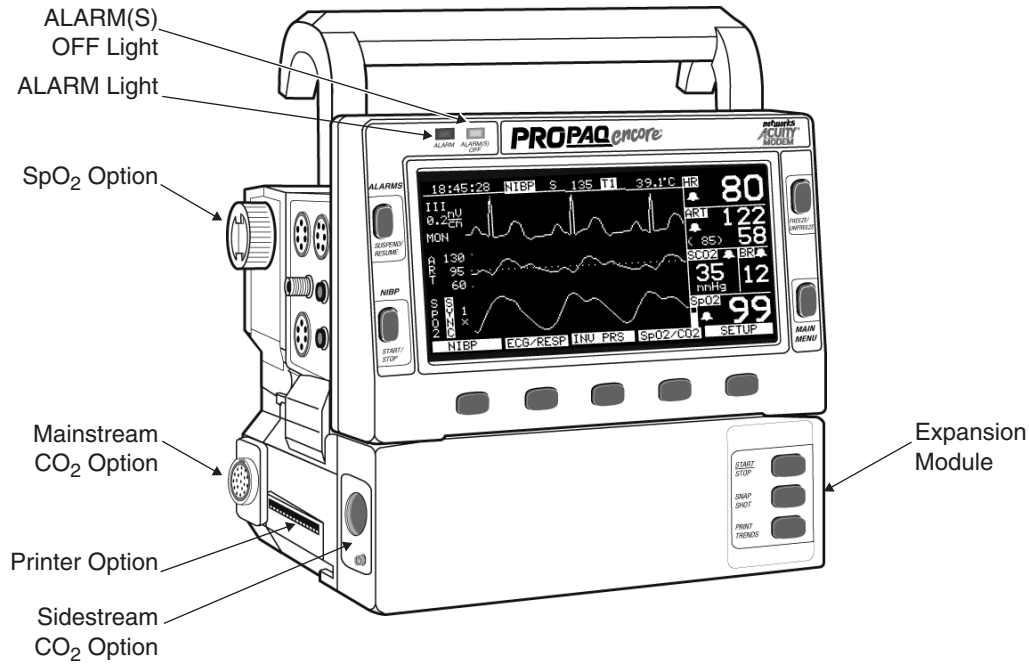
### Propaq Encore Models and Options

Three models of Propaq Encore monitors are available.

<b>Features common to all models</b>	ECG, 3-lead or 5-lead configurations, 0.05-40/0.5-40 Hz NIBP, with neonatal, pediatric and adult modes Temperature, 2 channels: YSI™ 400 and 700 series-compatible connectors Defibrillator Synchronization Real-time Analog output of ECG Electrocautery noise suppression on all channels except Impedance Pneumography
<b>206EL</b>	Two Invasive Pressure Channels
<b>204EL</b>	One Invasive Pressure Channel
<b>202EL</b>	No Invasive Pressure
<b>Options available for each model</b>	Pulse Oximetry (SpO <sub>2</sub> ) Capnography (CO <sub>2</sub> ) (available only with SpO <sub>2</sub> ): Mainstream Capnography (MCO <sub>2</sub> ) Sidestream Capnography (SCO <sub>2</sub> ) Dualstream Capnography (Both MCO <sub>2</sub> and SCO <sub>2</sub> ) Impedance Pneumography (RESP) (available only with SpO <sub>2</sub> ) Printer HP-compatible side panel

## Expansion Module

The Propaq Encore Expansion Module attaches to the monitor and houses additional capabilities. The Expansion Module can be fitted with the SpO<sub>2</sub>, CO<sub>2</sub>, and Printer options.



## Propaq Encore Pulse Oximetry Option (SpO<sub>2</sub>)

The Propaq Encore Pulse Oximetry option (SpO<sub>2</sub>) can be installed in the Expansion Module or in a smaller unit that attaches to the rear of the monitor:

- Masimo<sup>®</sup> Pulse Oximetry option (motion tolerant).
- Nellcor<sup>®</sup> Pulse Oximetry option (motion tolerant).
- Nellcor<sup>®</sup> Pulse Oximetry option (without motion tolerance).

## Capnography (CO<sub>2</sub>) Options

The Propaq Encore CO<sub>2</sub> options allow carbon-dioxide monitoring. The Mainstream CO<sub>2</sub> option and Sidestream CO<sub>2</sub> option allow CO<sub>2</sub> monitoring directly in the breathing circuit of a ventilator. The Sidestream CO<sub>2</sub> option also allows CO<sub>2</sub> monitoring of non-intubated patients through a cannula. The CO<sub>2</sub> options can be installed separately, or together as Dualstream CO<sub>2</sub> in the Expansion Module. These options require the Pulse Oximetry (SpO<sub>2</sub>) option.



## Impedance Pneumography (RESP) Option

The RESP option detects the rate or absence of respiratory effort, and is configured with the Pulse Oximetry option.

## Printer Option

The Expansion Module with Printer (EMP) provides a lightweight 3-channel recorder.

## Propaq-to-Acuity Option

This option allows communication between the Propaq Encore and the Acuity Central Monitoring System by means of an ethernet network system installed in your facility. The Acuity System operator can view the patient data and control most of the bedside Propaq functions. The Propaq Encore connects to the Acuity System through an Acuity network cable that plugs into the Propaq right side panel.

## Modem-Propaq Option

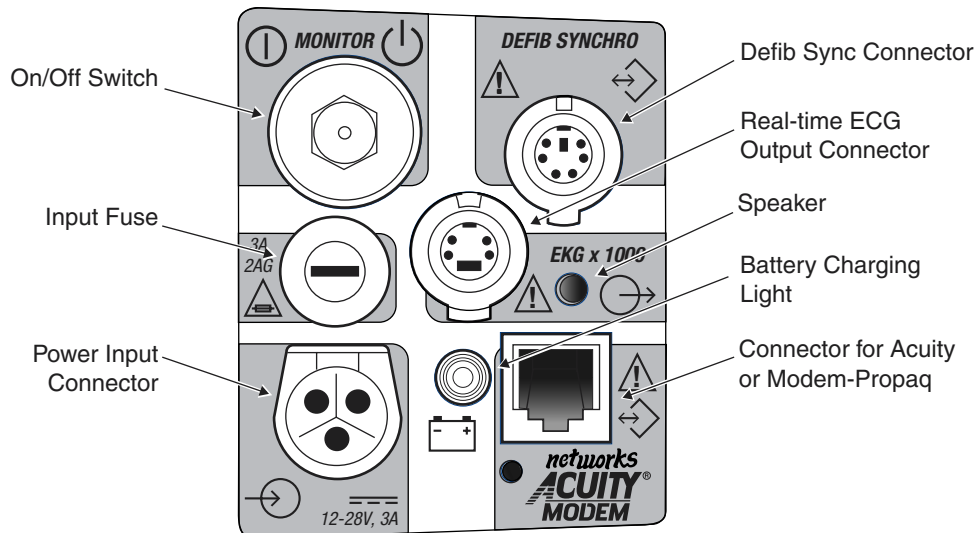
This option allows telecommunication between a Propaq Encore and the Acuity System by means of external modems. This option is configured with the Propaq-to-Acuity option. For more information refer to the *Modem-Propaq Reference Guide*.

## HP-compatible Side Panel Option

The HP connector-compatible option makes the Propaq Encore compatible with many Hewlett-Packard sensors and accessories used with the Hewlett Packard Component Monitoring System. This option replaces the standard Propaq Encore left side panel.

## Using the Propaq Encore

### System Controls (Right Side Panel)



**Warning** Safe interconnection between the Propaq Encore and other devices must comply with applicable medical systems safety standards such as IEC 601-1-1.

#### On/Off Switch

This switch turns the monitor on and off. The switch is recessed to prevent accidentally turning off the monitor, which would result in losing patient data.

#### Input Fuse

The input fuse, which protects the Propaq Encore against power surges, is a 3 Ampere fuse, externally replaceable by qualified service personnel. See page 124 for fuse replacement instructions.

#### Power Input Connector

This receptacle accepts the Welch Allyn ac power adapter, which must be used for ac mains operation and battery charging. The Propaq Encore is also designed to operate with other 12-28 volt, dc-only power sources, such as a vehicle battery system.

### ***Defib Sync Connector***

This connector allows connection with a LIFEPAK 5 or LIFEPAK 6s defibrillator for synchronized cardioversion. See page 131 for more information.

### ***Real-Time ECG Output Connector***

This connector provides a real-time analog ECG signal output.

### ***Battery Charging Light***

This green light turns on when a power source (ac power adapter or external dc source) is connected and the battery is charging. Although the monitor may be turned off, battery charging continues when an external power source is connected.

### ***Connector for Acuity or Modem-Propaq***

This connector allows either direct connection to an Acuity System, or connection to an external modem for telecommunication to an Acuity System. For more information about the Acuity System, see page 115. For more information about the Modem-Propaq, refer to the *Modem-Propaq Reference Guide*.

## **Alarm Lights**

Alarms and limits are described in detail beginning on page 91.

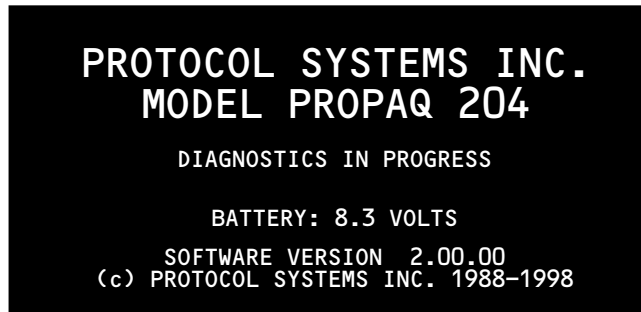
### ***ALARM Light***

When an alarm limit is violated, the red **ALARM** light turns on.

### ***ALARM(S) OFF Light***


When any alarm limit is turned off, the yellow **ALARM(S) OFF** light turns on.

## Power-up Screen



When you first turn on the monitor, the power-up screen displays information about the Propaq Encore and the monitor runs diagnostic tests to ensure proper functioning.

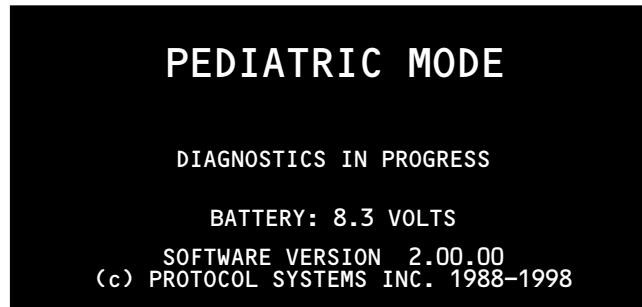A few seconds later, the top two lines of the screen are replaced with text indicating the current patient mode (adult, pediatric, or neonatal).



**Warning** Before you use a Propaq on a new patient, always turn off the Propaq for a few seconds, then turn it on again. This clears the prior patient's trend values, alarm limit settings, and NIBP cuff inflation target.

1. If the Propaq Encore has been used for a previous patient, switch the monitor off, then on again. The monitor will turn on in the powerup patient mode with the associated settings.



Verify that the powerup tone is produced. If the monitor has SpO<sub>2</sub>, verify two tones are produced to make sure that both speakers are working.

2. Verify the monitor is in the correct patient mode according to the patient's age. If the patient mode is not correct, change it (see page 29 to change the patient mode).



**Warning** Always check the patient mode when monitoring a new patient. The patient mode determines default alarm limits, maximum cuff inflation pressure, and internal algorithm settings.

3. Verify the battery voltage is sufficient for monitoring. If it is less than 7.4 V, connect to a power adapter (see page 119 for information about the power adapter).

## ***Powerup Equipment Alert: Program Fault, Settings Lost***

If a PROGRAM FAULT: SETTINGS LOST, TIME/DAY RESET equipment alert appears when you turn on the monitor, the monitor cannot recall the programmed custom settings and current time and date. This can occur if the battery is drained or after new software has been installed.

If this occurs, the monitor provides a special sequence of display windows to help you regain use of your monitor as quickly as possible. Do the following:

1. Connect an ac power adapter to recharge the battery (if the battery is drained).
2. Press any button below the equipment alert screen to acknowledge the alert. The monitor will display the Mode Setup window (shown on page 33).
3. Press these buttons to select one of the Factory patient modes for use:

<b>Factory Adult mode</b>	<b>POWERUP*, YES.</b>
<b>Factory Pediatric mode</b>	<b>NEXT, POWERUP*, YES.</b>
<b>Factory Neonatal mode</b>	<b>NEXT, NEXT, POWERUP*, YES.</b>

After you press **YES**, the monitor will display the Time/Day window.

4. Press **NEXT**, **UP**, and **DOWN** as needed to set the time and date. Then press **ENTER** to store the new time and date.



These display screens are only displayed in this order if the PROGRAM FAULT equipment alert occurs.

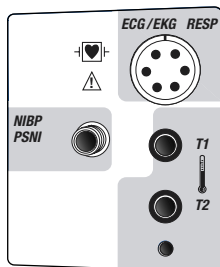
5. Turn off the monitor, then turn it on again so the settings will take effect.

The monitor is ready for use. If you want to store some customized patient mode program settings, refer to page 33.


If you follow these steps and the equipment alert reappears at powerup, the monitor may need to be serviced and the battery replaced. Contact a qualified service person.

## Patient Connections


The left side panels differ depending on the Propaq Encore model. All models have ECG, NIBP, and two temperature connectors. The Propaq Encore 204 left side panel includes one invasive-pressure connector, and the Propaq Encore 206 includes two invasive-pressure connectors.




Propaq Encore 202 EL

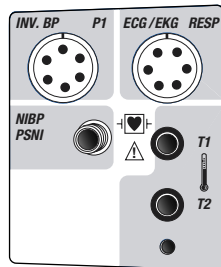


Propaq Encore 204 EL

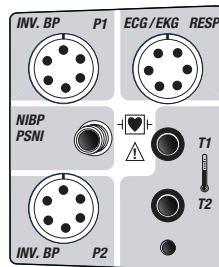


Propaq Encore 206 EL

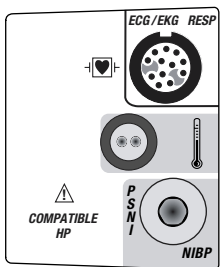

202 HP

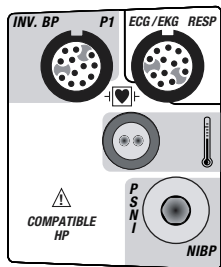


204 HP

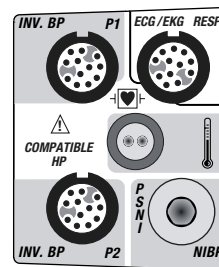


206 HP

On Propaqs with the Hewlett-Packard connector option, all models have only one temperature connector, the YSI 400 connector.


## Option Connectors

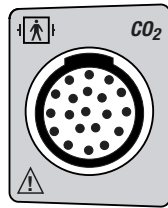

Mainstream CO<sub>2</sub> Connector

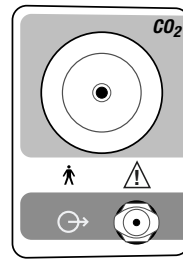


Sidestream CO<sub>2</sub> Connector

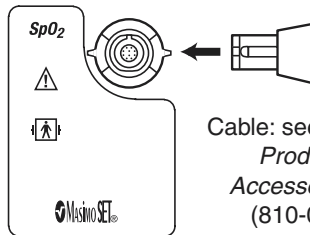


Cable: see Welch Allyn  
Products and  
Accessories Guide  
(810-0409-XX)

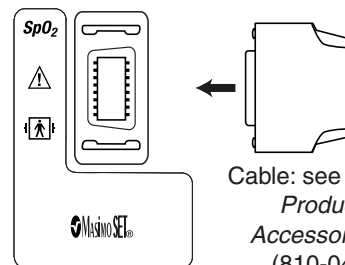Masimo SpO<sub>2</sub> Connector  
(motion tolerant)




Cable: see Welch Allyn  
Products and  
Accessories Guide  
(810-0409-XX)

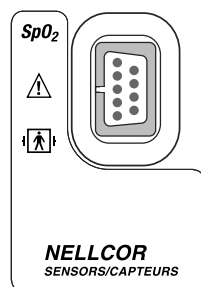Masimo SpO<sub>2</sub> Connector  
(motion tolerant)



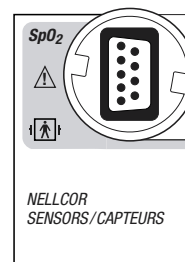Nellcor SpO<sub>2</sub> Connector  
(newer style, motion tolerant)



Nellcor SpO<sub>2</sub> Connector  
(older style, without motion  
tolerance)

## Propaq Encore Display

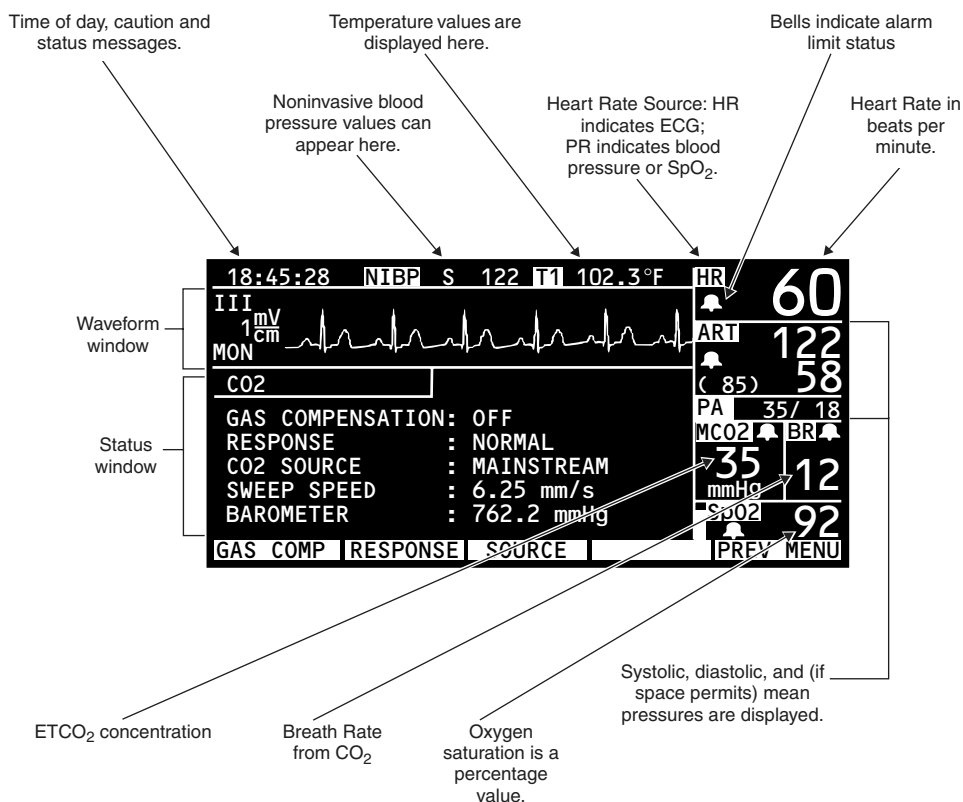
The display shows waveforms, vital sign numeric values, Propaq Encore status, and alarm information in different windows. Different vital sign numeric values (such as heart rate and blood pressures) have upper and lower range limits. If the Propaq Encore detects a vital sign value outside of the Propaq's measurable range, the monitor displays --- (below the range) or + + + (above the range) instead of the vital sign value.



**Warning** The Propaq Encore will show + + + for HR numerics between 301-350 beats per minute. Above 350 beats per minute, it may display incorrectly low heart rates, due to intermittent picking of R-waves.



Due to differences in software versions and standards required by different countries, the displays shown in this reference guide may be slightly different than the display on your Propaq Encore.




The screenspace is reallocated when vital signs are added or removed. By changing the size of the numeric windows below the heart rate, the Propaq Encore provides the best possible view of all numerics for vitals signs being monitored.



## Using the Propaq Encore

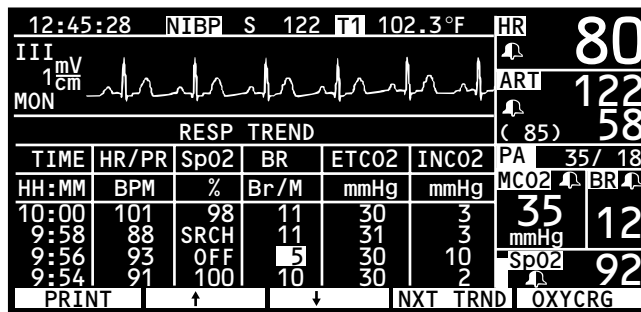You can select up to three waveforms to be shown on the Propaq. When only one waveform is selected, a trend window automatically appears below the waveform. While changing Propaq Encore settings, a status window may appear below the waveform.




Patient waveform and trend information can be simultaneously displayed, while numeric values are continuously updated.

## Propaq Encore Buttons

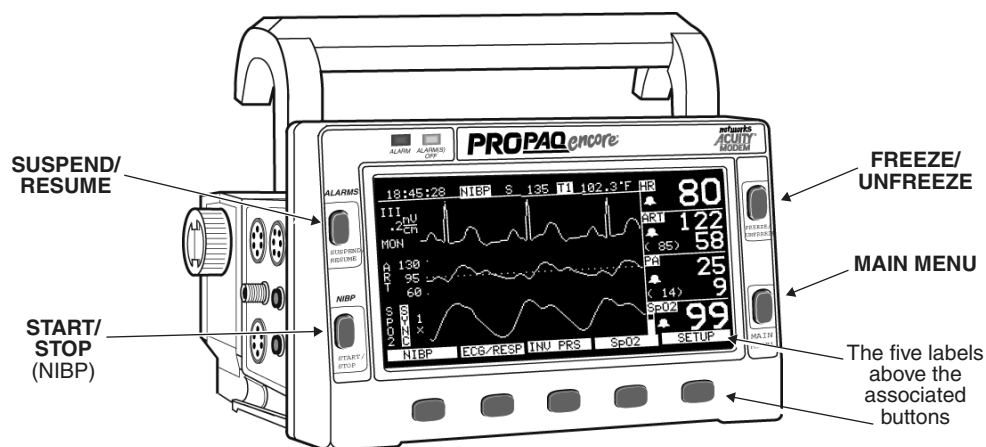The four buttons at the sides of the screen are reserved for the most commonly used functions.

- SUSPEND/RESUME** Suspends or resumes alarm tone.
- START/STOP** Starts and stops NIBP measurements. The **STOP** function will automatically vent the cuff.
- FREEZE/UNFREEZE** Freezes or “unfreezes” the waveforms. If only one or two waveforms are displayed and you press **FREEZE**, the frozen waveform(s) are shown along with an active waveform so you can continue to monitor the patient’s condition.
- MAIN MENU** Pressing **MAIN MENU** always returns the monitor to the top level menu.




The five buttons below the screen, and their associated labels located on the screen, provide access to the menus.


Later in this manual, the notation **A**, **B**, **C** is used as a shorter way to say “Press Button **A**, then **B**, then **C**.”

## Propaq Encore Menus

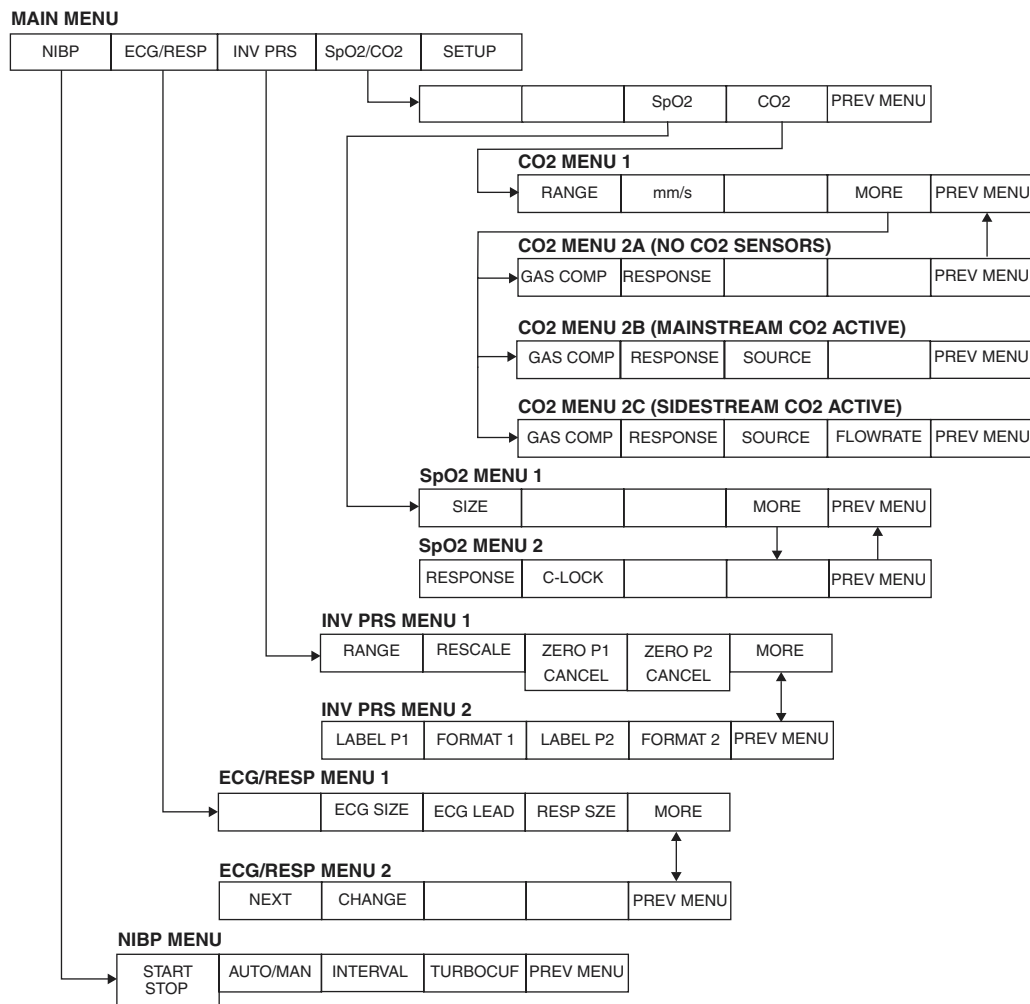
Menus for some patient vital signs are displayed only if that option is included in your Propaq.



## Key-press Route to Setup Menu 1

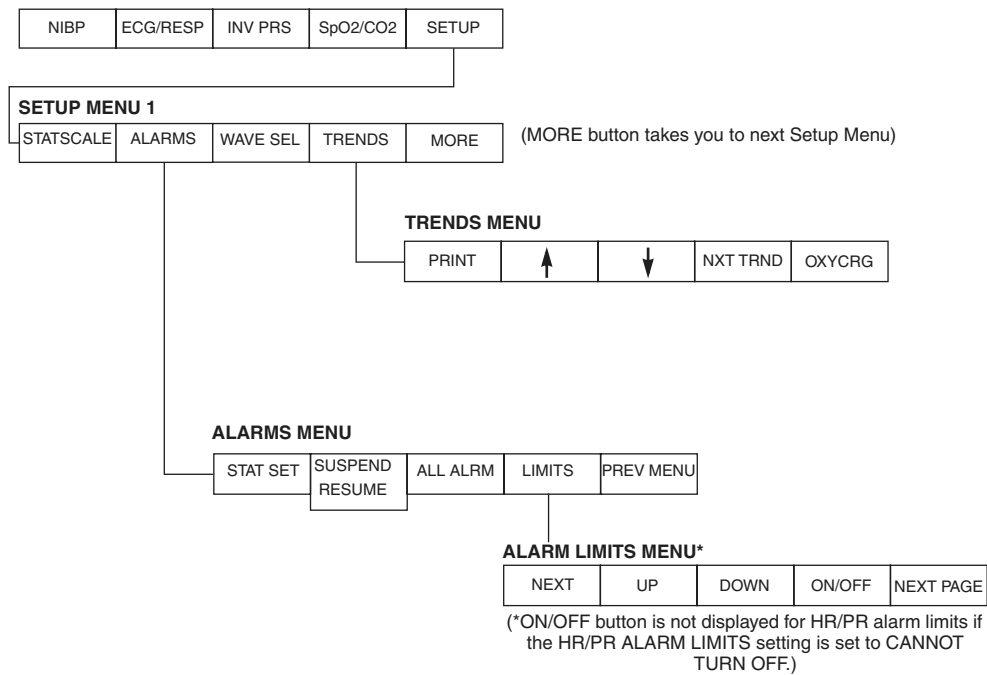

## Key-press Route to Setup Menu 2

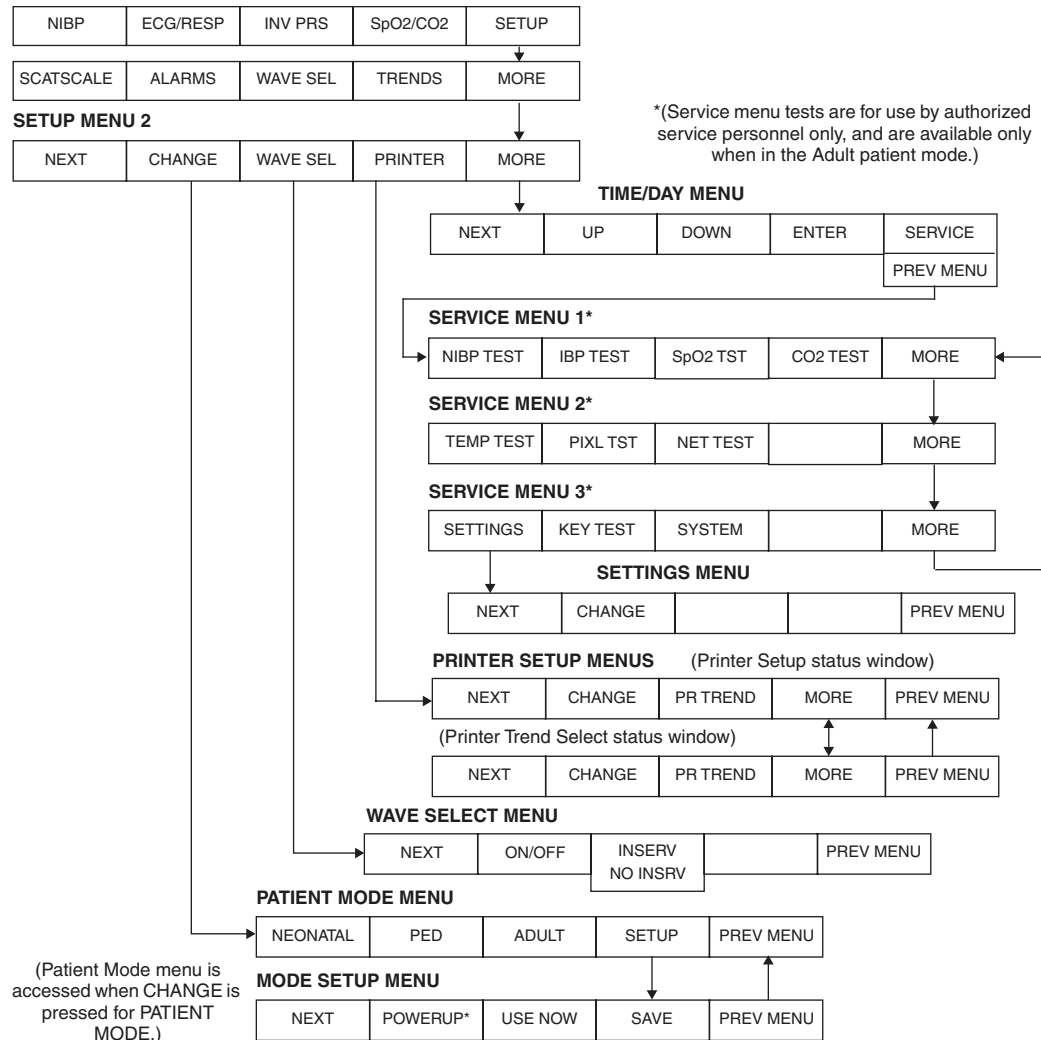

## Monitor Setup

Setup Menu 1 is accessed by pressing the **SETUP** button on the Main Menu.

STATSCALE	ALARMS	WAVE SEL	TRENDS	MORE
-----------	--------	----------	--------	------

<b>STATSCALE</b>	Automatically readjusts all waveform scales.
<b>ALARMS</b>	Allows access to the Alarms menu.
<b>WAVE SEL</b>	Allows you to turn on and off desired waveforms or NIBP numerics for display.
<b>TRENDS</b>	Allows access to the Trend settings and display.
<b>MORE</b>	Displays the next setup menu and the following status window:

SETUP	HR/PR	RR/BR	( 85)	58
CURRENT SOURCE : ECG	MC02		PA	35/ 18
SELECTED SOURCE : ART	---		MC02	BR
SWEEP (mm/s) : 12.5	3.13		35	12
ALARM TONE : HIGH			mmHg	
HR/PR TONE : LOW			SpO2	92
PATIENT MODE : ADULT				
NEXT	CHANGE	WAVE SEL	PRINTER	MORE

<b>NEXT</b>	Selects the next setting in the status window.
<b>CHANGE</b>	Changes the currently selected display setting. (Pressing <b>CHANGE</b> at PATIENT MODE allows you to choose between Adult, Pediatric, and Neonatal in a Patient Mode window.)
<b>PRINTER</b>	Allows access to the Printer Menu.
<b>MORE</b>	Allows access to the Time/Day window.
<b>CURRENT SOURCE</b>	When the selected HR/PR source is no longer available, the current source is the active source with highest priority. The RR/BR source cannot be manually selected. It will always be CO <sub>2</sub> if CO <sub>2</sub> is active. Otherwise, it will be ECG/RESP.
<b>SELECTED SOURCE</b>	The user-selected HR/PR source is displayed along with the HR/PR source currently being used by the monitor.
<b>SWEEP (mm/s)</b>	The selectable sweep speeds for HR/PR are 12.5, 25, and 50 mm/sec. The sweep speeds for RR/BR are 3.13, 6.25, and 12.5 mm/sec.
<b>ALARM TONE</b>	Sets the Alarm Tone volume to HIGH, MEDIUM, or LOW.

HR/PR TONE

Sets the Heart Tone volume to HIGH, MEDIUM, LOW, or OFF.

PATIENT MODE

Pressing **CHANGE** in this selection displays the following Patient Mode window:

PATIENT MODE		BATTERY: 9.2 VOLTS		( 85 )	58
SELECT PATIENT MODE BASED ON AGE:					
NEO	:	< 44	WEEKS	GEST. AGE	
PED	:	> 44	WEEKS	GEST. AGE,	
	:	< 9	YEARS		
ADULT:	:	> 9	YEARS		
NEONATAL		PED		ADULT	SETUP
				PREV	MENU

If you press **NEONATAL**, **PED**, or **ADULT**, a confirmation window appears, requiring you to confirm your selection:

PATIENT MODE				( 85 )	58
CHANGING SETTINGS AND MODE TO:					
CUSTOM: ADULT					
ADULT ALARM LIMITS WILL BE SET					
<b>ARE YOU SURE?</b>					
YES		NO			



Whenever you change the patient mode, the alarm limit settings are automatically changed to the defaults for that mode. If Custom settings have been set for that mode, the defaults are the Custom mode settings. If no Custom settings have been set, the defaults are the Factory Mode settings. See page 33 for more information about patient modes.

If you change the patient mode, the CO<sub>2</sub> alarm limits in the new mode might vary slightly from the originally-programmed CO<sub>2</sub> limits for the new mode. Check the CO<sub>2</sub> alarm limits.

If you press **SETUP** in the previous Patient Mode window, the Mode Setup window appears. This allows you to set custom patient modes and powerup defaults as described on page 33.

## Selecting Waveforms for Display

To select waveforms for display, press **SETUP, WAVE SEL**. Use the **NEXT** and **ON/OFF** buttons to turn on the desired waveforms in the wave select window:

SETUP		WAVE SELECT		( 85 )	58
ECG	: ON	RESP	: OFF	PA	35/ 18
ART	: <b>ON</b>	SpO2	: ON	MC02	BR
PA	: OFF	NIBP	: ON	35	12
CO2	: OFF			mmHg	
				SpO2	92
NEXT	ON/OFF	INSERV		PREV	MENU

### Display Priorities

You can turn on more than three waveforms, but only the first three waveforms listed in the wave select window that are monitored are displayed. The patient parameters being monitored are listed in the order they will be displayed if all are turned on.

Because of the critical nature of the ECG waveform, you cannot turn off ECG. However, if ECG is not monitored, another waveform will occupy its place.

The displayed waveforms are also the ones printed if a printer is attached.

## Setting the Time and Date

To set the time and date, from the Main Menu press **SETUP, MORE, MORE**. The monitor displays the Time/Day window:

TIME/DAY		BATTERY: 9.2 VOLTS		( 85 )	58
<u>TIME</u>		<u>DAY</u>		PA	35/ 18
H:MIN:S		MO/DA/YR		MC02	BR
07:45:32		06/12/97		35	12
				mmHg	
				SpO2	92
NEXT	UP	DOWN	ENTER	PREV	MENU

Press **NEXT**, **UP**, and **DOWN** as needed to set the time and date. Then press **ENTER** to store the new time and date.

## Time/Day Settings and Trends



**Warning** Changing the hour/minute/second setting for the monitor in the Time/Day window can cause the monitor to erase previously stored patient trend data.

When you change the hour/minute/second setting for the monitor in the Time/Day window, the monitor deletes any patient trend data that is older than five hours for non-NIBP trends or older than eight hours for NIBP trends according to the new clock setting.

However, if the monitor has not yet stored the full capacity of trends and you change the hour/minute/second setting to a time that is within the stored trend period, previously stored trends are not erased.

Changing the day, month, or year setting does not affect the stored patient trends.

## Changing the Date Format, Filter, and Units

To change the date format, ECG filter, or some measurement units, first make sure you are in the Adult patient mode. Then press **SETUP, MORE, MORE, SERVICE, YES** (to access the Service Menu), **MORE, MORE, SETTINGS**. The monitor displays the Settings window:

SETTINGS	BATTERY: 9.2 VOLTS	( 85)	58
DATE	: MO/DA/YR	PA	35/ 18
FILTER	: 60 Hz	MC02	BR
TEMP F/C	: °C	35	12
DECIMAL	: .	mmHg	
HR/PR ALARM LIMITS:	CAN TURN OFF	SpO2	92
CO2 UNITS:	mmHg		
NEXT	CHANGE		PREV MENU

NEXT	Selects the next setting in the status window.
CHANGE	Changes the currently selected display setting.
DATE	Sets the date format: Month/Day/Year, Day.Month.Year, or Year/Month/Day.
FILTER	Sets the ECG filter frequency. Make sure it is set to your ac mains frequency.
TEMP F/C	Sets the temperature display units: either degrees Fahrenheit or Celsius. If you change the units, the TEMP trends will <b>not</b> be cleared.



## Monitor Setup

DECIMAL	Sets the decimal character as either a period (.) or a comma (,).
HR/PR ALARM LIMITS	Allows or prohibits turning off the HR/PR alarm limits. If CANNOT TURN OFF is selected, the <b>ON/OFF</b> button is not displayed on the HR/PR Alarm Limits Menu.
CO <sub>2</sub> UNITS	Sets the CO <sub>2</sub> display units as mmHg, kPa, or percent (%). If you change the units, the CO <sub>2</sub> trends will be cleared and CO <sub>2</sub> alarm limit settings change to the factory default settings for the currently-used patient mode.



Any time you change the Date, Filter, Temp F/C, Decimal, HR/PR Alarm Limits (CAN or CANNOT TURN OFF) or CO<sub>2</sub> Units setting, the new setting also becomes the powerup default setting.

## Setting the Current, Custom, and Powerup Modes

The Propaq Encore has two sets of patient mode settings:

- Factory patient modes. The powerup settings and alarm limits for these patient modes are preset and cannot be changed. They are listed in Appendix B on page 158.
- Custom patient modes. You can customize the powerup settings and alarm limits for these patient modes. (See **SAVE** on page 34.)



**Warning** If any alarms are set to OFF and you select **SAVE** to store the settings for that CUSTOM patient mode, those alarms will be OFF whenever the Propaq powers up in that CUSTOM patient mode or when that CUSTOM patient mode is selected. Consider carefully before setting CUSTOM patient mode powerup alarms to OFF.



The alarm for apnea cannot be turned off at any time.

You can program the Propaq Encore to power up in any of the Factory patient modes or the Custom patient modes. You can also change the current patient mode during operation.



Whenever you change the patient mode, the alarm limit settings automatically change to the settings for that mode.

From the Main Menu, press **SETUP, MORE, CHANGE, SETUP**. The Mode Setup window appears:

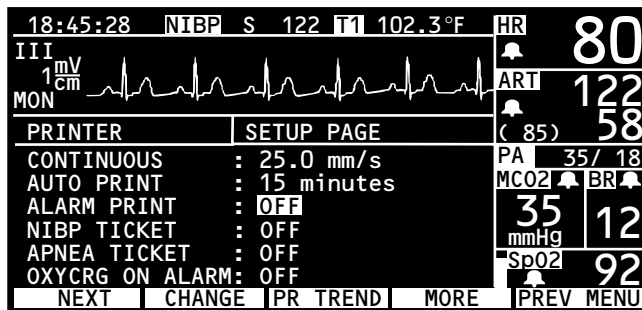
MODE SETUP		BATTERY: 9.3 VOLTS		ZEROED	
FACTORY:	*ADULT	PED	NEO	P2 NO ZERO	
CUSTOM :	ADULT	PED	NEO	MC02 BR	
				38	12
				mmHg	
				SpO2	97
* = PATIENT MODE ON POWERUP					
NEXT	POWERUP*	USE NOW	SAVE	PREV	MENU

The asterisk (\*) indicates which patient mode is currently selected for powerup.

- |                 |  |
|-----------------|--|
| <b>NEXT</b>     | Selects the next setting in the status window.   |
| <b>POWERUP*</b> | Selects the highlighted patient mode (and its associated settings) as the powerup mode. The selected powerup mode is marked by an asterisk (*). (This does not change the current patient mode.)   |
| <b>USE NOW</b>  | Selects the highlighted patient mode (and its associated settings) as the current patient mode. (This does not affect the powerup mode.)   |
| <b>SAVE</b>     | <p>Use this button to reprogram the settings of Custom patient mode:</p> <ol style="list-style-type: none"> <li>1. Make sure the patient mode you want to reprogram (<b>ADULT, PED</b> or <b>NEO</b>) is currently-used (to change patient modes, highlight the desired mode and press <b>USE NOW, YES</b>).</li> <li>2. Exit the Mode Setup window, then use other menus and buttons to set the monitor settings and alarm limits as desired.</li> <li>3. Re-enter the Mode Setup window, highlight the desired Custom mode, and press <b>SAVE, YES</b>.</li> </ol> |

## Printer Functions

Press **SETUP**, **MORE**, **PRINTER** to display the printer menu and setup window.



<b>NEXT</b>	Selects the next setting in the status window.
<b>CHANGE</b>	Changes the currently selected display setting.
<b>PR TREND</b>	Prints all trends turned on in the Printer Trend Select Window.
<b>MORE</b>	Pressing the <b>MORE</b> button displays another menu and status window.
<b>PREV MENU</b>	Returns you to the previous menu.
<b>CONTINUOUS</b>	Sets the print speed for real time (continuous) measurements to 6.25, 12.5, or 25 mm/sec. This sets the print speed for a printout obtained by pressing the <b>START/STOP</b> button on the printer.
<b>AUTO PRINT</b>	Automatically prints 8 seconds of patient information every 15 minutes, 30 minutes, 1 hour, 2 hours, or 4 hours. This is the latest patient information (real time). The print speed is automatically set to 25 mm/sec.
<b>ALARM PRINT</b>	Automatically prints upon an alarm. The Propaq Encore prints 20 seconds of patient information. The first 12 seconds contain information prior to the alarm. The print speed is automatically set to 25 mm/sec.
<b>NIBP TICKET</b>	Automatically prints an NIBP Ticket when the measurement is taken.
<b>APNEA TICKET</b>	When turned on, an Apnea Ticket is printed at the conclusion of an apnea alarm and at the one-minute clock interval if the apnea alarm does not cease.
<b>OXYCRG ON ALARM</b>	When turned on, an oxycardiogram will print if an HR/PR, SpO <sub>2</sub> , or RR/BR alarm occurs. For more information on OxyCRG, see page 112.

## PRINTER FAULT Messages

These PRINTER FAULT messages can appear in an equipment alert window.

### LOW BATTERY, PRINTER DISABLED

This message appears when the Propaq's battery voltage is less than 7.6 volts. To continue operation, plug the ac power adapter into the Propaq.

### CHECK DOOR

The door on the bottom of the printer is open. Close door to remove this message.

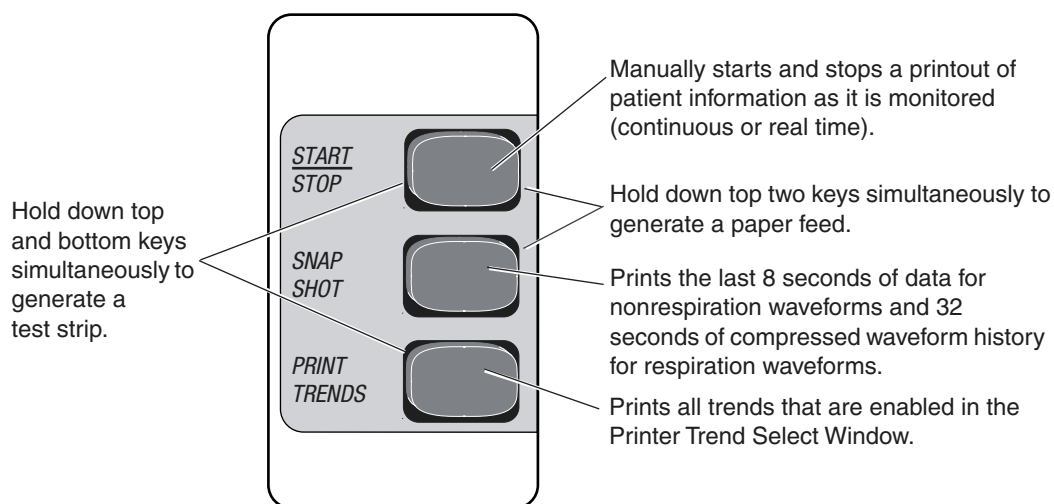
### PAPER OUT

To add printer paper, see page 128.

### OVERHEATING

The printer is overheating. Service may be required.

The front panel of the printer lets you control the basic printer functions.



If you press **FREEZE** prior to pressing **SNAPSHOT**, the printer prints the 8 seconds of patient information obtained prior to when you pressed **FREEZE**.

## *Learning the Propaq Encore*

### Using In-Service Mode

You can practice using the Propaq Encore without a patient simulator by using the Propaq's in-service mode of operation. The in-service mode cannot be activated while you are monitoring a patient. The message "SIMULATING" alternates with the time of day and patient mode on the display.

To begin practicing with your Propaq, disconnect all patient cables connected to the monitor. Leave the cuff connected so you can take NIBP measurements. If you have been monitoring a patient, turn off the Propaq Encore and turn it back on. From the Main Menu, press **SETUP**, **WAVE SEL**, **INSERV**.

The Propaq Encore has two sets of simulated patient information—an initial set and an alternate set. To change between them, press the **INSERV** button again.

If you connect a patient cable or set the NIBP channel to automatically take pressure measurements, the Propaq Encore stops simulating, goes through its powerup tests, and erases any simulated trend data it might have stored.

### What You Can Do With In-Service Mode

While using the in-service mode, you can press any of the Propaq Encore buttons, except for the **AUTO/MAN** button in the NIBP Menu, to change a function setting. You can also:

- change the ECG and RESP waveform sizes
- set alarm limits and cancel alarms
- STAT SET alarms
- customize the Propaq Encore settings
- change from °F to °C
- simulate invasive-pressure zeroing

## ***NIBP***

For noninvasive pressure measurements, keep the Propaq Encore in manual NIBP operating mode and take pressure measurements by pressing the **START** button. You can also press the NIBP Menu's **TURBOCUF** button to consecutively take pressure measurements for five minutes.

## ***Printer Message***

Simulated data can be printed on the Propaq Encore Printer. All printouts include the message "SIMULATED DATA" every four inches to prevent simulated data from being mistaken for actual patient data.

## **What You Cannot Do With In-Service Mode**

- You cannot use in-service mode to calibrate the monitor.
- You cannot set the Propaq Encore to take automatic noninvasive pressure measurements (except Turbocuf) while using in-service mode.
- You cannot use Defib Sync or Real-time ECG output while using in-service mode.
- You cannot activate in-service mode if you have been monitoring a patient.

# Patient Monitoring

---

## *ECG/RESP*

### Intended Use—Impedance Pneumography (RESP)

The Respiration channel is intended to detect the rate or absence of respiratory effort, deriving the signal by measuring the AC impedance between the selected terminals of the ECG electrodes. RESP displays a respiration rate and waveforms. Two respiration lead selections are available, Lead 1 (RA-LA) and Lead 2 (RA-LL).



**Warning** Impedance pneumography detects respiratory effort via changes in chest volume; therefore, impedance pneumography can be used to detect central apnea. However, apnea episodes with continued respiratory effort, such as obstructive apnea and mixed apnea, may go undetected. Always monitor and set alarms for SpO<sub>2</sub> when using impedance pneumography to monitor respiratory function.

With any monitor that detects respiratory effort via impedance pneumography, artifact due to patient motion, apnea mattress shaking, or electrocautery use may cause apnea episodes to go undetected. Always monitor and set alarms for SpO<sub>2</sub> when using impedance pneumography to monitor respiratory function.

The Propaq Encore automatically rejects cardiovascular artifact (CVA). This function is dependent upon accurate ECG R-wave detection. Therefore, always select the ECG lead with the most prominent QRS complex when monitoring respiration via impedance pneumography.

Don't place the Propaq Encore monitor with RESP in close proximity with another respiration monitor because the RESP measurement frequencies may interfere with one another.

Because pacemaker pulses in some instances may be falsely counted as breaths, impedance pneumography is not recommended for use on paced patients.



Impedance pneumography is not recommended for use with high frequency ventilation.

Since RESP is derived from the same leads as the ECG channel, the Propaq Encore determines which signals are cardiovascular artifact and which signals are a result of respiratory effort. If the breath rate is within five percent of the heart rate or a multiple or sub-multiple of the heart rate, the monitor may ignore breaths and trigger an apnea alarm.

## Intended Use–ECG

The Propaq Encore is intended for ECG monitoring of either a five-lead or three-lead configuration, including the Marriott configuration 1 (MCL1 requires all three electrodes). The five lead configuration can derive one of seven user-selected signals, Lead I, II, III, aVR, aVL, aVF, or V.

The monitor will automatically determine if only three lead wires are connected, and will automatically reduce the number of selectable leads to three (I, II, III). If four-wire ECG cables are used, they will be handled as if they were three-wire cables.



**Warning** The Propaq Encore monitor does not have automated arrhythmia analysis, therefore, some ventricular tachycardias and ventricular fibrillation may not be interpreted correctly and may display an inaccurate heart rate.

The Propaq Encore 200 series does not have automated ST segment monitoring, although with ECG set for extended bandwidth, ST segments may be accurately displayed and printed.

The Propaq Encore ECG's bandwidth is 0.5-40 Hz in Monitor Mode and 0.05-40 Hz in Extended Mode. Monitor Mode is useful to minimize baseline wander due to respiration or other artifact. However, in Monitor Mode, ST segments can be distorted, potentially causing underestimation of ST elevation and overestimation of ST depression. Always use Extended Mode when observing ST segment morphology on the display or printer.

The Propaq Encore can be used during procedures using electrosurgical machines and defibrillators. However, even though the ECG channel contains electrosurgical interference suppression (ESIS) circuitry, noise artifact may be displayed on the ECG trace while an electrosurgical device is in use. This will vary depending on ECG electrode placement and the operative site.





**Warning** High-intensity radio frequency (RF) energy from external sources, such as an improperly connected electrosurgical unit, can induce heat into electrodes and cables which can cause burns on the patient. Reading errors and damage to equipment may also result. This hazard can be reduced by (1) avoiding the use of small ECG electrodes, (2) selecting ECG electrode attachment points remote from the surgical site and from the electrosurgical return electrode, (3) using electrosurgical return electrodes with the largest practical contact area, and (4) assuring proper application of the electrosurgical return electrode to the patient.

Verify patient mode. Incorrect patient mode may result in inaccurate heart rates and inappropriate alarm settings.

Even though the Propaq Encore contains fully isolated patient-connected circuitry, it has not been specially designed for direct cardiac application.

The Propaq Encore can be used on patients with pacemakers. See page 47.

## ECG Connector and Applicable Accessories

Use only with accessories provided or recommended by Welch Allyn. Refer to the *Welch Allyn Products and Accessories* booklet.

To prevent injury, use the provided garment clips to route the ECG cables away from the patient's head.



**Warning** Use of ECG cables with loose or faulty detachable lead wires may cause erratic behavior of the ECG waveform, SpO<sub>2</sub>, C-Lock, and NIBP due to intermittent ECG lead wire connections.



**Caution** To protect the Propaq Encore from damage during defibrillation, for accurate ECG information, and for protection against noise and other interference, use only ECG electrodes and cables (namely, ones with internal current-limiting resistors) specified or supplied by Welch Allyn, and follow recommended application procedures.

## Preparation



**Warning** Use only ECG safety cables that are designed so that they cannot accidentally be plugged into an AC mains outlet or make contact with other hazardous electrical potentials including earth ground. To prevent damage during defibrillation, don't use ECG cables without 1K series resistors.

Preparing for ECG monitoring with the Propaq Encore requires you to prepare the monitor, prepare the patient, set up the ECG channel, and then set the ECG alarms.

### *Preparing the Monitor*

1. Inspect the ECG cable for wear, breakage, or fraying. Replace the cable if it shows signs of any of these. Plug the ECG cable into the ECG connector on the Propaq's left side panel.



**Warning** Before you use a Propaq on a new patient, always turn it off for a few seconds, then turn it on again. This clears the prior patient's trend values, alarm limit settings, and NIBP cuff inflation target.

2. If the monitor is off, press the OFF/ON switch to turn it on.
3. Select the patient mode appropriate for the patient (Neonatal, Pediatric, Adult). To change patient modes, see page 29.

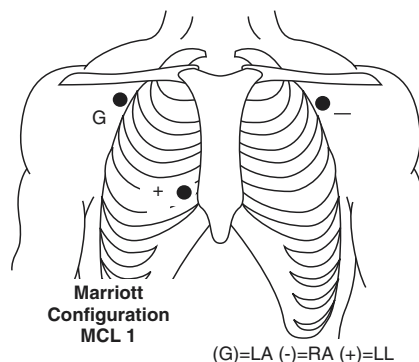
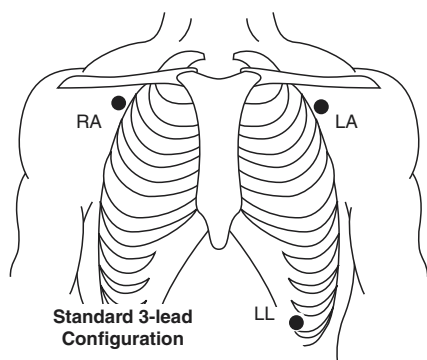
### *Preparing the Patient*

1. Thoroughly clean the skin areas where the electrodes will be attached. Attach lead wires to the electrodes before applying them to the patient.
2. If you are using pre-gelled electrodes, use only electrodes that have not expired. Make sure there is a generous amount of gel in the electrode and that it has not dried. For best results, use silver/silver chloride electrodes.
3. If you are using non-gelled electrodes, apply a  $\frac{1}{4}$  to  $\frac{1}{2}$  inch mound of gel over the electrode contact area.





Some electrodes may be subject to large offset potentials due to polarization. This effect is most likely when dissimilar metals are used for different electrodes, and may be severe enough to prevent obtaining an ECG trace. Furthermore, recovery time after application of defibrillator pulses may be compromised when using electrodes of dissimilar metals. Squeeze bulb electrodes, even if all of the same metal, are particularly vulnerable to this effect. Stainless steel needle electrodes are prone to having large erratic offset drifts, and are not recommended.

4. Apply the electrodes to the patient.



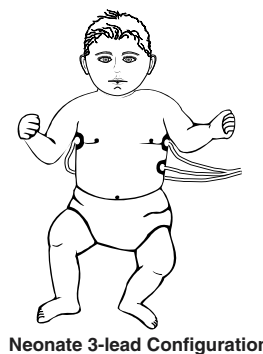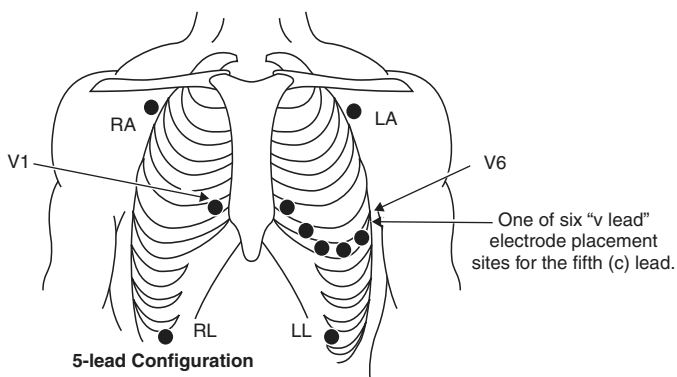If using MCL1, select LEAD II on the Propaq, and use all three electrodes.



5. Support the ECG cable so it does not stress the electrode wires, ECG cable connectors, or electrodes.



Two RESP leads are available. Choose the one that gives you the best signal. If neither signal is adequate, it may be necessary to experiment with nonstandard electrode placement such as placing the RA and LA electrodes on the respective mid-axillary lines just above the level of the nipples.

6. If an electrosurgical unit is going to be used, place the ECG cable and electrode wires as far as possible from the surgical site and from the electrosurgical return electrode and its cables. This will minimize interference.

By now there should be some kind of ECG waveform displayed on the monitor. A heart rate should be displayed to the right of the waveform. Depending on how the Propaq Encore is programmed, a beep tone may occur with each detected QRS event.

7. If there is no waveform, check the electrodes, wires, cable, and the monitor for a possible lead fault.

If an ECG electrode becomes disconnected or disrupted so that the Propaq Encore cannot receive the ECG signal, a message and tone are conveyed with an equipment alert.

### ***Setting Up the ECG/RESP Channel***

Press **ECG** or **ECG/RESP** (available with the Impedance Pneumography Option) to set the selections: **ECG SIZE**, **ECG LEAD**, **RESP SIZE** (available with Impedance Pneumography). The **MORE** button displays the second ECG/RESP menu and a status window with selections for HR/PR TONE, PACER DISPLAY, ECG BANDWIDTH, and RESP LEAD (available with Impedance Pneumography). If the patient being monitored has a pacemaker, you may want to turn on the Pacer indicator function.

### ***Setting ECG/RESP Alarms***

Set the alarm limits according to your hospital's standards.

Motion artifact or other factors can cause false HR/PR alarms. To help minimize false alarms, the Propaq delays or "holds off" triggering an HR/PR alarm for 3 seconds. During this holdoff period, if the Propaq detects that the patient's HR/PR vital sign has returned to acceptable limits, the Propaq cancels the alarm holdoff. The next time an HR/PR limit is violated, the Propaq starts a new 3-second HR/PR alarm holdoff period.

## How ECG/RESP is Displayed

Because of the critical nature of monitoring ECG, it is always displayed in the top part of the waveform display area. The ECG waveform is the only waveform that cannot be turned off using the Setup Wave Select Menu.

Respiration rate numerics are sourced from the CO<sub>2</sub> channel and displayed as BR when CO<sub>2</sub> is active. Otherwise, respiration rate numerics are sourced from the ECG/RESP channel and displayed as RR.

### *Patient Artifact*

Patient movement and other artifact might cause the waveform to move on the display. Most artifact such as this is automatically detected, and the waveform is adjusted so that it always remains centered in the waveform window.

Severe artifact and interference (such as interference from defibrillation) may cause the waveform to move off the display. The Propaq Encore will always automatically reposition the waveform in just a few seconds so you can see it again.

### *ECG/RESP Menus and Status Window*



<b>ECG SIZE</b>	Selects the ECG waveform size; sizes are shown in millivolts per centimeter (.2, .5, 1, 2, or 4 mV/cm) to the left of the waveform.
<b>ECG LEAD</b>	Selects the ECG lead. The available leads are lead I, II, III, aVR, aVL, aVF, or V. The Propaq's factory default lead setting is Lead II.

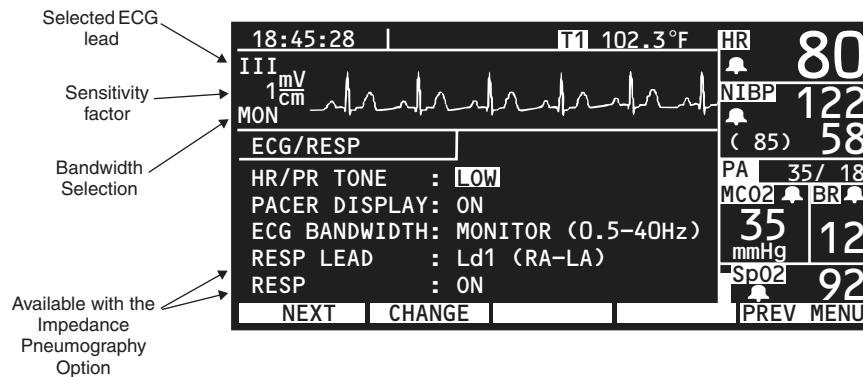
### *Impedance Pneumography Selections*

If your Propaq Encore includes the Impedance Pneumography Option, the following selection is also available:

<b>RESP SIZE</b>	Selects the RESP waveform size.
------------------	---------------------------------

The Size function “increases” and “decreases” the ECG or RESP waveform size. Each time you press a **SIZE** button, the waveform approximately doubles in height. When you reach the largest waveform size, the next press displays the smallest size.

When you press **MORE** in the first ECG/RESP Menu, a status window appears showing you the current ECG/RESP settings and additional selections.



- HR/PR TONE** Sets the heart tone loudness to LOW, MEDIUM, HIGH, or OFF. If SpO2 is monitored, the pitch of the tone varies with the SpO2 value.
- PACER DISPLAY** Turns on and off the pacer indicator in the ECG waveform.
- ECG BANDWIDTH** This selection allows you to determine the bandwidth for the data sent to the display and the printer. If the selection is Extended, the bandwidth is 0.05-40 Hz. If the selection is Monitor, the bandwidth is 0.5-40 Hz.

The Size function “increases” or “decreases” the ECG or RESP waveform size. Each time you press a **SIZE** button, the waveform approximately doubles in height. When you reach the largest waveform size, the next press displays the smallest size.



The QRS detector sensitivity threshold is not affected by changing the ECG display size. Likewise, the RESP breath detector threshold is not affected by changing the RESP display size.

## ***Impedance Pneumography Selections***

If your Propaq Encore includes the Impedance Pneumography Option, the following selections are also available.

- RESP LEAD** Selects the RESP lead. Choices are RA-LA and RA-LL, and choice is independent of ECG lead selection. Experiment with placement for best signals.
- RESP** Turns impedance pneumography (RESP channel) on or off.



The Propaq Encore counts as “breaths” respiratory efforts that are larger than two times background Cardiovascular Artifact (CVA).

Since RESP is derived from the same leads as the ECG channel, the Propaq Encore determines which signals are cardiovascular artifact and which signals are a result of respiratory effort. If the breath rate is within five percent of the heart rate or a multiple or sub-multiple of the heart rate, the monitor may ignore breaths and trigger an apnea alarm.

## Using the Propaq Encore With Pacemaker Patients



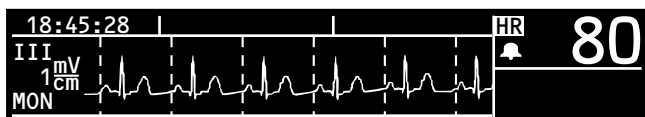
**Warning** Pacemaker signals can differ from one pacemaker to the next. The Association for Advancement of Medical Instrumentation (AAMI) cautions that “in some devices, rate meters may continue to count the pacemaker rate during occurrences of cardiac arrest or some arrhythmias. Do not rely entirely upon rate meter alarms. All pacemaker patients should be kept under close or constant observation.”

### *Pacemaker Signals*

If the patient being monitored has a pacemaker, the Propaq Encore detects and can indicate the occurrence of pacemaker signals. With the Propaq, pacemaker signals are not counted as heart beats as long as the pacemaker signal meets the pulse amplitude, pulse width, and overshoot/undershoot specifications listed in Appendix B.

### *Pacer Display*

On the Propaq Encore display, vertical dashed lines indicate each time a pacemaker signal is detected when the Propaq Encore PACER function is turned on. The waveform “spike” produced by the pacemaker will also be displayed if it contains sufficient energy. Whether the pacer is atrial, ventricular, or both, the indicator and the spike appear. If the PACER function is turned off, only the pacemaker spike is displayed:



### *Turn On and Off the Pacemaker Indicator*

The status of the pacer display is shown in the ECG status window. Turn PACER DISPLAY on or off using the **NEXT** and **CHANGE** buttons.

## Noise on the signal

Noise on the ECG signal may be detected as pacer signals, causing the pacer indicator to appear on the display. If you don't need to indicate pacemaker signals, you may want to turn off the pacemaker indicator for a better display of the ECG waveform.



**Warning** The presence of much pacer-like noise can cause the displayed heart rate to be erratic even though the ECG trace may look clean with the pacer indicator off. Fix the noise problem by using fresh ECG electrodes and an ECG cable whose lead wires make good connections.

## Using the Filter to Better Display a Waveform

The Propaq Encore includes a filter that reduces noise from the ac power signal and produces a much clearer ECG waveform. To set the filter, press **SETUP, MORE, MORE, SERVICE, YES** (to access the Service Menu), **MORE, MORE, SETTINGS** to display the Settings Menu. Select the **FILTER** function with the **NEXT** button and press the **CHANGE** button to change settings.

Set the filter to the ac mains frequency of your hospital (either 60 or 50 Hz).

## ECG Messages

If a lead fault occurs, an ECG equipment alert typically shows which lead failed. Otherwise, the Propaq Encore displays **MULTIPLE** in place of the failed leads.

**ECG LEAD CHANGED.** The Propaq Encore has automatically changed an ECG lead due to a lead wire or electrode problem.

**LEAD FAIL: REPLACE ELECTRODES.** There is excessive offset on the ECG channel.



## RESP Messages

**LEAD FAIL.** One or more electrodes are making very poor or no contact. Check for proper connection; replace electrodes if needed.

**INAPPROPRIATE ECG CABLE.** ECG cable appears not to contain 1K current limiting resistors. These resistors are required for RESP operation and to protect the Propaq Encore from damage during defibrillation. Replace cable with proper type.

**NOISY SIGNAL, CHECK ELECTRODES.** Electrodes are making poor contact and may be dried out. Replace electrodes.



**Warning** If a disconnected lead is in too close proximity to other electrical devices, it may cause false heart rate, a failure to detect apnea, or a failure to display a Lead Fail message.

## Invasive Pressure

*This section applies only to Propaq Encore Models 204EL and 206EL. If you don't have one of these models, you can skip this section.*

### Intended Use

The Propaq Encore invasive pressure channel is intended for measuring arterial, venous, and intracranial pressures using invasive transducers with 5  $\mu\text{V}/\text{V}/\text{mmHg}$  sensitivity. The Propaq Encore can be used with many types of transducers, including nondisposable, disposable dome, and fully disposable.

### Invasive Pressure Connectors and Transducers

Recommended transducers for use with the Propaq Encore are listed in the Welch Allyn *Products and Accessories* booklet. Do not use light-sensitive disposable transducers. Transducers must be used according to your hospital's standards and the manufacturer's recommendations. Always refer to the manufacturer's Directions for Use before using the transducer.



**Warning** If electrocautery is used, always avoid using any transducer with a conductive (metal) case that is electrically connected to its cable shield. Using a conductive transducer case with such a shield connection risks high-frequency burns at the ECG electrodes if the transducer case becomes earth grounded.

Although complete disconnections of invasive pressure transducers will be detected by the normal alarm functions, partial disconnection will not be detected, nor will the use of some incompatible transducers. The user must exercise reasonable measures to ensure that approved transducers are used and that pressure transducers are connected properly.

Before you use a Propaq on a new patient, always turn it off for a few seconds, then turn it on again. This clears the prior patient's trend values, alarm limit settings, and NIBP cuff inflation target.

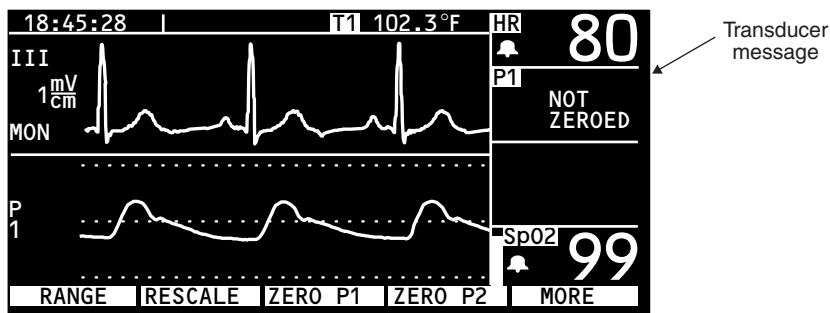
## Preparation

Preparing for invasive pressure monitoring with the Propaq Encore requires you to prepare the transducer, zero the transducer, set up the pressure channel, and set the invasive pressure alarm limits.

### Preparing the Transducer

1. Inspect the transducer cable for wear, breakage, or fraying. Replace it if the cable shows signs of any of these. Replace the transducer dome if necessary.
2. Apply the transducer according to your hospital's procedures. Always refer to the transducer manufacturer's Directions for Use.
3. If the transducer is a disposable unit with separate cable, connect the transducer to the transducer cable. Plug the transducer cable into an invasive pressure connector on the left side panel.

The message NOT ZEROED (or NO ZERO, depending on the zone) immediately appears in the blood pressure numerics window for the invasive pressure channel being used.



### Zeroing the Transducer

1. To zero the transducer, open the transducer's stopcock to atmospheric air. Allow a few seconds for the transducer to settle.
2. If the ZERO menu is not displayed, press the following Propaq Encore buttons: **INV PRS**, then **ZERO P1** (or **ZERO P2**). The word ZEROING appears in the numerics window during zeroing. The button label changes to **CANCEL** to allow you to cancel the zeroing process if necessary.

3. Wait for a tone to briefly sound and the word ZEROED to appear in the blood pressure numerics window. You will then see the pressure scale to the left of the waveform, and the pressure numerics appear.
4. Close the transducer's stopcock.
5. If the transducer will not zero, the words ZERO REJECTED (or NO ZERO, depending on zone) will appear in the numerics window. Press **CANCEL** and try zeroing again. You won't see the pressure values and the scales until an acceptable zero reference is established.
6. Check that the transducer is open to atmospheric air and that it is properly connected to the Propaq Encore, then try zeroing again. The Propaq Encore will not allow zeroing to occur if the pressure waveform is pulsatile, if there is too much noise in the signal or if the transducer's offset is too great. Once the channel is zeroed, the pressure scale appears next to the waveform.

If the transducer still does not zero, try another transducer or another cable.

### ***Setting up the Pressure Channel***

Press **INV PRS** to set the invasive pressure channel selections: **RANGE**, **RESCALE**, and **ZERO P1/ZERO P2**. Press **MORE** to set **LABEL P1/LABEL P2** and **FORMAT 1/FORMAT 2**.

### ***Setting the Invasive Pressure Alarms***

Set the alarm limits according to your hospital's standards.

### ***Rezeroing a Transducer***

You can rezero a transducer at any time, after again opening the transducer stopcock to atmospheric air. If the transducer has already produced pressure readings, rezeroing provides a new zero reference for the Propaq.

If the zero value is not accepted, the Propaq Encore continues to use the previous zero reference and displays the pressure values and waveforms based on that value. If the new zero value is accepted, the new pressure values based on the new zero value are displayed, and the waveform is adjusted according to the new scale.



**Warning** If a **ZERO** button is pressed after an invasive pressure channel has been successfully zeroed and is currently monitoring a pressure waveform, the message **ZERO REJECTED** will display in the invasive pressure numerics window. This message will preempt the valid invasive pressure numerics until the **CANCEL** button in the Invasive Pressure Menu is pressed.

If the invasive pressure channel enters an alarm condition while the **ZERO REJECTED** message is overriding the invasive pressure numerics, no invasive pressure numerics will flash to indicate invasive pressure is in alarm.

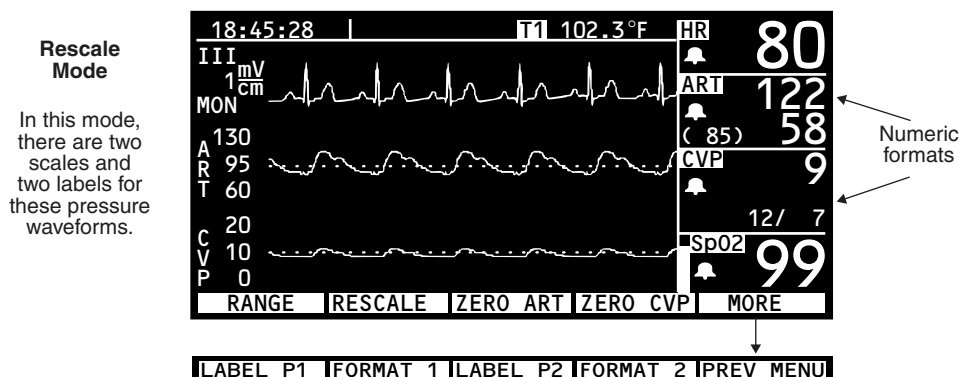
To remove the **ZERO REJECTED** message and to restore the invasive pressure numerics during an invasive pressure alarm, you must return to the invasive pressure menu and press **CANCEL**. This will restore the invasive pressure numerics.

## How Invasive Pressure is Displayed

From the invasive pressure signal, the Propaq Encore displays both a pressure waveform and pressure numeric values (systolic, diastolic, and mean). The waveform is displayed in a waveform window (if the waveform is turned on in the wave select window). The numerics are displayed in the blood pressure numerics windows.

The Propaq Encore allows you to identify the pressure measurement with a selectable label, and the numerics can be displayed in different formats.


The pressure waveform scales are not displayed until you zero the transducer. Once the zero reference has been established the scales automatically appear.



<b>RANGE</b>	Sets the display to Range Mode. All invasive pressure waveforms monitored are displayed against the same scale. You can select one of five Propaq Encore pressure scales. If two waveforms have a great difference in their pressures, the higher pressure waveform may not be visible if it is out of range of the scale. Press <b>RANGE</b> until the desired scale appears.
<b>RESCALE</b>	Sets the display to Rescale Mode. Each invasive pressure waveform is displayed against its own scale. Each time you press the button, the scale is automatically selected based on the highest and lowest pressure levels of each pressure waveform.
<b>ZERO</b>	Zeroes the selected pressure channel, or cancels zeroing in process.
<b>CANCEL</b>	The <b>ZERO</b> button changes to <b>CANCEL</b> while zeroing.
<b>LABEL</b>	Selects a label for the pressure channel. The selectable labels are: ART—arterial, PA—pulmonary artery, CVP—central venous pressure, ICP—intracranial pressure, UA—umbilical artery, and UV—umbilical vein.  You can still use the generic Propaq Encore pressure label, P1 or P2.
<b>FORMAT</b>	The Propaq Encore displays the invasive pressure values in two different numeric formats in the pressure numerics window. You can select which pressure value(s) are most prominently displayed.

**Range Mode**

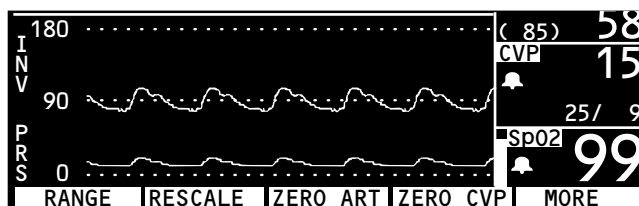In this mode, there is one scale and one label for both pressure waveforms.



## Invasive Pressure Messages

The following messages can appear in the numerics window.

**NOT ZEROED (or NO ZERO).** No zero reference has been established. The monitor displays the pressure waveform, but to protect against erroneous readings, the pressure waveform scale is not displayed. To remove this message, zero the transducer.

**ZEROING.** This message briefly appears as the transducer is being zeroed.

**ZEROED.** This message appears after the zero value has been accepted. It remains for eight seconds and is replaced by the current pressure values.

**ZERO REJECTED (or REJECT).** Unable to establish a zero reference value. The message remains until the **CANCEL** button is pressed.

**CANCELED.** This message appears if **CANCEL** is pressed while the channel is zeroing.

These equipment messages can appear in an equipment alert window.

**TRANSDUCER NOT DETECTED.** The transducer connection is broken.

**TRANSDUCER SHORT CIRCUIT.** This message appears when the Propaq Encore senses a short in the transducer. The transducer should be replaced.

**INCOMPATIBLE TRANSDUCER.** Check the compatible transducers listed in the Welch Allyn *Products and Accessories* booklet to ensure you are using a compatible transducer.

## NIBP

### Intended Use

The Propaq Encore noninvasive blood pressure channel (NIBP) indirectly measures arterial pressures using an inflatable cuff. If ECG is also monitored, the Propaq Encore synchronizes the NIBP measurement process to the occurrences of the R-wave, increasing accuracy in cases of extreme artifact, diminished pulses, or some dysrhythmias.



**Warning** The patient's limb should be periodically observed to ensure that the circulation is not impaired for a prolonged period of time.

The Propaq Encore should never be used to monitor NIBP on one patient while simultaneously monitoring ECG on another patient.

If a noninvasive blood pressure measurement is suspect, repeat the measurement. If you are still uncertain about the reading, use another method.

Do not attempt to take NIBP pressures on patients during cardiopulmonary bypass.

Some or all NIBP safety functions are disabled in the NIBP TEST screen in the Service Menu. Do not attempt to conduct NIBP TEST when the cuff is attached to a patient.

The Propaq Encore NIBP channel has been calibrated to agree with a central invasive blood pressure. Diastolic pressures may be 5 to 10 mmHg lower than the auscultatory equivalent. Systolic pressures may be lower than radial invasive equivalent.

**Neonatal Mode** is intended for use on infants of up to about 44 weeks gestational age in neonatal care settings. The Neonatal Mode provides the lowest cuff pressure and shortest inflation time limits to ensure patient safety and comfort.

**Pediatric Mode** is intended to be used on larger infants and small children up to nine years old in pediatric care settings. This mode supports the widest range of cuff sizes and a higher range of patient numerics for the hypertensive infant or child while still restricting the cuff pressure and inflation times to limits lower than those allowed for adults.

**Adult Mode** provides the full range of patient numerics and cuff pressures but limits the cuff sizes available to the standard child cuff and larger.



## NIBP

NIBP measurements are affected by normal physiological pressure variations from reading to reading. Normal respiration may affect pressure by as much as 10 to 20 mmHg. Patient's emotional state, body position, and cuff fit may also adversely affect NIBP measurements. In some individuals, the act of taking blood pressure readings may alter the blood pressure. Successive readings on the same patient may vary for the above reasons.

The static accuracy of the Propaq's internal manometer can be verified by a qualified biomedical engineer using a mercury column manometer (refer to the *Propaq Encore Service Manual*). The accuracy of the Propaq's determination of systolic, diastolic, and mean pressures in a clinical setting can only be assessed by careful statistical analysis of controlled clinical trials of representative patient populations.

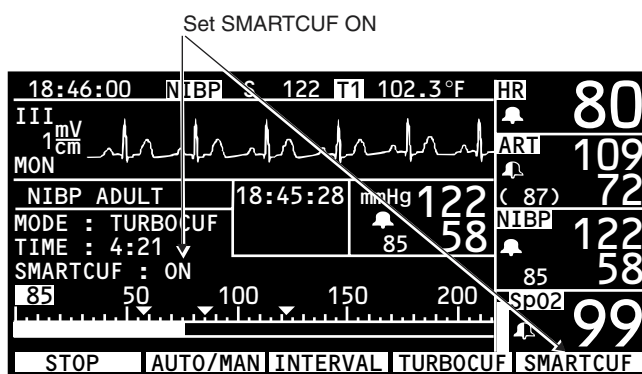
## Improve NIBP Accuracy with Smartcuf™

NIBP measurements can be adversely affected by many factor such as cardiac arrhythmias, sudden changes in blood pressure, body motions such as convulsions or shivering, bumping the cuff, vibration, vehicle motion, or weak pulses.


The patented Smartcuf software filtering technology greatly increases NIBP measurement accuracy in the presence of motion artifact or diminished pulses. Smartcuf synchronizes the NIBP reading with the R-wave of the patient's ECG to eliminate noise created by external stimuli such as patient motion or vibration. The monitor must perform ECG monitoring while using Smartcuf.

To enable the Smartcuf filter:

- Connect the ECG leads to the patient and perform ECG monitoring during NIBP.
- From the Main Menu, press **NIBP** to display the NIBP Menu and set Smartcuf to ON.



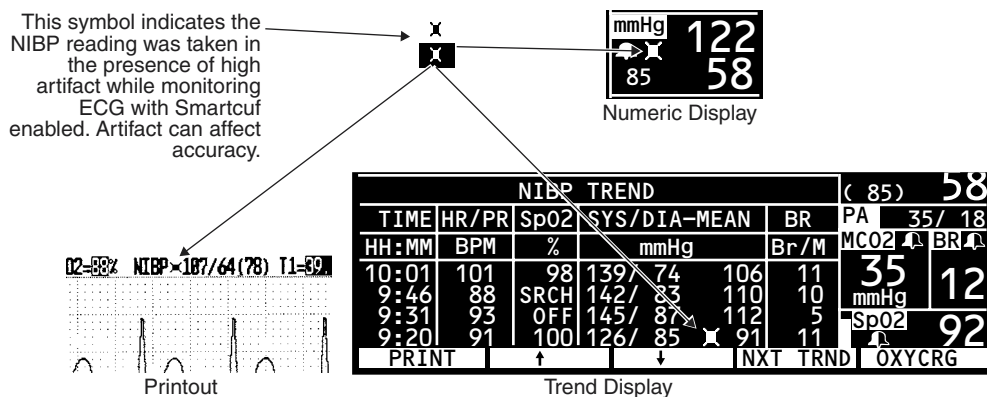If artifact is so severe while Smartcuf is enabled that it affects the accuracy of an NIBP measurement, that measurement is marked with a special symbol on the display and on printouts



There may be some situations where it is desirable to disable Smartcuf. This may include situations with very extreme motion artifact, certain types of arrhythmias, or other situations where it is not possible to obtain a good ECG signal. NIBP measurements can still be performed when Smartcuf is disabled.

To disable Smartcuf, from the Main Menu press **NIBP** to display the NIBP Menu and set Smartcuf to OFF.

## NIBP Connector and Cuffs

Use Propaqs with approved hoses and cuffs listed in the Welch Allyn *Products and Accessories* booklet. The Propaq Encore uses a single-hose cuff. Cuffs that conform to AAMI or AHA guidelines should be used. Select the proper size of cuff based on the limb circumference.

	Neonate	Pediatric	Adult
Hoses	Neonate/Infant	Adult	Adult
Cuffs (typical cuff labeling)	Neonate #1 to #5 (disposable); newborn, infant (reusable)	Newborn, infant, small child, child, small adult	child, small adult, adult, large adult, thigh
Recommended limb circumference	up to 15 cm	7.7 to 25 cm	greater than 15 cm



**Warning** When monitoring NIBP, match the Propaq patient mode to the style of the cuff. For neonates, set the monitor to Neonatal Mode unless the circumference of the limb is too large for the cuff. In that case, use the Pediatric Mode. In the Pediatric Mode, the maximum cuff inflation pressure can exceed 150 mmHg, and two retries are allowed.

For information about patient mode specifications, see “NIBP Specifications” on page 142.

## Preparation

Setting up for noninvasive blood pressure monitoring requires three steps: place the cuff on the patient and connect the cuff to the monitor, set up the NIBP channel, and set the NIBP alarm limits.



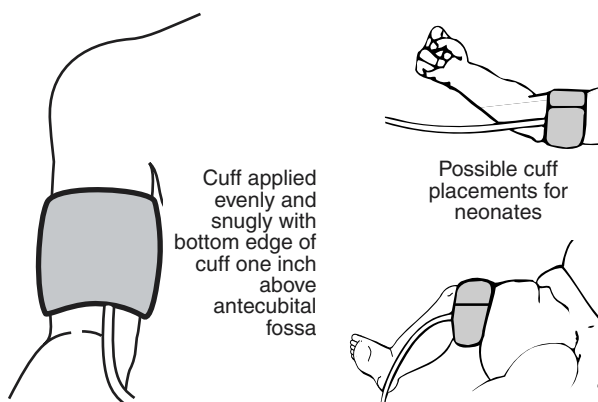
**Warning** Before you use a Propaq on a new patient, always turn it off for a few seconds, then turn it on again. This clears the prior patient’s NIBP cuff inflation target, trend values, and alarm limit settings.

At powerup, the Propaq has an NIBP default inflation pressure (cuff inflation target) based on the patient mode (see page 142 for the values). After each NIBP measurement, the Propaq adjusts the target inflation pressure to optimize the next NIBP measurement. To avoid possible patient discomfort, be sure to turn the monitor off and then on between different patients to reset the cuff inflation target to the default value.

### Place and Connect the Cuff

1. Squeeze as much air from the cuff as you can before placing it on the patient.
2. Place the cuff on the limb.

When you place the cuff, it should ideally be placed at the same level as the heart. If above the heart, add 1.9 mmHg to the NIBP measurement for every inch the cuff is above the heart. If below the heart, subtract 1.9 mmHg for every inch.



The cuff should fit snugly, but not be uncomfortable. The hose must not be kinked or pinched.

Ensure that the cuff tubing is centered over the brachial artery.



A cuff that is not properly connected to the patient may result in a false reading if the patient and cuff are moved by motion artifact or clinical personnel during the NIBP measurement. Always verify the cuff is properly placed on the patient.


3. Screw the hose connector onto the NIBP connector on the monitor's left side.
4. If motion artifact such as shivering, coughing, or vehicle motion interferes with NIBP readings, do the following:

Position the patient's limb away from the body so the applied cuff is not in contact with the patient's body or any other object such as a bed rail. Try to keep the cuff at the same level as the heart.

Connect the ECG leads to the patient and perform ECG monitoring during NIBP.

### Set Up the NIBP Channel

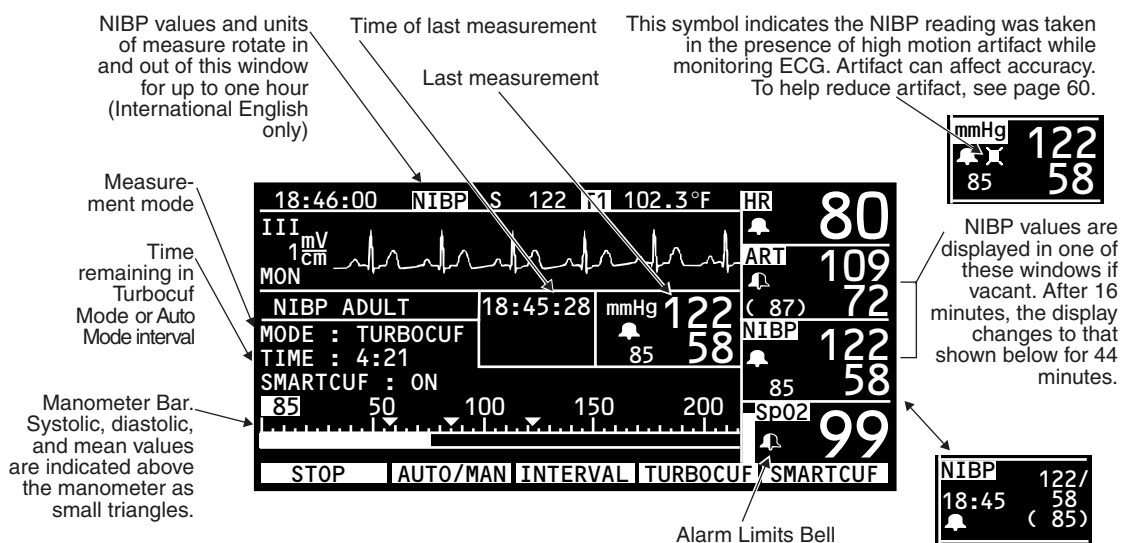Press the **NIBP** button to display the status window and menu.



**START/STOP** Starts and stops NIBP measurements. Any time the Propaq Encore is taking a noninvasive pressure measurement, the **START** button changes to **STOP** so you can stop the measurement in progress. This button initiates the same action as the **START/STOP** button at the left side of the screen. Pressing **STOP** will automatically vent the cuff.

## NIBP

**AUTO/MAN** This button switches the mode between Automatic or Manual Mode. The Manual Mode is the default unless you change it by reprogramming your Propaq. Measurements can be taken at intervals of 1, 2, 3, 5, 10, 15, 30, and 60 minutes. Press **START** to initiate a measurement.

**INTERVAL** Selects the interval at which NIBP measurements are automatically taken. The interval you select, ranging from one minute to 60 minutes, is shown on the display next to the word TIME.

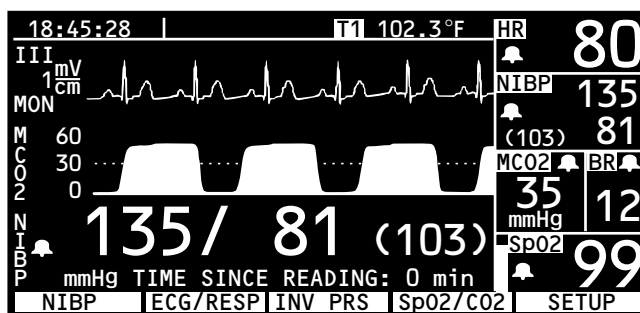
**TURBOCUF** Automatically starts NIBP measurements and continues to take as many measurements as possible within five minutes.

**SMARTCUF** Enables or disables the Smartcuf motion artifact filter. NIBP measurements can still be taken when Smartcuf is off. Artifact may interfere with the accuracy of NIBP measurements with Smartcuf off.


### NIBP Displayed in Waveform Window

By turning on NIBP in the wave select window, the NIBP numerics can be displayed in a waveform window.


The numerics are shown in large characters for 16 minutes for each new measurement taken...



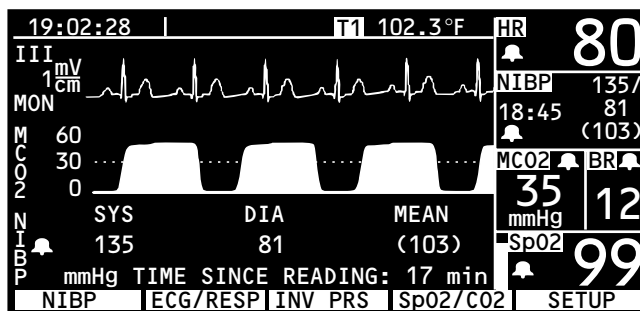... and then they change to the smaller characters for 44 minutes. The numerics are removed after 60 minutes.



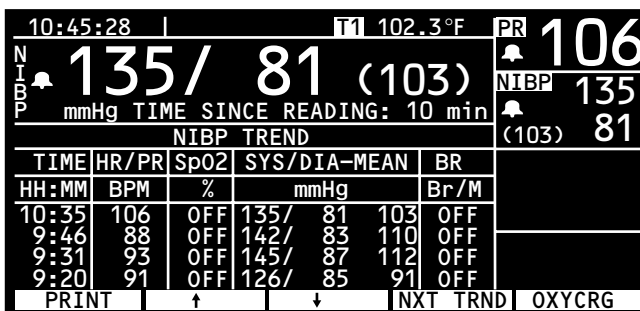If NIBP is the only vital sign being monitored, the numerics are displayed in a waveform window above a trend window.



## Set the NIBP Alarm Limits

Set the alarm limits according to your hospital's standards.

## NIBP Display Default Settings

To select which vital sign waveforms are displayed, press **MAIN MENU**, **SETUP**, **WAVE SEL** to display the Wave Select window.

SETUP		WAVE SELECT		( 85) 58	
ECG	: ON	RESP	: OFF	PA	35/ 18
ART	: <b>ON</b>	SpO2	: ON	MC02	BR
PA	: OFF	NIBP	: ON	35	12
CO2	: OFF			mmHg	
				SpO2	92
NEXT	ON/OFF	INSERT		PREV	MENU

To enable the display of a waveform or large NIBP numerics, select ON (ECG cannot be set to OFF). The monitor displays the first three active waveforms set to ON in the order of priority listed in the Wave Select window. If NIBP is ON and only one or two other waveforms are ON and active, the monitor displays large NIBP numerics in a waveform window.

## Important Information About Automatic Measurements

A blood pressure measurement will begin when the minute of the time of day clock is evenly divisible by the interval. For example, if the interval is set to 10 (minutes), measurements will begin at the hour and at 10, 20, 30, 40, and 50 minutes past the hour. Note, however, that for intervals 1, 2, or 3 (minutes), measurements begin 1, 2, or 3 minutes after the interval is set. For example, if the 1 minute interval is selected at 10:45:20, the next measurement starts at 10:46:20.

The start time may be delayed if the previous measurement ended within 30 seconds of the scheduled start time, because the monitor requires that the cuff pressure be below 5 mmHg for a minimum of 30 seconds between measurements to allow time to restore blood flow to the limb.

## NIBP Messages

The following NIBP messages can appear in the equipment alert window. An NIBP caution message also appears in the numerics window. If an error number (ERR# x) is listed in an NIBP trend printout or display, it indicates that the corresponding NIBP equipment alert occurred.

**AIR LEAK, CHECK HOSE (ERR# 1).** The Propaq Encore could not properly inflate cuff. Check the hose and cuff for obvious leaks, such as O-rings in the hose connections.

**CUFF NOT DETECTED (ERR# 2).** During cuff inflation the detected pressure did not sufficiently rise. Check that the cuff connection is tight and take the measurement again.

**KINKED HOSE, CHECK HOSE (ERR# 3).** The Propaq Encore could not properly inflate cuff. Check for a kinked hose between the monitor and the patient.

**OVERPRESSURE CONDITION (ERR# 4).** The pressure in the cuff exceeded the acceptable limits for patient mode. Check the hose and try taking another measurement.

**WEAK PULSES, CAN'T FIND SYS/DIA (ERR# 5).** There are not enough pulses to determine the systolic or diastolic pressures, but a mean pressure is available. Try reapplying the cuff after squeezing as much air from it as you can.

**ARTIFACT, CAN'T FIND SYS/DIA (ERR# 6).** The systolic or diastolic pressures are unreliable due to artifact, but a mean pressure is available. May be caused by patient motion.

**NO PULSES DETECTED (ERR# 7).** The cuff may not be properly applied to the patient, or the patient may not have detectable pulses due to shock or arrhythmias.



**Warning** The Propaq Encore cannot differentiate between physiologic and cuff application causes of the NO PULSES DETECTED message. Always evaluate the patient for presence of life threatening conditions whenever this message occurs.

**CONNECT ECG TO REDUCE NIBP ARTIFACT (ERR# 8).** NIBP artifact prevents a valid reading. Connect ECG electrodes to improve NIBP measurements.

**NO VALID BLOOD PRESSURE FOUND (ERR# 9).** This message can occur due to motion artifact, the Propaq Encore being set in the wrong patient mode, or the wrong hose or cuff being used in relation to the patient mode.

**CALIBRATING, PLEASE WAIT (ERR# 10).** The Propaq Encore periodically recalibrates the NIBP channel to ensure it can properly make NIBP determinations. Normal monitor operation continues while the NIBP channel is calibrating. If the NIBP channel has not updated its calibration in 15 minutes, the channel will briefly deactivate until a new calibration has occurred.

**LOW BATTERY, NIBP DISABLED (ERR# 11).** The battery lacks sufficient voltage to be able to operate the NIBP channel. Connect the Encore to the ac power adapter.

**SERVICE REQUIRED, NIBP DISABLED (ERR# 12).** Have the monitor serviced.

**CUFF TOO LARGE FOR PATIENT MODE (ERR# 13).** The monitor detects a cuff too large for the current patient mode. First, verify the patient mode. If the patient mode is correct, make sure the cuff fits snugly. If this alert occurs in Neonatal Mode, change the patient mode to Pediatric Mode and check alarm limits. If the alert occurs in Pediatric Mode, change to Adult Mode and check the alarm limits. Note that different pressures and retries are used for each mode as stated in “NIBP Specifications” on page 142.

**KINKED OR NEONATE HOSE (ERR# 14).** This message occurs when the neonate hose is detected in adult patient mode. Change the hose or the patient mode selection.

**ARTIFACT PRESENT, MINIMIZE ARTIFACT (ERR# 15).** The monitor has detected too much artifact to allow accurate readings. Take steps to reduce artifact. Position the patient’s limb away from the body so the applied cuff is not in contact with the patient’s body or any other object such as a bed rail. If the Smartcuf motion artifact filter is on, make sure that the ECG leads are properly connected to perform ECG monitoring during NIBP. If the Smartcuf motion artifact filter is off, consider turning it on (and connect ECG if not already connected).

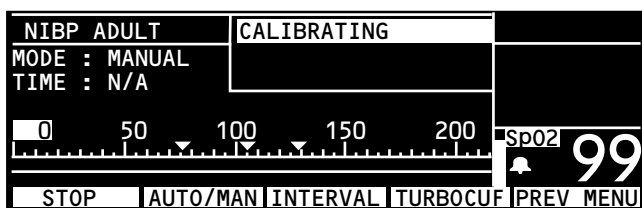
The following messages can appear in the NIBP status window.

**CALIBRATING.** The NIBP channel is running an internal calibration.

**DISABLED, LOW BATT.** See LOW BATTERY, NIBP DISABLED above.

**NIBP DISABLED, SERVICE REQUIRED.** See SERVICE REQUIRED, NIBP DISABLED above.

**RETRY.** Since the Propaq Encore did not receive a valid NIBP reading, it will automatically attempt to take another reading.





***NIBP IN PROGRESS Message***

The monitor displays this message when noise or artifact such as vehicle motion causes a delay while measuring NIBP. To remove the message, press any button below the screen. To cancel the NIBP measurement, press **STOP**.

# Temperature

## Intended Use

Propaq Encore monitors provide two temperature channels (except for the HP-side panel option). When both channels are active, the difference temperature ( $\Delta T$ ) is also displayed. You can select °C or °F.

## Temperature Connectors and Probes

Propaqs are intended to be used with approved temperature probes listed in the Welch Allyn *Products and Accessories* booklet. Other temperature probes that do not match the performance specifications of these approved probes may produce incorrect temperature readings.

## Preparation

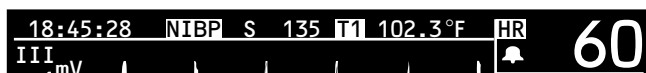


**Warning** Application and use of metal-jacketed temperature probes that come in contact with conductive objects or clinical personnel during electrocautery may cause burns at the patient-probe/electrode contact points.

1. Place the probe on the patient, and plug it into one of the connectors on the Propaq's side panel. Within a few seconds, the Propaq Encore will display the temperature.
2. To select the temperature units (°C or °F), press **SETUP, MORE, MORE, SERVICE, YES** (to access the Service Menu), **MORE, MORE, SETTINGS**. Use the **NEXT** and **CHANGE** buttons to select and set the temperature units as desired. The Propaq Encore automatically updates the temperature display to show the newly selected units. Changing units does not clear Temperature trends.
3. Set the alarm limits according to your hospital's standards.

## How Temperature is Displayed

Temperature is displayed as a numeric only, in a window at the top of the Propaq Encore screen, in °C or °F. This area displays all temperature measurements (T1, T2,  $\Delta T$ ), one at a time.



## Temperature Messages

The following messages can appear in an equipment alert window. A temperature caution message will also appear in the temperature numeric window when one of these messages appears (except **PROBE NOT DETECTED**).

**PROBE NOT DETECTED.** This message occurs when the Propaq Encore has successfully measured temperature and a probe is then disconnected. Reconnect the probe or acknowledge the equipment alert by pressing any menu button.

**PROBE SHORT.** Verify that the probe is properly inserted in the left side panel. If so, replace probe.

**CALIBRATION ERROR, TEMP DISABLED.** This message appears when the Propaq Encore has detected that it cannot accurately measure the temperature. The monitor should be serviced.

Malfunction of the temperature probes may result in inaccurate readings. Confirm suspect readings.

## Pulse Oximetry (SpO<sub>2</sub>)



**Warning** Oxygen saturation measurements using pulse oximetry are highly dependent on proper placement of the sensor and patient conditions. Patient conditions such as shivering and smoke inhalation may result in erroneous oxygen saturation readings. If pulse oximetry measurements are suspect, verify the reading using another clinically accepted measurement method, such as arterial blood gas measurements on a co-oximeter.

Tissue damage can be caused by incorrect application or use of a sensor (e.g., wrapping the sensor too tightly, applying supplemental tape, failing to periodically inspect the sensor site, leaving a sensor on too long in one place). Refer to the Directions for Use provided with each sensor for specific instructions on application and use, and for description, warnings, cautions, and specifications.

Sensors exposed to ambient light while not applied to a patient can exhibit semi-normal saturation readings. Be sure the sensor is securely placed on the patient and check its application often to ensure accurate readings.

Inaccurate measurements may be caused by venous pulsations.

The pulse oximetry option can be used during defibrillation, but the readings may be inaccurate for a short time.

**Interfering Substances:** Carboxyhemoglobin may erroneously increase readings; the level of increase is approximately equal to the amount of carboxyhemoglobin present. Methemoglobin may also cause erroneous readings. Dyes, or any substances containing dyes, that change usual arterial pigmentation may cause erroneous readings.

Before you use a Propaq monitor on a new patient, always turn off the monitor for a few seconds, then turn it on again. This clears the prior patient's trend values, alarm limit settings, and NIBP cuff inflation target.

Each SpO<sub>2</sub> sensor is designed for application to a specific site on the patient within a certain size range. To obtain optimal performance, use an appropriate sensor and apply it as described in the sensor's directions for use.

If excessive ambient light is present, cover the sensor site with opaque material to block the light. Failure to do so may result in inaccurate measurements. Light sources that can affect performance include surgical lights (especially those with a xenon light source), bilirubin lamps, fluorescent lights, infrared heating lamps, and direct sunlight.

## Pulse Oximetry (SpO<sub>2</sub>)

If NIBP will be monitored while using SpO<sub>2</sub>, place the NIBP cuff on a different limb than the SpO<sub>2</sub> sensor to help reduce unnecessary SpO<sub>2</sub> alarms. For optimal measurements, avoid placing the SpO<sub>2</sub> sensor on the same limb as an arterial catheter or intravascular line.

Loss of pulse signal can occur if the sensor is too tight, there is excessive ambient light, an NIBP cuff is inflated on the same limb as the sensor, there is arterial occlusion proximal to the sensor, the patient is in cardiac arrest or shock, or the patient has hypotension, severe vasoconstriction, severe anemia, or hypothermia.

## Perform SpO<sub>2</sub> Monitoring with Masimo Option

1. Attach the sensor to the patient according to the sensor manufacturer's instructions, observing all warnings and cautions.



**Warning** Use only Masimo accessories and sensors with the monitor with Masimo SpO<sub>2</sub> option as listed in the Welch Allyn *Products and Accessories* booklet (810-0409-XX).

2. Inspect the Masimo SpO<sub>2</sub> cable. Replace it if it shows any signs of wear, breakage, or fraying. Plug the sensor into the cable and plug the cable into the Propaq monitor.

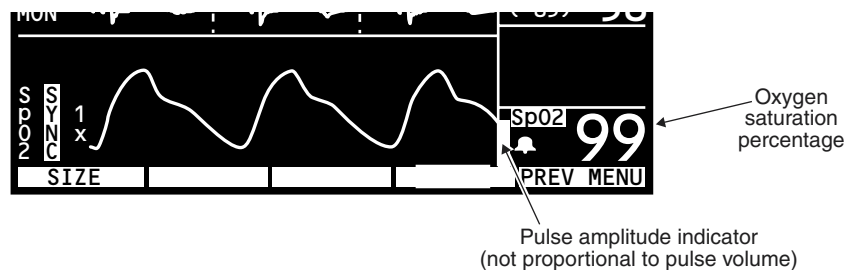


The monitor displays STANDBY in the SpO<sub>2</sub> numeric window until it measures and displays the SpO<sub>2</sub> value.

As oxygen saturation increases and decreases, the pitch of the heart tone rises and falls.

The monitor self-calibrates the SpO<sub>2</sub> channel whenever the monitor is first turned on or a sensor is first connected to the SpO<sub>2</sub> channel.

3. From the Main Menu, press **SpO<sub>2</sub>** (or **SpO<sub>2</sub>/CO<sub>2</sub>**, then **SpO<sub>2</sub>**) to display the SpO<sub>2</sub> menu similar to the following:



4. Press **SIZE** to adjust the waveform size for best viewing (1x, 2x, 4x, or 8x).
5. Adjust the placement of the sensor until a good SpO<sub>2</sub> waveform is displayed. A waveform with artifact may cause erroneous oxygen saturation readings.
6. Set alarm limits according to your hospital's standards.



To help minimize false alarms, the Propaq monitor briefly delays or "holds off" triggering both audible and visual alarms for limit violations for SpO<sub>2</sub>% and Pulse Rate for 10 seconds. After the alarm hold-off period begins, if the monitor detects that the patient's vital sign has returned to acceptable limits, the monitor cancels the alarm hold-off. The next time a vital sign limit is violated, the monitor starts a new hold-off period.

The "averaging time" for SpO<sub>2</sub> measurements is fixed at eight seconds.

7. If patient movement interferes with measurements, consider the following possible solutions:
  - be sure the sensor is secure and properly applied
  - use a new sensor with fresh adhesive backing
  - select a different type of sensor
  - move the sensor to a less active site

## Perform SpO<sub>2</sub> Monitoring with Nellcor Option

1. Attach the sensor to the patient according to the sensor manufacturer's instructions, observing all warnings and cautions.



**Warning** Use only Nellcor accessories and sensors with the monitor with Nellcor SpO<sub>2</sub> option as listed in the Welch Allyn *Products and Accessories* booklet (810-0409-XX).

2. When using a sensor extension cable, inspect the cable before use. Replace it if it shows any signs of wear, breakage, or fraying. Plug the sensor into the cable and plug the cable into the Propaq monitor, or plug the sensor directly into the monitor.
3. If the monitor SpO<sub>2</sub> receptacle has a locking ring, lock the connector in place by turning the locking ring clockwise until it stops. For other connectors, make sure the plug is all the way in.



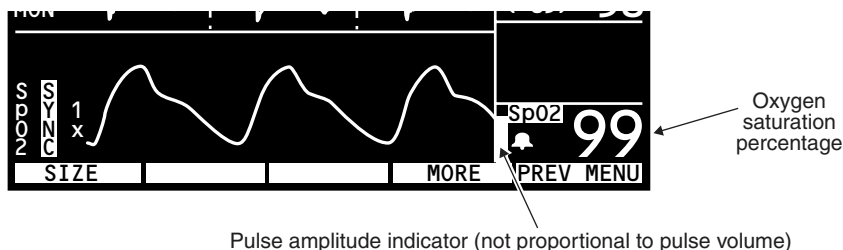
**Caution** If you see the error message **DEFECTIVE SpO<sub>2</sub> SENSOR**, either the sensor is not compatible with the monitor or the sensor is not working properly. Check the Welch Allyn *Products and Accessories* booklet (810-0409-XX) to be sure the sensor is compatible. If compatibility is not a problem, try another sensor.



The monitor displays STANDBY in the SpO<sub>2</sub> numeric window until it measures and displays the SpO<sub>2</sub> value. As oxygen saturation increases and decreases, the pitch of the heart tone rises and falls.

The Nellcor SpO<sub>2</sub> option periodically performs an internal adjustment which causes the SpO<sub>2</sub> waveform to appear flat for a brief period.

- From the Main Menu, press **SpO<sub>2</sub>** (or **SpO<sub>2</sub>/CO<sub>2</sub>**, then **SpO<sub>2</sub>**) to display the first SpO<sub>2</sub> menu similar to the following:



- Press **SIZE** to adjust the waveform size for best viewing (1x, 2x, 4x, or 8x).



At high magnification (4x, 8x), some waveforms may appear truncated. To view these waveforms, reduce the size until the complete waveform appears.

- Adjust the placement of the sensor until a good SpO<sub>2</sub> waveform is displayed. A waveform with artifact may cause erroneous oxygen saturation readings.
- Press **MORE** to display the second SpO<sub>2</sub> menu:



- Press **RESPONSE** to select the appropriate time required to measure SpO<sub>2</sub>:

Response	Indications for Use
NORMAL: 5-7 seconds	Use for relatively stable patients.
FAST: 2-3 seconds	Use when patient movement or other artifact is not present.
SLOW: 10-15 seconds <sup>a</sup>	Use when patients exhibiting movement are preventing accurate measurement at NORMAL setting.

a. SLOW setting is not applicable to the Nellcor SpO<sub>2</sub> option with motion tolerance.

9. If the C-LOCK function is desired, press **C-LOCK** to set it to ON.



C-LOCK synchronizes the pulse oximeter's systole determination to the R-wave to reduce the effects artifact may have on SpO<sub>2</sub> measurements. Under some conditions you may find more stable SpO<sub>2</sub> readings with C-LOCK set to ON. SYNC appears next to the waveform when synchronization to the ECG has been obtained. Synchronization takes a few seconds to establish the first time. If C-LOCK is on and the HR source is SpO<sub>2</sub>, the heart rate source is automatically changed to ECG. An ECG signal must be present or C-LOCK does not activate.

If you get false SpO<sub>2</sub> alarms with patients with low perfusion states or multiple arrhythmias, try turning off C-LOCK.

10. Set alarm limits according to your hospital's standards.



To help minimize false alarms, the Propaq monitor briefly delays or "holds off" triggering both audible and visual alarms for limit violations for SpO<sub>2</sub>% and Pulse Rate for 10 seconds. After the alarm hold-off period begins, if the monitor detects that the patient's vital sign has returned to acceptable limits, the monitor cancels the alarm hold-off. The next time a vital sign limit is violated, the monitor starts a new hold-off period.

11. If patient movement interferes with measurements, consider the following possible solutions:

- be sure the sensor is secure and properly applied
- use a new sensor with fresh adhesive backing
- select a different type of sensor
- move the sensor to a less active site



## Perform SpO<sub>2</sub> “Spot-Check” Monitoring

The SpO<sub>2</sub> Standby Mode allows you to remove the SpO<sub>2</sub> sensor from a patient without having to disable all alarms or disconnect the SpO<sub>2</sub> sensor cable from the Propaq monitor. You can therefore perform intermittent or “spot-check” SpO<sub>2</sub> monitoring.

1. While monitoring SpO<sub>2</sub>, remove the SpO<sub>2</sub> sensor from the patient, but leave it connected to the monitor. When the monitor detects the lack of a pulsatile waveform, it sounds a patient alarm and displays this menu:



2. Press **STANDBY** to place SpO<sub>2</sub> into the Standby Mode.



The monitor suspends the SpO<sub>2</sub> alarm tone indefinitely and displays **STANDBY** in place of SpO<sub>2</sub> numerics. SpO<sub>2</sub> remains in the Standby Mode until the SpO<sub>2</sub> sensor is reapplied to a patient. Other vital sign monitoring is not restricted. By contrast, if you press **SUSPEND** instead of **STANDBY**, the monitor temporarily suspends all alarm tones; however, the alarm tone resumes after 90 seconds if the SpO<sub>2</sub> sensor is still disconnected from the patient.

3. To resume SpO<sub>2</sub> monitoring, reapply the SpO<sub>2</sub> sensor to a patient.



The monitor exits the Standby Mode and resumes SpO<sub>2</sub> monitoring

The message **STBY** on the SpO<sub>2</sub> trend display and trend printouts indicates the monitor was in the SpO<sub>2</sub> Standby Mode.

## SpO<sub>2</sub> Messages

The following status message can appear in the equipment alert window or the SpO<sub>2</sub> numeric display area:

**NO SENSOR DETECTED** appears in the equipment alert window and indicates a probe has been disconnected from the monitor after being plugged in for more than a few seconds.

**SEARCH:** during the search time, the SpO<sub>2</sub> channel tries to detect blood pulsing through the measurement site. Once the measurement has been established, the oxygen saturation value is displayed in the numeric window.

**STANDBY** is displayed in the numeric window when the SpO<sub>2</sub> sensor is disconnected from the patient, an alarm occurs, and you press the **STANDBY** button. **STANDBY** is also displayed if you first plug the SpO<sub>2</sub> sensor cable into the monitor connector before attaching the SpO<sub>2</sub> sensor to the patient.

**DEFECTIVE SpO<sub>2</sub> SENSOR.** If you see the error message **DEFECTIVE SpO<sub>2</sub> SENSOR**, either the sensor is not compatible with the monitor or the sensor is not working properly. Check the Welch Allyn Products and Accessories booklet (810-0409-XX) to be sure the sensor is compatible. If compatibility is not a problem, try another sensor.

## Capnography (CO<sub>2</sub>)

### Intended Use

The Propaq's Capnography (CO<sub>2</sub>) option is intended to noninvasively measure the following vital signs or events: End-tidal CO<sub>2</sub> (ETCO<sub>2</sub>), Inspired CO<sub>2</sub> (INCO<sub>2</sub>), Breath Rate, and Apnea.

The CO<sub>2</sub> option is available as Mainstream CO<sub>2</sub>, Sidestream CO<sub>2</sub>, or Dualstream CO<sub>2</sub>. Although Dualstream CO<sub>2</sub> provides both Mainstream CO<sub>2</sub> and Sidestream CO<sub>2</sub> monitoring, only one method can be used at a time. The v option is required for any Propaq equipped with the CO<sub>2</sub> option.

CO<sub>2</sub> reading accuracy is affected by the presence of interfering gases and vapors. If the CO<sub>2</sub> option is used on patients who are being administered oxygen (O<sub>2</sub>) or nitrous oxide (N<sub>2</sub>O), be sure to set the appropriate compensation setting using the **GAS COMP** button.



**Warning** Before you use a Propaq on a new patient, always turn it off for a few seconds, then turn it on again. This clears the prior patient's trend values, alarm limit settings, and NIBP cuff inflation target.

Avoid exposing a Propaq with the CO<sub>2</sub> option to non-patient sources of CO<sub>2</sub> such as vehicle engine exhaust or smoke. When such exposure is possible, avoid opening the printer door. Exposure to these CO<sub>2</sub> sources can temporarily trap v within the monitor or Mainstream v sensor housing, even when monitor power is off. This can temporarily cause an erroneous elevated CO<sub>2</sub> measurement baseline until the trapped CO<sub>2</sub> leaks out and the baseline returns to zero (which can require as long as 3-24 hours).



**Caution** The Mainstream v and Sidestream CO<sub>2</sub> options are not recommended for use during magnetic resonance imaging (MRI) procedures. The magnetic fields involved will permanently damage the CO<sub>2</sub> sensor.



CO<sub>2</sub> monitoring outside the specified operating temperature range can cause inaccurate CO<sub>2</sub> readings. The operating temperature range for the CO<sub>2</sub> option is different than the range of 0° to 40° C for other Propaq Encore functions:

Mainstream CO<sub>2</sub> operating temperature: 10° to 40° C

Sidestream CO<sub>2</sub> operating temperature: 5° to 40° C

## Mainstream CO<sub>2</sub> Option

The Mainstream CO<sub>2</sub> option measures the carbon dioxide content of a patient's inhaled and exhaled breath. A mainstream sensor is attached to an airway adapter in series with a ventilator's patient breathing circuit.

Patients using Mainstream CO<sub>2</sub> must either be intubated or breathing through a well-fitting face mask connected to a breathing system such as an anesthesia circle system.

## Sidestream CO<sub>2</sub> Option

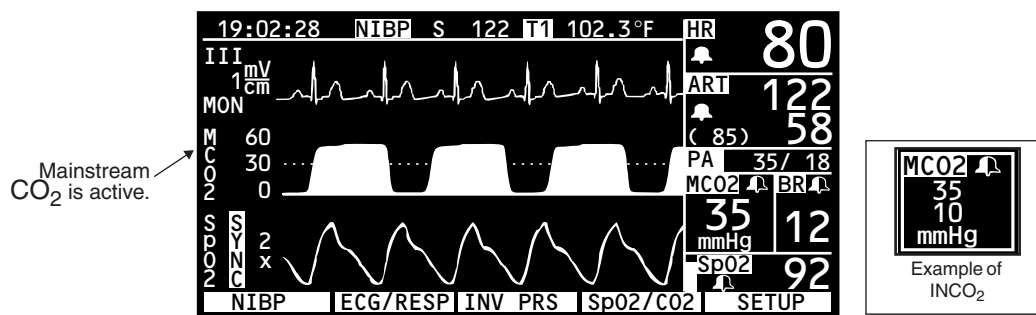
The Sidestream CO<sub>2</sub> option measures the carbon dioxide content of a patient's inhaled and exhaled breath. A sidestream sensor is located within the Propaq monitor. The patient's expired gas is aspirated from the airway and sent through a sampling line and sidestream watertrap to the internal sensor.

Patients using Sidestream CO<sub>2</sub> can either be intubated or non-intubated using a CO<sub>2</sub> Sampling cannula or a combination CO<sub>2</sub> Sampling/Oxygen Delivery nasal cannula.


## CO<sub>2</sub> Measurements and Display

The measured CO<sub>2</sub> levels are normally displayed as a waveform and an ETCO<sub>2</sub> numeric value. (The CO<sub>2</sub> waveform can be viewed when it is selected for display and the higher-priority IBP waveforms are not displayed.) If the INCO<sub>2</sub> numeric value is at an alarm level greater than 7.5 mmHg (or 1 kPa or 1%), it is also displayed.

If Mainstream CO<sub>2</sub> is active, the Propaq displays **MCO2**. If Sidestream CO<sub>2</sub> is active, the Propaq displays **SCO2**. The Propaq displays **CO2** if either or both the mainstream sensor and sidestream watertrap are installed in the Propaq but neither are active.



## Capnography (CO<sub>2</sub>)

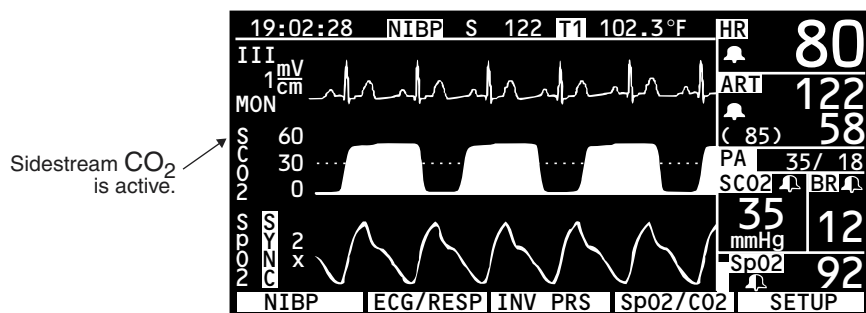

Displayed values of ETCO<sub>2</sub> and INCO<sub>2</sub> are the highest and lowest values (respectively) of CO<sub>2</sub> measured during the time interval set by the RESPONSE setting on the CO<sub>2</sub> Menu.

You can set upper and lower alarm limits for ETCO<sub>2</sub>, and an upper alarm limit for INCO<sub>2</sub>. The monitor only displays the numeric value for inspired CO<sub>2</sub> (INCO<sub>2</sub>) if it is in alarm or if it is greater than or equal to 7.5 mmHg (or ≥1 kPa or 1%). Refer to the CO<sub>2</sub> specifications in the Propaq Encore Reference Guide for more information.

### Breath Rate Measurements

Breath Rate (BR) is determined from the CO<sub>2</sub> sensor. The Propaq displays a numeric BR value next to the CO<sub>2</sub> values. You can set upper and lower alarm limits for BR.

### Apnea Events

In the Adult and Pediatric Mode, you can set the apnea delay to 6, 10, 15, 20, 25, or 30 seconds. In the Neonatal Mode, you can set the apnea delay to 6, 10, 15, or 20 seconds. The Propaq initiates an alarm in response to each apnea event longer than the apnea delay setting.

When an apnea event is detected, the BR numeric automatically goes to 0 and an apnea alarm occurs. After the alarm ceases, the Propaq prints an Apnea Ticket if the Apnea Ticket setting in the Printer Setup window is set to ON.

### Numeric Area Status Messages

The following status messages can appear in the numeric display area:

**OFF** indicates no CO<sub>2</sub> source is selected.

**SRCH** indicates the MCO<sub>2</sub> or SCO<sub>2</sub> sensor is preparing for a measurement.

**UNCAL** indicates the monitor has detected a problem such as a lack of calibration, an obstruction, or a low battery.

**WARM UP** indicates Mainstream CO<sub>2</sub> has been activated and is preparing for operation. This typically requires 30 seconds at room temperature.

**START UP** indicates Sidestream CO<sub>2</sub> has been activated and is preparing for operation. This typically requires 30 seconds at room temperature.

## CO<sub>2</sub> Display Menus and Status Window

To access the first CO<sub>2</sub> menu, press **SpO<sub>2</sub>/CO<sub>2</sub>, CO<sub>2</sub>** from the Main Menu

RANGE	mm/s		MORE	PREV MENU
-------	------	--	------	-----------

**RANGE** Selects the CO<sub>2</sub> waveform scale (range)

**mm/s** Sets the display sweep speed for CO<sub>2</sub> and RESP

### Mainstream CO<sub>2</sub> Menu and Status Window

When Mainstream CO<sub>2</sub> is active, press **MORE** from the first CO<sub>2</sub> menu to access the MCO<sub>2</sub> menu and the MCO<sub>2</sub> status window:

CO <sub>2</sub>			( 85)	58
GAS COMPENSATION: OFF			PA	35/ 18
RESPONSE : NORMAL			MCO <sub>2</sub>	BR
CO <sub>2</sub> SOURCE : MAINSTREAM			35	12
SWEEP SPEED : 6.25 mm/s			mmHg	
BAROMETER : 762.0 mmHg			SpO <sub>2</sub>	92
GAS COMP	RESPONSE	SOURCE	PREV MENU	

**GAS COMP** Selects the measurement compensation for CO<sub>2</sub> measurements

**RESPONSE** Sets the response time for CO<sub>2</sub> measurement (NORMAL, FAST, OR SLOW)

**SOURCE** Changes between Mainstream CO<sub>2</sub> and Sidestream CO<sub>2</sub> monitoring (if both options are installed), or disables CO<sub>2</sub> monitoring

Setting CO<sub>2</sub> SOURCE to OFF allows you to disable CO<sub>2</sub> monitoring without removing the sensor. When the CO<sub>2</sub> SOURCE is set to OFF, the Propaq displays OFF for CO<sub>2</sub> numerics.

## Capnography (CO<sub>2</sub>)

SWEEP SPEED	The CO <sub>2</sub> and RESP display sweep speed (set by the <b>mm/s</b> button)
BAROMETER	The ambient barometric pressure



To change CO<sub>2</sub> units, use the Settings Menu under the Service Menus. Changing these units will clear CO<sub>2</sub> trends.

## Sidestream CO<sub>2</sub> Menu and Status Window

When Sidestream CO<sub>2</sub> is active, press **MORE** from the first CO<sub>2</sub> menu to access the SCO<sub>2</sub> menu and the SCO<sub>2</sub> status window:

CO2	( 85 )	58
GAS COMPENSATION:	OFF	PA 35/ 18
RESPONSE	: NORMAL	SCO2 35 BR 12
CO2 SOURCE	: SIDESTREAM	mmHg
SWEEP SPEED	: 6.25 mm/s	SpO2 92
BAROMETER	: 762.0 mmHg	
FLOWRATE	: 175 mL/min	
GAS COMP	RESPONSE	SOURCE FLOWRATE PREV MENU

<b>GAS COMP</b>	Selects the measurement compensation for CO <sub>2</sub> measurements
<b>RESPONSE</b>	Sets the response time for CO <sub>2</sub> measurement (NORMAL, FAST, OR SLOW)
<b>SOURCE</b>	Changes between Mainstream CO <sub>2</sub> and Sidestream CO <sub>2</sub> monitoring (if both options are installed), or disables CO <sub>2</sub> monitoring  Setting CO2 SOURCE to OFF allows you to disable CO <sub>2</sub> monitoring without removing the sensor. When the CO2 SOURCE is set to OFF, the Propaq displays OFF for CO <sub>2</sub> numerics.
<b>SWEEP SPEED</b>	The CO <sub>2</sub> and RESP display sweep speed (set by the <b>mm/s</b> button)
<b>BAROMETER</b>	The ambient barometric pressure
<b>FLOWRATE</b>	Sets the sampling flow rate to either 90 or 175 ml/min. You can change the flow rate while Sidestream CO <sub>2</sub> is active



To change CO<sub>2</sub> units, use the Settings Menu under the Service Menus. Changing these units will clear CO<sub>2</sub> trends.

## Mainstream CO<sub>2</sub> Monitoring

1. Select the appropriate airway adapter.

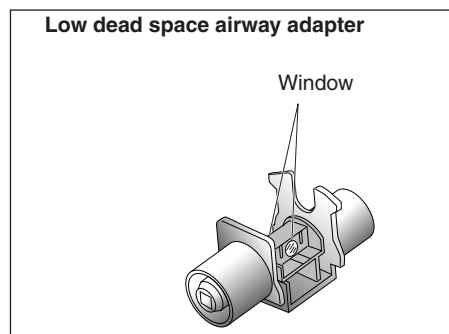
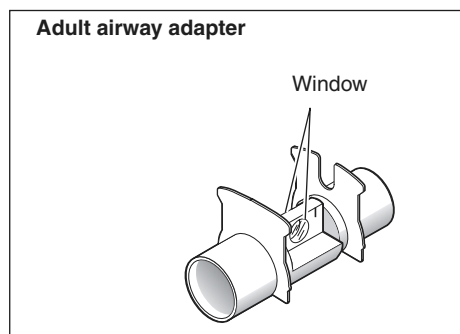
Use only accessories provided or recommended in the Welch Allyn *Products and Accessories* booklet. Always refer to the manufacturer's Directions for Use for instructions about operation, cleaning, and replacement.



**Warning** Do not attempt to verify operation of the CO<sub>2</sub> sensor by blowing through it directly. Always blow through an attached airway adapter. Otherwise, a small amount of CO<sub>2</sub> from your breath may enter the CO<sub>2</sub> sensor housing and cause a small shift in the measured CO<sub>2</sub> values. It may take 3-24 hours for the sensor to return to proper calibration.

Do not clean and/or reuse a single-patient-use airway adapter. When a single-patient-use airway adapter becomes occluded, replace it.

2. Connect the adapter, ventilator circuit, and CO<sub>2</sub> sensor according to the manufacturer's instructions.



**Warning** Prior to using an airway adapter, always look through the window lumen and inspect the adapter for inadvertently lodged obstructions and for window integrity.

If the sensor does not easily slide onto the adapter, do not attempt to force these components together. They fit together in only one way. Take care not to damage the glass windows.

After attaching the sensor to the adapter, check the adapter again for proper placement of the sensor and for window integrity.





**Warning** When attaching the airway adapter, position the adapter so the sensor is on top to avoid fluid collection in the sensor airway slot. Any concentration of fluids here can cause inaccurate CO<sub>2</sub> readings.

Always check to make sure there are no leaks in the breathing circuit. Check all of the connections.

3. Plug in the CO<sub>2</sub> sensor cable to the Mainstream CO<sub>2</sub> connector on the Propaq left side panel.
4. Set up the CO<sub>2</sub> channel and alarms. Follow the steps described on page 81.



**Warning** When disconnecting the CO<sub>2</sub> sensor from the tracheal or endotracheal tube, check the sensor to determine how hot it is. If it is too hot for patient comfort, do not allow it to come into contact with the patient.



When disconnecting the airway adapter from the ventilator circuit, always detach the CO<sub>2</sub> sensor from the airway adapter before removing the airway adapter from the ventilator circuit.

## How to Set Up the CO<sub>2</sub> Channel and Alarm Limits

After you connect the Mainstream CO<sub>2</sub> sensor or Sidestream CO<sub>2</sub> watertrap, follow these steps to set up the CO<sub>2</sub> channel and alarm limits.



After you connect a Mainstream CO<sub>2</sub> sensor or Sidestream CO<sub>2</sub> watertrap, the Propaq displays the waveform briefly without a value range. It displays WARM UP (for Mainstream) or START UP (for Sidestream) in the CO<sub>2</sub> numerics window. After about 20 seconds, the Propaq displays the CO<sub>2</sub> measurement and waveform range.

1. Press **SpO<sub>2</sub>/CO<sub>2</sub>, CO<sub>2</sub>** to display the first CO<sub>2</sub> Menu.
2. Press **RANGE** until you see the desired waveform scale range on the Propaq screen. The range choices are shown below.

Units	Waveform Scale Range		
mmHg	0-100	0-60 (default setting)	0-30
kPa	0-14	0-8	0-4
%	0-14	0-8	0-4



If an inspired value is displayed indicating patient rebreathing (non-zero INCO<sub>2</sub>), check the patient breathing circuit for proper function. Verify the sensor calibration against room air. If the Propaq continues to display inspired values, return the sensor to Welch Allyn for service.

3. Press **mm/s** to select either 3.13, 6.25 or 12.5 mm/sec. The default setting is 6.25 mm/sec. This sweep speed also applies to RESP.
4. Press **MORE** to view the CO<sub>2</sub> status window.
5. If either O<sub>2</sub> or N<sub>2</sub>O is being administered to a patient, press **GAS COMP** to set the proper gas compensation as listed below. If no other gas is being administered, set GAS COMPENSATION to OFF. (OFF is the default setting.)

Gas Administration Level/GAS COMP Setting	ETCO <sub>2</sub> or INCO <sub>2</sub> Value
OFF	CO <sub>2</sub> value = actual CO <sub>2</sub> value
O <sub>2</sub> > 50%, No N <sub>2</sub> O	CO <sub>2</sub> value = actual CO <sub>2</sub> value x 1.03
N <sub>2</sub> O > 50%	CO <sub>2</sub> value = actual CO <sub>2</sub> value x 0.952



If the ETCO<sub>2</sub> value is displayed as + + +, verify calibration against a known reference gas. If the sensor calibration is not accurate, return it to Welch Allyn for service.

6. Press **RESPONSE** to select either NORMAL, SLOW, or FAST.

The FAST setting is recommended where a sudden step change in ETCO<sub>2</sub> is of concern, such as that induced by an air embolus in certain neurosurgical procedures.

A SLOW response will decrease ETCO<sub>2</sub> false alarms when breath morphology varies considerably from one breath to the next.

The default setting is NORMAL.

Response Time Setting	Sampling Period	Typical Indications for Use
<b>FAST</b>	15 seconds	During neuroanesthesia
NORMAL	30 seconds	During routine use
<b>SLOW</b>	45 seconds	To decrease ETCO <sub>2</sub> false alarms

7. Set alarm limits for ETCO<sub>2</sub>, INCO<sub>2</sub>, and Breath Rate.



**Warning** For patient safety, it is recommended that the Breath Rate alarm limits always be turned on and set appropriately.

Motion artifact or other factors can cause false RR/BR alarms. To help minimize false alarms, the Propaq delays or “holds off” triggering a RR/BR alarm for 5 seconds. After this holdoff period begins, if the Propaq detects that the patient’s RR/BR vital sign has returned to acceptable limits, the Propaq cancels the alarm holdoff. The next time a RR/BR limit is violated, the Propaq starts a new 5-second RR/BR alarm holdoff period.

8. Set the alarm limit for Apnea Delay—this is the maximum time allowed between two successive breaths before an Apnea alarm occurs.

After the first breath has been detected, the Apnea alarm limit setting will be automatically turned on for as long as the CO<sub>2</sub> channel is active. The RR/BR and Apnea Alarm Limit window is shown below. STAT SET does not affect the Apnea limit.

ALARM LIMITS		UPPER	LOWER	( 85 )	58
RR/BR	Br/M:	30	5	PA	35/ 18
APNEA	DELAY:	15	secs	MC02	45
				mmHg	12
				SpO2	92
NEXT	UP	DOWN	ON/OFF	NEXT	PAGE

## Sidestream CO<sub>2</sub> Monitoring



**Warning** If the Sidestream CO<sub>2</sub> option is connected to a ventilatory circuit, be sure to adjust appropriate ventilator or anesthesia system settings to compensate for the sampling flow volume (90 or 175 ml/min) that is aspirated from the ventilatory circuit by the Sidestream CO<sub>2</sub> option.

Do not use Sidestream CO<sub>2</sub> if flammable anesthetic gases are in use.

Setting up for Sidestream CO<sub>2</sub> monitoring requires these major steps:

1. Connect the watertrap.
2. Set up the CO<sub>2</sub> channel and alarm limits.

3. Connect to a non-intubated patient.

OR

4. Connect to an intubated patient.

Use only accessories provided or recommended in the Welch Allyn *Products and Accessories* booklet. Always refer to the manufacturer's Directions for Use.



Breath rates greater than 50 breaths/minute may reduce the reported ET<sub>CO<sub>2</sub></sub> values. Select the 175 ml/min flow rate to minimize errors at higher breath rates.

The 175 ml/min flow rate is recommended for intubated adult patients.

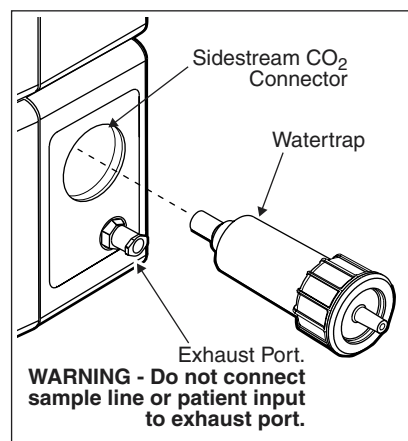
When monitoring a small child with a rapid respiratory rate, Mainstream CO<sub>2</sub> can provide a more accurate representation of the expired CO<sub>2</sub> waveform.

### Connect the Watertrap

1. Firmly insert the Sidestream CO<sub>2</sub> watertrap into the Sidestream CO<sub>2</sub> connector on the Propaq left side panel.



**Warning** The watertrap is disposable and should only be used for a single patient. Do not reuse the watertrap for another patient.



### Set Up the CO<sub>2</sub> Channel and Alarm Limits

1. Follow the steps described on page 81.

## Connect to a Non-Intubated Patient

1. Position the cannula on the patient according to the manufacturer's instructions.



**Warning** The cannula is disposable and should only be used for a single patient. Do not reuse the cannula for another patient.

If oxygen is being delivered while using Sidestream CO<sub>2</sub>, be sure to use a CO<sub>2</sub> Sampling and O<sub>2</sub> Delivery Cannula. Using a different type of cannula could obstruct oxygen delivery.

The exhaust port for Sidestream CO<sub>2</sub> is an output for the expired gases from the patient and any connected breathing apparatus. The exhaust port is intended only for connection to gas collection equipment such as gas scavenger devices. **Do not allow any other connection to the exhaust port.**

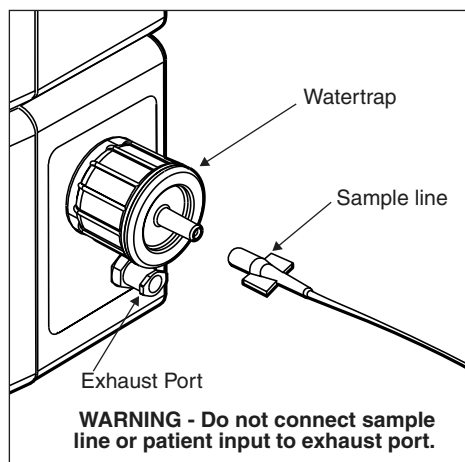


If you use a gas scavenging system with Sidestream CO<sub>2</sub>, be sure to install it according to the manufacturer's instructions. The scavenging system should comply with ISO 8835-3:1997 (E).

2. Connect the sample line to the cannula and the watertrap. Make sure that the sample line is firmly connected.



Sidestream CO<sub>2</sub> accuracy decreases if additional tubing is connected to the sample line. Avoid connecting additional tubing to the standard sample line.



## Connect to an Intubated Patient

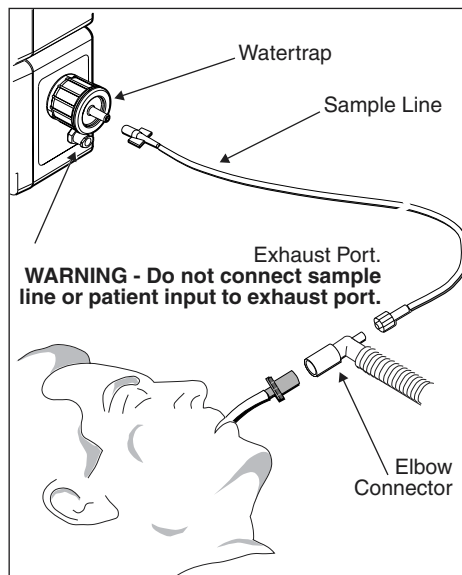
1. Connect the gas sampling elbow and elbow connector into the patient's breathing circuit according to the manufacturer's instructions.



**Warning** The exhaust port for Sidestream CO<sub>2</sub> is an output for the expired gases from the patient and any connected breathing apparatus. The exhaust port is intended only for connection to gas collection equipment such as gas scavenger devices. **Do not allow any other connection to the exhaust port.**



If you use a gas scavenging system with Sidestream CO<sub>2</sub>, be sure to install it according to the manufacturer's instructions. The scavenging system should comply with ISO 8835-3:1997 (E).



2. Connect the sample line to the elbow connector and the watertrap. Make sure the sample line is firmly connected.



Sidestream CO<sub>2</sub> accuracy decreases if additional tubing is connected to the sample line. Avoid connecting additional tubing to the standard sample line.

## CO<sub>2</sub> Messages

### Mainstream CO<sub>2</sub> Messages

Equipment messages for the Mainstream CO<sub>2</sub> option can appear on the display in an equipment alert window and in numeric zones.

If a sensor is damaged, contact Welch Allyn' Technical Services Department for information on sensor service options.

**ALTIMETER FAILURE - RANGE.** The Propaq is operating at an altitude outside the Mainstream CO<sub>2</sub> option's operating altitude range of -2,000 to 15,000 feet. Returning the monitor to within this range automatically cancels this message and restores operation.

**ALTIMETER FAILURE - RATE.** The altimeter has detected that the ambient pressure is changing at a rate greater than 100 mmHg/minute. When the rate of change is back within the 100 mmHg/minute range, disconnect and reconnect the CO<sub>2</sub> sensor to the Propaq.

**DEGRADED WAVEFORM, CHECK ADAPTER (UNCAL appears in the numerics area).** The Mainstream CO<sub>2</sub> adapter is obstructed or the CO<sub>2</sub> sensor has failed. The CO<sub>2</sub> waveform is displayed without range values. Replace the adapter or replace the sensor.

**LACK OF WAVE, CHECK ADAPTER, SENSOR.** Either the airway adapter is obstructed or the CO<sub>2</sub> sensor has failed. Replace the airway adapter if it is obstructed. The sensor must be unplugged and plugged in again.

**LOW BATTERY, HEATER DISABLED (UNCAL appears in the numerics area).** The Propaq's battery voltage is too low. The CO<sub>2</sub> waveform is displayed without range values. To continue operation, plug the ac power adapter into the Propaq.

**NO MAINSTREAM SENSOR DETECTED (SRCH appears in the numerics area).** The Mainstream CO<sub>2</sub> sensor has been disconnected from the Propaq after providing CO<sub>2</sub> values. Disconnect and reconnect the sensor to the Propaq if necessary.

**NON-PROTOCOL SENSOR (UNCAL appears in the numerics area).** A CO<sub>2</sub> sensor has been connected that does not match Welch Allyn's specifications. The CO<sub>2</sub> waveform is displayed without range values. Replace the sensor with a Welch Allyn CO<sub>2</sub> sensor.

**SENSOR FAILURE, CALIBRATION ERROR.** A sensor is defective or out of calibration. The sensor will be disabled. Replace the sensor.

**SENSOR FAILURE - EEPROM.** The sensor has failed. Replace the sensor.

**SENSOR FAILURE - HEATER.** The sensor's temperature control circuit or the Propaq's CO<sub>2</sub> circuitry has failed. Try replacing the sensor. If the message reappears, have the Propaq serviced.

**SENSOR FAILURE - MOTOR DRIVE.** The sensor's motor drive (in the sensor head) has failed. Replace the sensor.

**SENSOR TEMPERATURE TOO HIGH.** The sensor's temperature is too high. The sensor's ambient operating range is 10° to 46° C. When the ambient temperature returns to this range, this message is automatically removed and operation is restored.

**WARM UP or WARM (appears in the numerics area).** The sensor heater is warming up. Wait 20 to 30 seconds for the sensor to heat. Values should appear in the numerics area when the sensor is sufficiently warm.

## ***Sidestream CO<sub>2</sub> Messages***

Equipment messages for the Sidestream CO<sub>2</sub> option can appear on the display in an equipment alert window.

**ALTIMETER FAILURE - RANGE.** The Propaq is operating at an altitude outside the Sidestream CO<sub>2</sub> option's operating altitude range of -2,000 to 15,000 feet. Returning the monitor to within this range automatically cancels this message and restores operation.

**ALTIMETER FAILURE - RATE.** The altimeter has detected that the ambient pressure is changing at a rate greater than 100 mmHg/minute. When the rate of change is back within the 100 mmHg/minute range, disconnect and reconnect the CO<sub>2</sub> sensor to the Propaq.

**ALTIMETER NOT CALIBRATED - EEPROM -** The Sidestream CO<sub>2</sub> option has not been calibrated. Refer the Propaq to a Biomedical Engineer for calibration.

**AMBIENT TEMPERATURE TOO HIGH.** The sensor temperature is too high. The Sidestream CO<sub>2</sub> option is disabled until the ambient temperature is within the operating range specifications.

**AMBIENT TEMPERATURE TOO LOW.** The sensor temperature is too low. The Sidestream CO<sub>2</sub> option is disabled until the ambient temperature is within the operating range specifications.

**CALIBRATION ERROR, SERVICE REQUIRED.** Send the Propaq to a Biomedical Engineer for service.

**DEGRADED WAVEFORM, SERVICE REQUIRED.** Send the Propaq to a Biomedical Engineer for service.

**LACK OF WAVEFORM, SERVICE REQUIRED.** Send the Propaq to a Biomedical Engineer for service.

**MOTOR FAILURE, SERVICE REQUIRED.** The sensor hardware has failed. Send the Propaq to a Biomedical Engineer for service.

**NO WATERTRAP DETECTED.** There is no Sidestream CO<sub>2</sub> watertrap installed. Install a watertrap.

**OCCCLUSION - CHECK EXHAUST PORT/TUBING.** Blockage has been detected on the pneumatic exhaust port. Check the exhaust port and related tubing for occlusions. Make sure that the sampling line and any inputs to the patient breathing apparatus are not connected to the exhaust port.



**OCCLUSION - CHECK WATERTRAP/TUBING.** Blockage has been detected on the Sidestream CO<sub>2</sub> input. Check the watertrap, sample line, and any connected tubing for occlusion.

**PUMP FAILURE, SERVICE REQUIRED.** The pump is not able to maintain the target flow rate. Send the Propaq to a Biomedical Engineer for service.

**SIDESTREAM STICK EEPROM FAILURE.** Send the Propaq to a Biomedical Engineer for service.

**SSP BOARD EEPROM FAILURE.** Send the Propaq to a Biomedical Engineer for service.



# Alarms and Limits

## Alarms Window and Menu

You can set the alarms in normal operating mode or in-service mode. You can also have the Propaq Encore automatically print the patient's vital sign information when an alarm occurs.

## Alarm Status Window

This status window shows, with a bell symbol, the alarm status of each patient channel.

Bells only appear when at least one limit is turned on and the vital sign parameter is being monitored. For NIBP, the bell always appears whenever at least one limit is turned on, even when an NIBP reading is not being taken.

- The full bell shows you that all alarm limits are turned on.
- The half bell indicates that at least one alarm limit is turned off.
- The absence of a bell shows that no alarm limits are turned on.

These bells also appear in the numerics windows.

## Alarms Menu

The Alarms Menu shown below the status window lets you access other alarm functions to individually set alarm limits or automatically set them.

ALARM STATUS		( 85)	58
APNEA :	🔔	RR/BR :	🔔
HR/PR :	🔔	SpO2 :	🔔
ART :	🔔	NIBP :	🔔
PA :	🔔	TEMP :	🔔
CO2 :	🔔		
STAT SET	SUSPEND	ALL ALRM	LIMITS
		PREV MENU	

## Adjusting the Alarm Tone Volume

You can adjust the alarm tone to one of three volumes. Except for temporarily suspending the alarm tone using the **SUSPEND** button, you cannot turn off the tone completely.

To adjust the alarm tone volume, press **SETUP**, then **MORE**. The status window appears. Use the **NEXT** button to select ALARM TONE, and the **CHANGE** button to change the setting.

SETUP	HR/PR	RR/BR	( 85)	58
CURRENT SOURCE :	ECG	MC02	PA	35/ 18
SELECTED SOURCE :	ART	---	MC02	BR
SWEEP (mm/s) :	12.5	3.13	35	12
ALARM TONE :	HIGH		mmHg	
HR/PR TONE :	LOW		SpO2	92
PATIENT MODE :	ADULT			
NEXT	CHANGE	WAVE SEL	PRINTER	MORE



**Warning** At the highest volume alarm level, the sound pressure level will not exceed safe limits (OSHA HSM 73-1101, 1972). However, additional precautions may be required in patients under treatment with ototoxic medications.

## Types of Propaq Alarms

### Patient Alarms: Definitions and Indications

The Propaq Encore can alert you to changing patient conditions through its programmable alarm functions. After you set alarm limits, any vital sign that violates any of its limits results in both audible and visual alarm indications.

The Propaq Encore also shows you when any alarm limit is turned off by illuminating the amber **ALARM(S) OFF** light.

Patient and Alarm Limit Status	Red ALARM Light	Yellow ALARM(S) OFF Light	Tone
Patient in alarm condition, and all alarm limits are on	FLASH	OFF	ON
Patient in alarm condition, and at least one alarm limit is off	FLASH	ON	ON
Patient alarms suspended (whether in alarm condition or not)	ON	ON or OFF	OFF
Patient not in alarm condition, and at least one alarm limit is off.	OFF	ON	OFF
Equipment alert, patient not in alarm condition	OFF	FLASH	ON for 1 second, OFF for 4 seconds
Equipment alert, patient alarms suspended	ON	FLASH	OFF

A steady, high-pitched alarm tone sounds whenever a patient alarm limit is violated.

The alarm tone continues until:

- the patient condition changes and no longer violates the limit
- you *suspend the alarm tone* by pressing the **SUSPEND** button
- you *adjust the alarm limit* so the vital sign does not violate it
- you *turn off* the violated alarm limit

The alarm tones are summarized in the table:

Alarm Condition	Tone
Patient alarm	Continuous ON
Apnea alarm	ON for 1 second, OFF for 1 second
Equipment alert	ON for 1 second, OFF for 4 seconds

## Life-Threatening Alarms

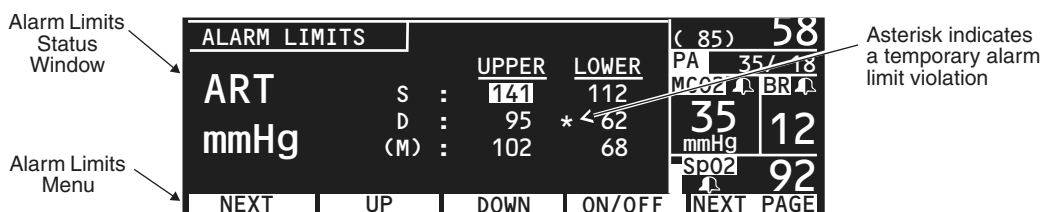
Life-threatening alarms are automatically turned on when the channel is activated, and cannot be turned off by the user. A life-threatening alarm is the highest priority patient alarm. You are notified of any life-threatening alarm as soon as it is detected. An apnea alarm is a life-threatening alarm.

## Temporary Alarms

A temporary alarm occurs when a patient's condition changes, causes an alarm to occur, changes again, and no longer violates a limit.

After a temporary alarm, review each vital sign's alarm limit window until you find the limit with an asterisk (\*) next to it. This limit caused the alarm. The asterisk(s) are cleared when you leave that limits window.

1. Press **ALARMS** and then **LIMITS** to display the first limits window.
2. Press **NEXT PAGE** (in the event of multiple alarm violations) until all asterisks are located.



## Responding to Patient Alarms

### *Suspend an Alarm Tone*

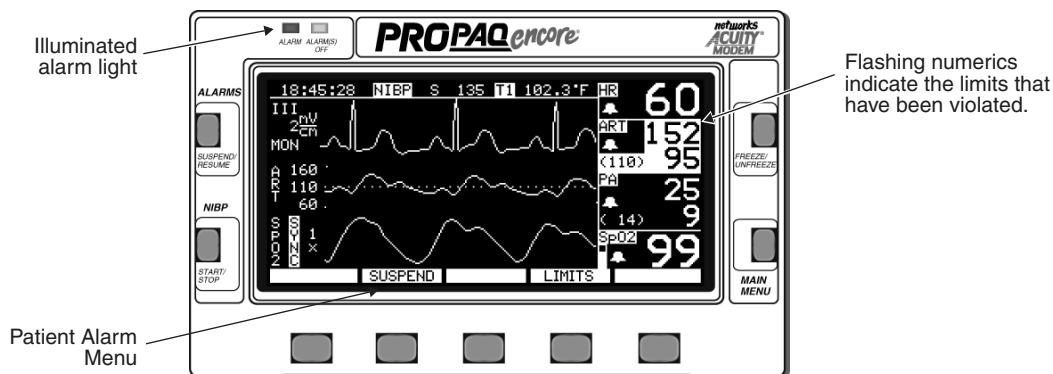
You can temporarily turn off an alarm tone by pressing the **SUSPEND/RESUME** key to the left of the screen or the **SUSPEND** button in the Patient Alarm Menu. The tone is suspended for 90 seconds. During that period, the visual alarm indications continue. When you suspend the alarm, the red ALARM light is on.

You can “unsuspend” the alarm before 90 seconds has elapsed by pressing the **SUSPEND/RESUME** key to the left of the screen, or the **RESUME** button in the Alarms Menu. If an alarm condition still exists, the tone will again sound. For NIBP (both pressure values and pulse rates derived from NIBP), pressing **RESUME** will not resume the NIBP alarm because NIBP is not continuously measured.



**Warning** Suspending an alarm suspends **all** alarm tones for 90 seconds or until the **RESUME** button is pressed.

If the **LIMITS** button is available, you can make Propaq Encore alarm limit adjustments by pressing **LIMITS**.



### *Turn Off or Change Alarm Limits*

Possible adjustments in order to cancel the alarm are:

- Turn off the alarm limits only for the violating alarm limit values by pressing **ON/OFF** in the Limits Menu after you have selected the appropriate vital sign limits window (**LIMITS, NEXT PAGE, ON/OFF**)



The **ON/OFF** button is not available for HR/PR alarm limits if HR/PR ALARM LIMITS in the Settings window is set to CANNOT TURN OFF.

- Manually change the violated alarm limit value by selecting the violated limit value and adjusting it
- Turn off all alarm limits by pressing **MAIN MENU, SETUP, ALARMS, ALL ALARM**

## Responding to Life-Threatening Alarms

During a life-threatening alarm, you can suspend the alarm tone and adjust the Alarm Limits.

## Alarm Prints

If you have the Printer option, you can set the Propaq Encore to print when a patient alarm occurs. See page 107 for more information about printing.



## Alarm Holdoffs

### Propaq Audio and Visual Alarm Holdoffs

Motion artifact or other factors can cause false vital sign alarms. To help minimize false alarms, the Propaq briefly delays or “holds off” triggering alarms for certain vital sign limit violations.

After the alarm holdoff period begins, if the monitor detects that the patient’s vital sign has returned to acceptable limits, the Propaq cancels the alarm holdoff. The next time a vital sign limit is violated, the Propaq starts a new alarm holdoff period. The following table summarizes the alarm holdoff periods.

Monitoring Function	Alarm Holdoff Time Period
HR/PR	3 seconds (except NIBP PR)
SpO <sub>2</sub>	10 seconds
RR/BR	5 seconds

When the Propaq is connected to the Acuity Central Monitoring System, Acuity does not receive any alarm indication from the monitor until this holdoff period expires.

### Propaq Audio Alarm Holdoff with Acuity

When a Propaq Encore in Adult or Pediatric Mode is connected to an Acuity System, the audio alarms at the bedside Propaq can be delayed up to 4 minutes and 15 seconds. The delay time is selected in Acuity software at the time of Acuity installation. Visual alarm indications are not delayed.



This audio alarm holdoff does *not* occur in the Neonatal mode.

## Setting Alarm Limits

### Adjust Limits with STAT SET

When you want to quickly set all alarm limits, the Propaq Encore can calculate new alarm limits using the patient's current values. Press the **STAT SET** button in the Alarms Menu (**SETUP, ALARMS, STAT SET**). The Propaq Encore activates all alarms and calculates the limits for all monitored vital signs, except apnea delay. The monitor performs a mathematical function (addition, subtraction, or multiplication) on the current value of the vital signs to arrive at the new limits. The formulas for these calculations are shown in the table.

Vital Sign	If the Patient's Vital Sign Value is	Then Calculated New Lower Limit is	Then Calculated <sup>a</sup> New Upper Limit is
Heart Rate	HR ≤ 99 100-250 HR ≥ 251	HR x 0.8 HR - 20 Unchanged	HR x 1.2 HR + 20 250
Pulse Rate	PR ≤ 99 PR ≥ 100	PR x 0.8 PR - 20	PR x 1.2 PR + 20
Invasive Pressure	Inv Prs ≤ 25 26 - 99 Inv Prs ≥ 100	Inv. Pressure - 5 Inv. Pressure x 0.8 Inv. Pressure - 20	Inv. Pressure + 5 Inv. Pressure x 1.2 Inv. Pressure + 20
NIBP	NIBP ≤ 25 26 - 99 NIBP ≥ 100	NIBP - 5 NIBP x 0.8 NIBP - 20	NIBP + 5 NIBP x 1.2 NIBP + 20
Respiration Rate/Breath Rate	RR/BR ≤ 25 26 - 99 RR/BR ≥ 100	RR/BR - 5 RR/BR x 0.8 RR/BR - 20	RR/BR + 5 RR/BR x 1.2 RR/BR + 20
Temperature	Temp ≥ 0°C	Temp - 0.5	Temp + 0.5
SpO <sub>2</sub>	SpO <sub>2</sub> ≥ 0%	SpO <sub>2</sub> - 5 (min. limit 50%)	100% (adult and pediatric mode) SpO <sub>2</sub> + 5 (neonate mode)
ETCO <sub>2</sub>	ETCO <sub>2</sub> ≥ 0 mmHg  ETCO <sub>2</sub> ≥ 2.0 (% or kPa)	ETCO <sub>2</sub> - 5 mmHg (min. 15 mmHg) ETCO <sub>2</sub> - 0.7 (% or kPa) (min. 2.0% or 2.0 kPa)	ETCO <sub>2</sub> + 10 mmHg  ETCO <sub>2</sub> + 1.4 (% or kPa)
INCO <sub>2</sub>	INCO <sub>2</sub> ≥ 0 mmHg INCO <sub>2</sub> ≥ 0 (% or kPa)	Not affected by STAT SET	INCO <sub>2</sub> + 5 mmHg INCO <sub>2</sub> + 0.7 (% or kPa)
Apnea Delay	Not affected by STAT SET		

a. New alarm limits calculated by STAT SET cannot be outside the allowable alarm limit range. If a new limit is calculated to be above or below the allowable alarm limit range, it defaults to the maximum or minimum alarm limit allowed for that vital sign.



**Warning** If a patient's vital sign value falls outside of the upper or lower alarm range limit, **STAT SET** turns off the alarm and the alarm limit except for the following:

1. The lower alarm limits for SpO<sub>2</sub> and ETCO<sub>2</sub> are not turned off by **STAT SET**.
2. If HR/PR ALARM LIMITS in the Settings window is set to CANNOT TURN OFF, **STAT SET** affects HR/PR alarm limits as follows: for an overrange HR/PR patient value (displayed as +++), the upper alarm limit is set to the maximum alarm limit, and the lower alarm limit is unchanged; for an underrange HR/PR patient value (displayed as - - -), the lower alarm limit is set to the minimum alarm limit and the upper alarm limit is unchanged; for an indeterminate HR/PR patient value (displayed as ???), the alarm limits are unchanged.

## Turning On and Off All Limits

If you want to turn all limits on or off without changing their values, press **ALL ALARM** in the Alarms Menu (a confirmation window appears requiring you to confirm your selection). The Alarms Status Window displays bells when alarms are on.

## Changing Individual Limits

Setting each limit requires selecting the desired vital sign parameter window, selecting the limit to adjust, and changing the limit.

1. From the Main Menu, press **SETUP, ALARMS**.
2. Press **LIMITS** to display the alarm limits window and the Limits Menu.
3. Press **NEXT PAGE** to change to the desired alarm limit window.
4. Press the **NEXT** button to move the cursor.
5. Press **UP** or **DOWN**, or **ON/OFF** (if **ON/OFF** is available) to set the limit to the desired limit value.
6. When the limit is set, select the next limit with the **NEXT** button. Or, to select another vital sign, press **NEXT PAGE**.

Repeat steps 4 through 6 to adjust the desired limits. Continue the process until you've set all the limits you want.

## Customizing Propaq Encore Alarm Limits

Although adjusting the Propaq's alarm limits can be automatic with **STAT SET**, you can also customize the limits you want to take effect every time you turn on the monitor. See “Setting the Current, Custom, and Power-up Modes” on page 33.

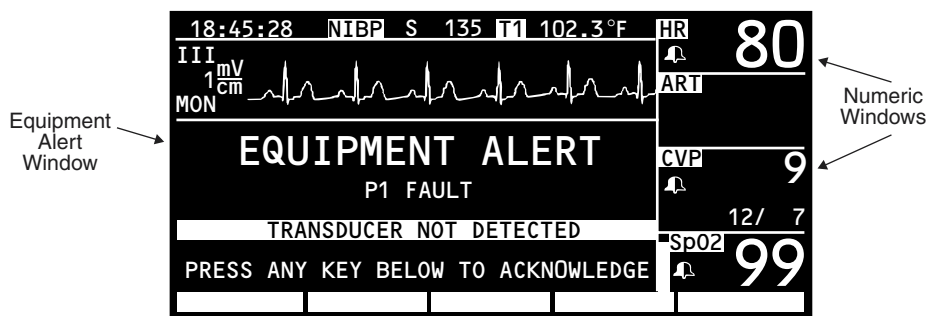
## Equipment Alerts

### Definitions and Indications

The Propaq Encore can alert you to changing equipment conditions, such as disconnected or faulty sensors, low battery voltage, lost programmed settings, and other conditions that can affect patient monitoring.

If an equipment alert condition is detected, a high-pitched alarm tone sounds for one second at five-second intervals. This alert tone pattern repeats until you respond to the equipment alert by pressing any button located at the bottom of the screen or until the equipment condition is corrected. You can also press **SUSPEND** and suspend the audio alarm tone for 90 seconds.

An equipment alert window also appears on the display identifying the condition. Certain equipment conditions also result in a caution message. Some of these messages alternate with the time of day, others appear in corresponding numeric windows.



Some equipment conditions resulting in an equipment alert also interrupt normal patient monitoring and cause a patient alarm. Patient alarms always take precedence over an equipment alert.

### Responding to Equipment Alerts

The equipment alert tone continues to sound and the window remains on the display until:

- the equipment condition is corrected
- you press any button located at the bottom of the screen

For some equipment alerts, all you need to do is correct the condition and the Propaq Encore resumes normal monitoring.

If the equipment condition also caused a patient alarm, you will need to first suspend the alarm tone by pressing **SUSPEND**, then take the required action.

The Propaq Encore returns the channel to its settings prior to the equipment alert. If a patient alarm occurred and you turned off any alarm limits, you will need to turn them back on. You can see which alarm limits are turned on by pressing **SETUP, ALARMS**.

## **Power-up Equipment Alert: Program Fault, Settings Lost**

If a PROGRAM FAULT: SETTINGS LOST, TIME/DAY RESET equipment alert appears when you turn on the monitor, the monitor cannot recall the programmed custom settings and current time and date. This can occur if the battery is drained or after new software has been installed.

If this occurs, the monitor provides a special sequence of display windows to help you regain use of your monitor as quickly as possible. Refer to page 21 and perform the described steps.

If you follow these steps and the equipment alert reappears at powerup, the monitor may need to be serviced and the battery replaced. Contact a qualified service person.

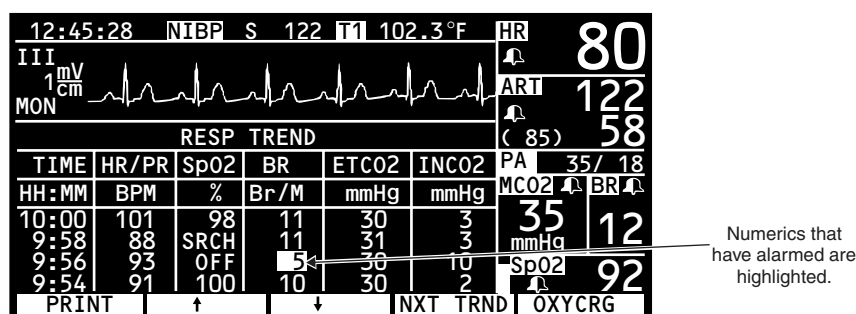
# Trends

Every two minutes, the Propaq Encore collects the monitored vital sign numerics and stores them in its trend memory, which can save the last five hours of trend information. (See “NIBP Trends” on the next page for exceptions.) All this information can be printed and viewed as a trend print.

## The Trend Status Window and Menu

There are five trends: NIBP, RESP, P1, P2, and TEMP. Except for NIBP, all vital signs are continuously monitored from the time monitoring begins to the time it ends.

Each trend shows the time of the reading, the HR/PR measurements, the SpO<sub>2</sub> value (if configured), and other values. The three blood pressure trends show systolic, diastolic, and mean pressures, and all but the TEMP trend show respiration rate, if it's available. Columns on the trend table show the word 'OFF' for the vital signs not being monitored.



The Trend Menu allows you to select trended data for display and print the displayed trend if a printer is attached.

- PRINT** Prints the displayed trend.
- ↑ ↓** The Up Arrow button allows you to scroll up to the most current reading and the Down Arrow lets you scroll down to the oldest reading, four readings at a time.
- NXT TRND** Allows you to cycle through the current display of each trend group.
- OXYCRG** Prints an oxycardiogram. For more information on OxyCRG, see page 112.

## How Trends are Accumulated

For each of the four continuous trends, numerics are sampled every two minutes to a maximum of 150 samples (up to five hours). When the maximum is reached, the trending continues but the older data is no longer stored. The readings are displayed in descending order, most recent first.

The following is programmable: which trend window comes up first. The NIBP trend window is the factory default.

To clear trend data before connecting a new patient, the monitor power must be cycled. This will prevent the trend data of a previous patient from being attributed to the new patient.

## NIBP Trends

A maximum of 128 NIBP readings are collected (up to 8 hours). NIBP is not measured continuously like other vital sign parameters. The numerics on this trend are captured at the time of the NIBP reading.

NIBP TREND							( 85 )	58
TIME	HR/PR	SpO2	SYS/DIA-MEAN		BR	PA	35/ 18	
HH:MM	BPM	%	mmHg		Br/M	MC02	BR	
10:01	101	98	139/ 74	106	11	35	12	
9:46	88	SRCH	142/ 83	110	10	mmHg		
9:31	93	OFF	145/ 87	112	5	SpO2	92	
9:20	91	100	126/ 85	91	11			
PRINT		↑	↓	↔	NXT TRND	OXYCRG		

When the SEARCH message appears in an NIBP TREND display or printout, it indicates that the monitor was not able to complete an NIBP measurement during that time period.

This symbol indicates that the NIBP reading was taken in the presence of high motion artifact while monitoring ECG. Artifact can affect accuracy. To help reduce artifact, see page 59.

If an error number (ERR# x) is listed in an NIBP TREND printout or display, it indicates that an NIBP equipment alert occurred. See page 62 for NIBP alert error numbers and definitions.

## Displaying Trends

Displayed trends show the last five hours of data. Trends are displayed if you are at the main menu and have only one waveform turned on in the wave selection window. Trends are also displayed when you press the **TRENDS** button in the first Setup Menu.



## Selecting a Trend (NXT TRND)

Select the trend you want displayed by pressing the **NXT TRND** button. The trend is identified by a label at the top of the table.

P1 TREND	P2 TREND	NIBP TREND	RESP TREND	TEMP TREND
Time	Time	Time	Time	Time
HR/PR	HR/PR	HR/PR	HR/PR	HR/PR
SpO <sub>2</sub>	SpO <sub>2</sub>	SpO <sub>2</sub>	SpO <sub>2</sub>	SpO <sub>2</sub>
SYS/DIA-Mean	SYS/DIA-Mean	SYS/DIA-Mean	RR or BR	T1
RR or BR	RR or BR	RR or BR	ETCO <sub>2</sub>	T2
			INCO <sub>2</sub>	ΔT



# Printing

---

## *Printing Patient Data*

### **Printing Waveforms**

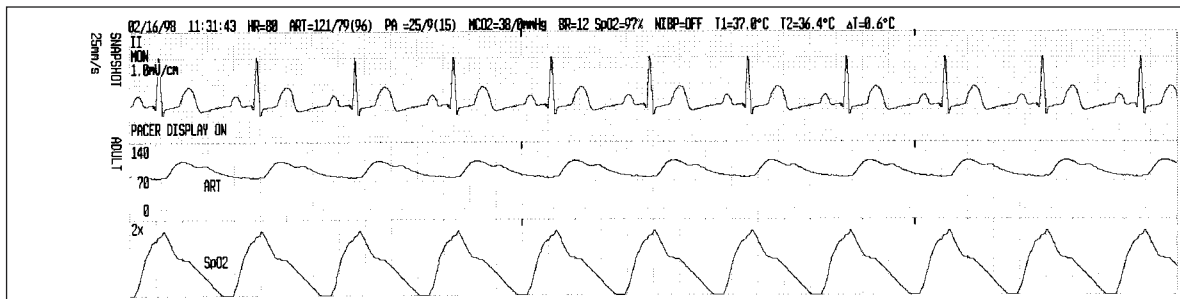
If you have an expansion module set up with your monitor, you can print any waveform by pressing either the **SNAPSHOT** or **START/STOP** button on the expansion module when the waveform is displayed on the monitor. You can also set the monitor to print automatically.

The number of seconds of data shown on the printout depends on the print speed set in the Printer Setup window.

Patient vital sign numerics are printed above the waveforms.

### *The ECG Waveform*

ECG is printed on a grid with major divisions (dotted lines) every 5 mm and minor divisions (single dots) every 1 mm. The ECG waveform is always printed if ECG is monitored.



### *The Invasive-pressure Waveforms*

Pressure waveforms are printed on a grid with major divisions (vertical dotted lines) every 5 mm, and the pressure scale grids are printed horizontally.

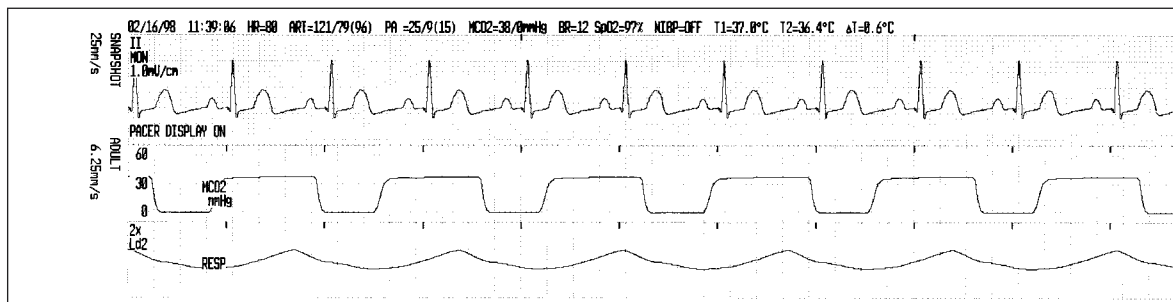
### *The SpO<sub>2</sub> Plethysmograph*

The plethysmograph is printed without a horizontal grid. The size is printed on the printout. Although no vertical scale is displayed for the SpO<sub>2</sub> waveform, a size indicator is displayed to show the relative gain of the waveform.

## CO<sub>2</sub> and RESP Waveforms

CO<sub>2</sub> waveforms are printed on a grid with major divisions (vertical dotted lines) every 5 mm, and the pressure scale grids are printed horizontally. The waveforms and numerics are labeled as **MCO<sub>2</sub>** for Mainstream CO<sub>2</sub> and **SCO<sub>2</sub>** for Sidestream CO<sub>2</sub>. If neither CO<sub>2</sub> option is active, the label is simply **CO<sub>2</sub>**.

On the **SNAPSHOT** command only, the CO<sub>2</sub> and RESP waveforms will be printed with different sweep speeds than the other waveforms.



## Printing NIBP Measurements

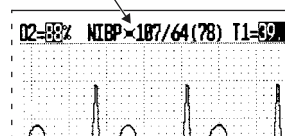
You can print the results of an NIBP measurement each time one occurs. This printout is called the NIBP Ticket. The NIBP Ticket must be turned on in the Printer Setup window.

1. To turn on NIBP TICKET, press **SETUP, MORE, PRINTER**. The printer setup window appears.
2. Use the **NEXT** and **CHANGE** buttons to select and turn on the NIBP TICKET.

NIBP TICKET			02/16/98				
TIME	HR/PR	SpO2	SYS	/ DIA	- MEAN	BR	
HH:MM	BPM	%	mmHg		Br/M		
11:46	80	97	120	/ 63	80	12	
ADULT							

This symbol indicates that the NIBP reading was taken in the presence of high motion artifact while monitoring ECG. Artifact can affect accuracy. To help reduce artifact, see page 59.

NIBP TICKET			02/16/98		
TIME	HR/PR	SpO2	SYS / DIA	- MEAN	BR
HH:MM	BPM	%	mmHg		Br/M
11:46	80	97	120 / 63	80	12



## Printing the Apnea Ticket

The Apnea Ticket documents the length of each apnea episode.

The Apnea Ticket is printed  
after the patient resumes  
breathing . . .

. . . and/or 60 seconds after the last  
breath was detected if the patient has  
not resumed breathing

APNEA TICKET (BR) 02/16/98		
TIME	HR/PR	SpO2
H:MIN:S	BPM	%
LAST BREATH:		
15:45:53	80	OFF
RESUMED BREATHING:		
15:47:02	80	OFF
ELAPSED TIME:		
00:01:09		
ADULT		

APNEA TICKET (BR) 02/16/98		
TIME	HR/PR	SpO2
H:MIN:S	BPM	%
LAST BREATH:		
15:45:53	80	OFF
BREATHING NOT RESUMED:		
15:46:53	80	OFF
ELAPSED TIME:		
00:01:00		
ADULT		

To set the Propaq to print an Apnea Ticket after an apnea event, follow these steps.

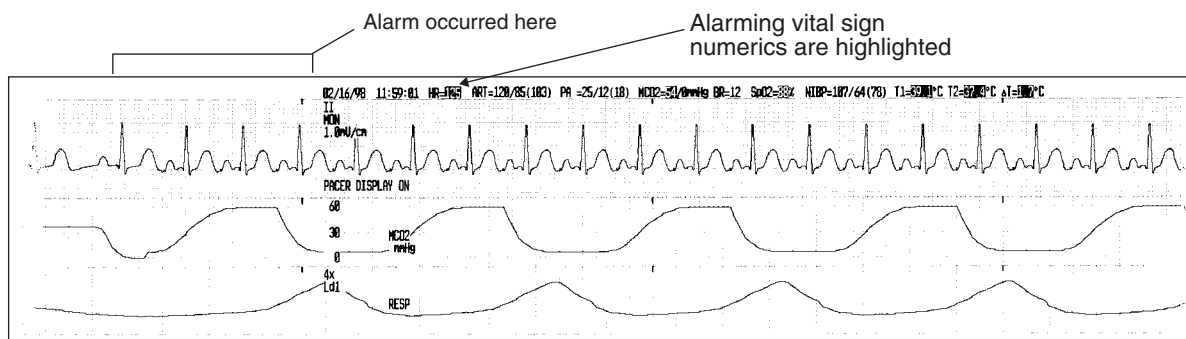
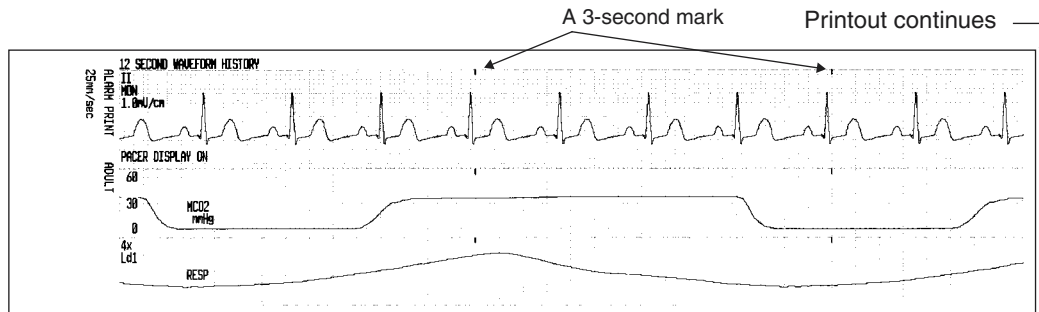
1. From the main menu, press the following buttons: **SETUP**, **MORE**, **PRINTER**.
2. Press the **NEXT** button until APNEA TICKET is highlighted in the Printer Setup window.
3. Press the **CHANGE** button until APNEA TICKET is set to ON.

## Printing When a Patient Alarm Occurs

To set up the printer to print on a patient alarm, follow these steps.


1. From the Main Menu, press the following buttons: **SETUP, MORE, PRINTER.**
2. Press the **NEXT** button to select ALARM PRINT in the printer setup window.
3. Press the **CHANGE** button until ALARM PRINT is set to ON.

PRINTER	SETUP PAGE	( 85 ) 58
CONTINUOUS	: 25.0 mm/s	PA 35/ 18
AUTO PRINT	: 15 minutes	MC02 35 BR 12
ALARM PRINT	: OFF	mmHg
NIBP TICKET	: OFF	SpO2 92
APNEA TICKET	: OFF	
OXYCRG ON ALARM	: OFF	
NEXT	CHANGE	PR TREND MORE PREV MENU



OxyCRG

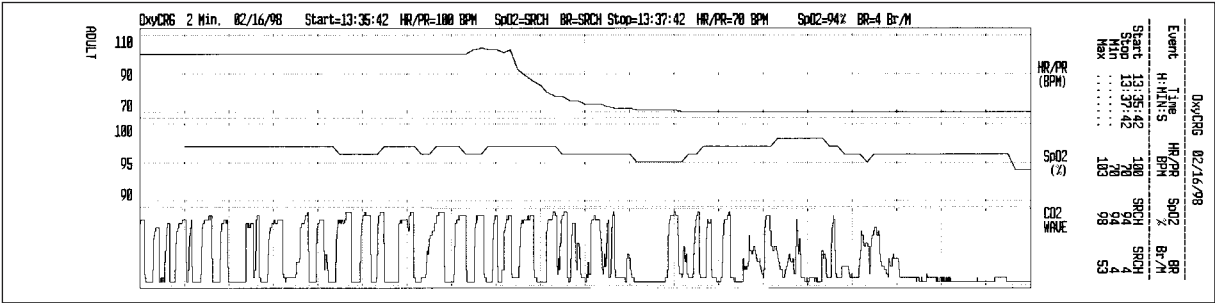The OxyCRG is a graphical printout of two minutes of continuous HR/PR and SpO<sub>2</sub> numerics, and a condensed respiratory waveform. If any of the parameters have been completely inactive for the two minutes prior to the initiation of the print, the associated band will be empty.



OxyCRG on Alarm

When an alarm condition is detected, a print will be initiated if OXYCRG ON ALARM is turned on. If an SpO<sub>2</sub> or HR/PR alarm condition is detected, an OxyCRG will be queued to print 60 seconds after the alarm is detected. If an Apnea or RR/BR alarm condition is detected, an OxyCRG will be queued to print 75 seconds after the alarm is detected.

The parameters which have alarmed in the two-minute period are indicated by the highlight of the corresponding labels.



## Printing Trends

Printed trends are useful for reviewing the patient's vital signs over the last several minutes to the last five hours. The Propaq enables you to print one trend or several trends with a press of a button, or automatically at 4-hour intervals.

### Printing a Single Trend

The best way to print just one trend is with the **PRINT** button in the Trends Menu. When you press the **PRINT** button, the displayed trend is printed. If you want to print a trend different from the one displayed, press **NXT TRND** until the desired trend is shown.

### Printing Several Trends

The best way to print several trends at one time is to set up the printer to print the trends you want and then press the **PRINT TRENDS** button on the expansion module, or the **PR TREND** button in the Printer Menu, whenever you want the trends printed.

1. From the Main Menu, press **SETUP, MORE, PRINTER, MORE**. The printer trend select window appears.
2. Using the **NEXT** and **CHANGE** buttons, select each of the trends you want printed and turn them on. Turn off all other trends.
3. Now, each time you want to print the selected trends, press **PRINT TRENDS**.

NIBP TREND		02/16/98					
TIME	HR/PR	SpO2	SYS	/ DIA	- MEAN	BR	
HH:MM	BPM	%	mmHg			Br/M	
10:46	80	97	114	/ 67	81	12	
10:45	80	97	116	/ 67	82	12	
10:44	80	97	118	/ 68	81	12	
ADULT							

This symbol indicates that the NIBP reading was taken in the presence of high motion artifact while monitoring ECG. Artifact can affect accuracy. To help reduce artifact, see page 59.

## Automatic Trend Prints

To automatically print trends at 4-hour intervals, activate AUTO TREND and select the print times.

Use the **CHANGE** button to set the trend print times according to the start time (clock hour) of each shift or to OFF.

Once the Propaq is set up, it will print all the selected trends at each 4-hour interval. You can place the printouts in the patient's record at the end of each work shift.

Starting times are selectable for the 4-hour interval printouts

PRINTER	TREND SELECT PAGE	( 85 )	58
AUTO TREND: 03 07 11 15 19 23	PA	35/ 18	
NIBP : OFF	MC02	BR	
RESP : ON	35	12	
ART : OFF	mmHg		
PA : ON	SpO2	92	
TEMP : OFF			
NEXT	CHANGE	PR TREND	MORE
		PREV MENU	

# Acuity Central Monitoring System

---

## ***Intended Use***



**Warning** Use of equipment, accessories, and parts not recommended or supplied by Welch Allyn could result in inaccurate patient information or damage to the system.

When Acuity is inoperable, be sure to keep Acuity patients under close surveillance, especially those prone to arrhythmias. Use Acuity only in conjunction with close surveillance by trained clinicians.

Connect the Propaq Encore to an Acuity system only. Connecting to other networks could damage the monitor or injure the patient. If in doubt about the network jacks or devices, consult your facility's Biomedical Engineering Department.

Make sure the Acuity network cable is not damaged. The Acuity network cable is the sole link between the Propaq Encore and the Acuity Central Monitor.

The Acuity system is used as a central monitoring system for Propaq Encore monitors configured to interface with Acuity.

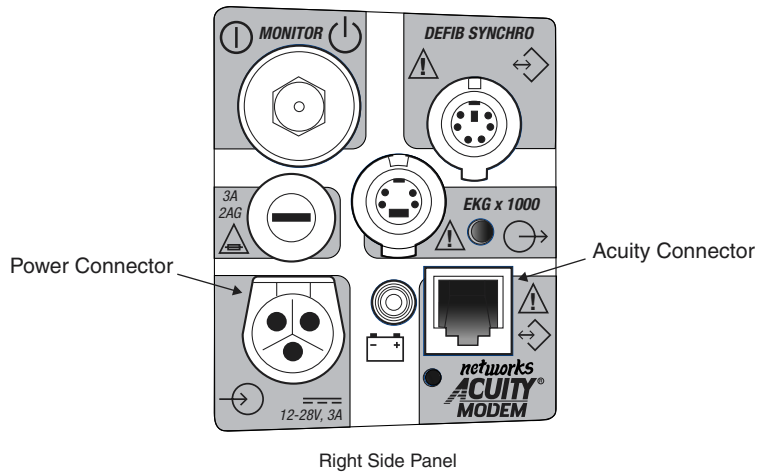
The Acuity system may be used to monitor all patients. For neonatal patients, use all Acuity features except the Welch Allyn Cordless Acuity and the arrhythmia detection option. The Acuity system is intended to be used in compliance with the instructions in this *Propaq Encore Reference Guide*, the *Acuity System Reference Guide*, and accepted hospital and clinical protocols.

## ***Connecting to the Acuity System***

1. If the Propaq Encore has already been connected to the patient, save the patient's Trends and Alarm Limit settings by keeping the monitor turned on. (The Propaq Encore transmits up to five hours of trend information when you connect it to the Acuity network.)

If the Propaq Encore has not been connected to the patient, clear any prior patient's trends and alarm limit settings by turning off the Propaq Encore and after a few seconds, turning it on again.

2. If the Propaq Encore is not already connected to the patient, attach leads and sensors to the patient as described in Chapter 2 of this reference guide.
3. Plug in the Acuity network cable to the Acuity network jack on the Propaq Encore side panel shown in the following figure. Plug in the other end of the cable to the bedside Acuity network jack.



**Caution** Accessory equipment connected to the analog and digital interfaces must be certified according to the respective IEC standards (e.g., EN 60950 for data processing equipment and EN 60601-1 for medical equipment). Furthermore, all configurations shall comply with the system standard IEC 601-1-1. Anyone connecting additional equipment to the signal input or output connectors is configuring a medical system, and is therefore responsible that the system complies with the requirements of the system standard IEC 601-1-1. If in doubt, consult your Biomedical Engineering Department.

4. Connect the AC adapter to the Propaq Encore and the wall outlet to charge the battery. Check to see that the battery charging light on the monitor's right side panel is on.
5. Confirm the patient identification at the bedside or enter the patient information at the Acuity Central Monitor using the Patient ID Setup Window.
6. If alarm limits have not been set, do so at the Propaq Encore or at the Acuity Central Monitor using the Alarms Setup Window.

## Press NET OFF to Disconnect from Acuity



**Warning** If you don't set alarm limits, Acuity uses preset settings (for arrhythmia limits), and the powerup default settings for the Propaq Encore.



When a Propaq Encore in Adult or Pediatric Mode is connected to an Acuity System, the audio alarms at the bedside Propaq can be delayed up to 4 minutes and 15 seconds. The delay time is selected in Acuity software at the time of Acuity installation. Visual alarm indications are not delayed.

## Key-press Route to Acuity Menu

### MAIN MENU

NIBP	ECG/RESP	INV PRS	SpO2/CO2	SETUP
------	----------	---------	----------	-------

### SETUP MENU 1

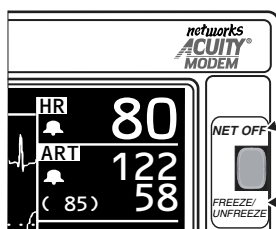
STATSCALE	ALARMS	WAVE SEL	TRENDS	MORE
-----------	--------	----------	--------	------

### ACUITY MENU

		NET OFF	SNAPSHOT	PREV MENU
			FRZ PRNT	

## Press NET OFF to Disconnect from Acuity

When you want to disconnect the monitor from Acuity, be sure to use the **NET OFF** button as described below. The **NET OFF** button is a safety feature to help make sure the patient is not disconnected accidentally. If you simply disconnect the Acuity network cable from the monitor without using the **NET OFF** button, the monitor and Acuity both generate equipment alert messages that must be acknowledged by an operator.



### NET OFF

When the monitor is connected to Acuity, pressing this button disconnects the monitor from the Acuity network.

### FREEZE/UNFREEZE

When the monitor is NOT connected to Acuity, pressing this button freezes or unfreezes the displayed waveforms.

1. To disconnect the Propaq Encore from the Acuity network, press the front panel **NET OFF** button.

2. Within 15 seconds, disconnect the Acuity network cable from either the Propaq Encore side panel or the bedside jack. If the patient will no longer be monitored with the Propaq Encore, turn off the monitor to erase trend information.

If you do not disconnect the Acuity network cable within about 30 seconds, the monitor and Acuity attempt to reconnect and prompt you to confirm the patient ID.

Whenever the monitor is connected to Acuity, the **NET OFF** function overrides the **FREEZE/UNFREEZE** function for the **NET OFF** front panel button. Therefore, you cannot freeze the displayed waveforms when the monitor is connected to Acuity. Also, you cannot initiate a Freeze Print at Acuity when connected to Acuity. When the monitor is not connected to Acuity, pressing the **FREEZE/UNFREEZE** button freezes or unfreezes the displayed waveforms.

## ***Printing at Acuity***

You can print various waveforms from the Propaq Encore to the Acuity system printer. To print a waveform shown on the display screen, press **SETUP, ACUITY, SNAPSHOT**. If you press the **FREEZE** button on the front of the Propaq Encore, the button changes from **SNAPSHOT** to **FRZ PRNT**.

## ***Network Alert Message***

When the Propaq Encore is connected to the Acuity system, it constantly exchanges information with Acuity. If the Propaq Encore detects an interruption in this flow of information, it displays an alert message: **NETWORK FAULT, CHECK ACUITY/DATA COMM CONNECTION**.

If the Propaq Encore displays this message, check the Acuity network cable to be sure it is plugged in to the side panel and to the bedside jack. If the cable is damaged, replace the cable.

If the cable appears undamaged and the Acuity system is operating normally, ask your service personnel to check the network and the Propaq Encore Acuity connector.

# Power Sources

---

For in-hospital operation and recharging from ac mains, an ac power adapter plugs into the monitor. Use only a Welch Allyn ac power adapter to ensure protection against risk (leakage) current hazards.

The Propaq Encore can also be powered and recharged from a dc source (isolated from ac mains) capable of supplying 12-28 Vdc and continuously supplying 25W.



**Caution** When a transport vehicle's battery system is used to provide input power to the Propaq Encore, surges caused by a defect in the vehicle's power system may blow a fuse in the Propaq's side panel or cause further damage to the Propaq.

## ***Power Adapter Intended Use***

Welch Allyn power adapters are intended to be used only with Propaq monitors, and Propaq Encore monitors are intended to be recharged using only a Welch Allyn power adapter with a mating plug, and rated for your ac mains.

The power adapter contains symbols on its labeling. For definitions of these symbols, see page 11.



**Warning** Place the power adapter where it cannot fall and harm someone.



**Caution** Use of other than Welch Allyn power adapters with the plug rated for your ac mains can damage or compromise the safety of the Propaq Encore monitor and may require fuse replacement in the power adapter. Verify that the Power Adapter is set for the proper mains voltage prior to plugging it into the Propaq.

Do not autoclave the power adapter. Do not operate the power adapter with a damaged case, mains power cord, or plug.

## Verifying Proper Power Adapter Configuration

Prior to using the power adapter, check it for proper voltage selection by looking in the small indicator window on the front end (by the power switch). If the number in the window does not match your ac mains source voltage (100-120V or 200-240V), the adapter should be reconfigured. See the illustration on page 121.

Your biomedical technician can change the voltage setting and fuses on the Power Adapter and can verify that your facility is using the correct power cord.

### ***Power Adapter Configurations***

Power adapter part numbers are listed with their rated input, fuse, output, and applicable usage. Check to be sure you are using the correct power adapter for your mains power source by comparing the part number on the power adapter to the table below. Always replace fuses with the fuses rated for the power adapter.

Part Number/ Connector Style	Rated Input Voltage	Rated Serviceable Fuses	Rated Output	Application
With power switch: 503-0054-00 Without power switch: 503-0093-XX	100V-120V ac, 500 mA, 50/60 Hz	T800 mA/250V Time-Delay 5 x 20 mm	16-24V dc 25 VA	25 Watt requirement in countries with 100V-120V power systems.
With power switch: 503-0054-01 Without power switch: 503-0092-XX	200V-240V ac 250 mA, 50/60 Hz	T400 mA/250V Time-Delay 5 x 20 mm	16-24V dc 25 VA	25 Watt requirement in countries with 200V-240V power systems.

## Replacing the Power Adapter Fuses

The Power Adapter contains two fuses that can easily be replaced by service personnel if necessary. The adapter can contain spare fuses.



**Caution** Replace each fuse only with the specified type.



## Procedure

1. Unplug the power adapter's removable cord from the ac mains outlet and the power adapter.
2. Turn the power adapter so you can see the window that indicates the voltage setting.
3. Using a small, flat-blade screwdriver, carefully pry the fuse module from the power adapter.
4. Replace both fuses with the specified type.

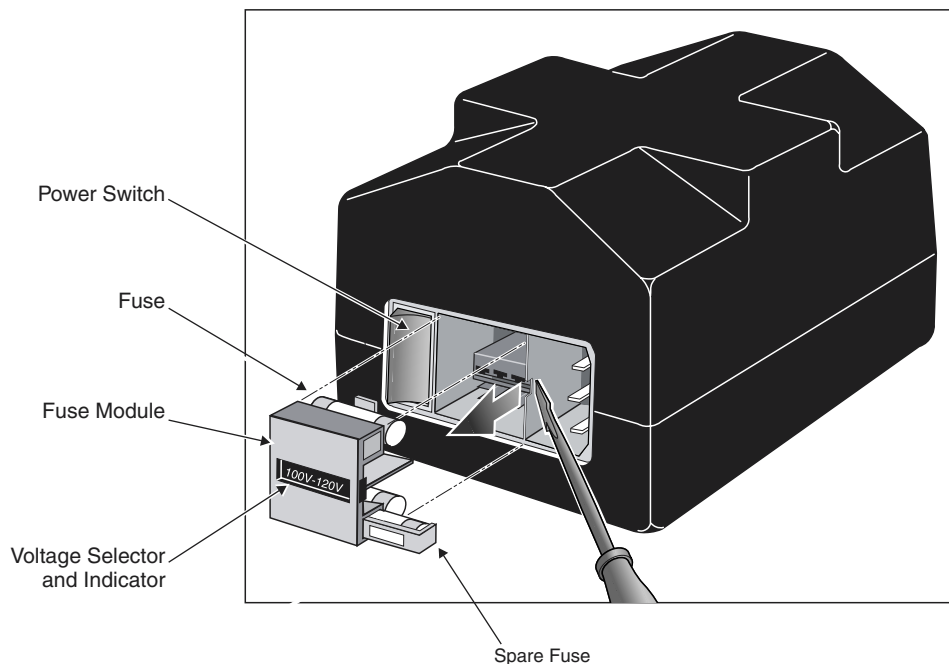


Both fuses should be replaced at the same time, even if only one fuse has opened due to an overcurrent situation. The unopened fuse may have been stressed and could become unreliable.



**Caution** Spare fuses are contained in housings next to the fuses in the fuse module as shown in the following picture. Between the fuses is a small printed-circuit board (PCB) that sets the power adapter to the desired ac mains voltage. When handling the fuse module, the PCB may slide out.

If the small PCB between the fuses has slipped out of place, slide it back into place in the fuse module, and verify that the voltage setting indicated in the window on the fuse module is correct. If the voltage setting is incorrect, simply slide the PCB out of the fuse module, rotate it 180° and slide it back into place.





**Caution** Make sure the voltage selector indicates the proper ac input voltage. If you change the adapter voltage setting, you must replace all fuses to match the appropriate type specified on the bottom of the power adapter. The only fuses contained in the power adapter when shipped from the factory are fuses specified for the original adapter input voltage setting.

## ***Battery Care***



**Caution** Leaving the monitor's lead-acid batteries in a completely discharged state may result in permanent battery damage. The batteries should be kept fully charged.

## **Recharging Time**

The battery charges to full capacity within eight hours (if the monitor remains off).

## **Monitor Functions Resumed**

In most instances, most monitor functions are usable immediately after plugging in the ac power adapter and cycling the power switch. More charging time may be required before the NIBP, CO<sub>2</sub>, and printer can be operated.

## **Operating Times Using Battery Power**

The amount of time you can operate the Propaq Encore on each battery charge depends upon many variables including active options, frequency of NIBP measurements, frequency and length of print strips, ambient temperature, battery age and condition, and what information is displayed.

For monitors without the Expansion Module or SpO<sub>2</sub> option, typical monitor operating time is about 2 hours at 25° C for a new, fully-charged battery. This is when all patient channels are active and measurements are taken every 15 minutes.

For monitors without the Expansion Module but with the SpO<sub>2</sub> option, operating time is about 5 hours.

For monitors with the Expansion Module and printer, SpO<sub>2</sub> and CO<sub>2</sub> options, typical operating time is about 3 hours under the above conditions when print strips are generated every 15 minutes.

## Monitor Functions Based on Battery Voltage

As battery voltage drops during extended monitor battery operation, error messages are displayed and monitor functions are discontinued in order of priority.

## Checking Battery Voltage

The Propaq's battery voltage is displayed on the initial powerup screen. The battery voltage is also displayed with the Settings window.

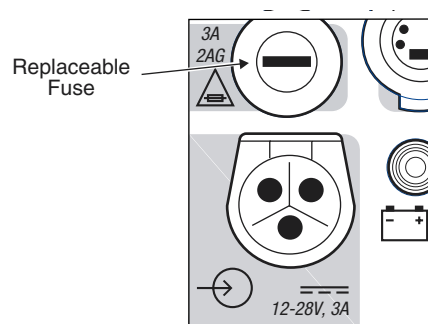
## Replacing the Fuse

The Propaq Encore is protected against power surges by a 3 Ampere fuse, which can easily be replaced in the right side panel. Fuse replacement should only be performed by a qualified service person.



If the green Battery Charging lamp does not light when an ac adapter is connected, this fuse may be blown.

1. Disconnect the Propaq Encore from the patient.
2. Disconnect the ac power adapter from the Propaq Encore's power input connector.
3. Using a small screwdriver or similar device, unscrew the fuse carrier by turning it counterclockwise.
4. Remove the fuse holder and replace the fuse with a 3 Ampere, 2AG, 250V (fast or slow acting) fuse. This fuse can be ordered from Welch Allyn or its service centers using part number 503-0058-00.



# Care and Maintenance

## *Inspect and Clean the Monitor and Accessories*

Before cleaning, thoroughly inspect the monitor and all accessories for any signs of damage, cracks, or improper mechanical function of keypads, switches, connectors, and printer paper door. While gently bending and flexing cables and tubing, inspect for damage, cracks, cuts, abrasions, extreme wear, exposed wires or bent connectors. Confirm connectors securely engage. Report damage or improper function to your service department.

Equipment	Cleaning Instructions	Approved Cleaning Solutions	
Propaq CS Monitor <sup>a</sup> <sup>b</sup>	<ul style="list-style-type: none"> <li>• Wipe with a nearly-dry cloth moistened with cleaning solution.</li> <li>• Thoroughly wipe off any excess cleaning solution. Do not let water or cleaning solution run into connector openings or crevices.<sup>c</sup></li> </ul>	Warm water Liquid soap T.B.Q. <sup>®b</sup> Hydrogen peroxide solution	Coverage <sup>®</sup> Wex-cide <sup>®b</sup> Windex <sup>®</sup> Ovation <sup>®</sup>
NIBP cuff	<ul style="list-style-type: none"> <li>• Wipe gently with cloth dampened with cleaning solution.</li> <li>• Thoroughly wipe off excess cleaning solution. To avoid harming cuff function, do not let water or cleaning solution enter cuff tubing.</li> </ul>	Common hospital disinfectants, including Clorox <sup>®</sup> liquid bleach (1:10 solution of Clorox/water), isopropyl alcohol, Lysol <sup>®</sup> solution, PhisoHex <sup>®</sup> , Quatricide <sup>®</sup> , Virex <sup>®</sup> and Vesphene <sup>®</sup>	
Cables, tubing, CO <sub>2</sub> sensor <sup>d</sup>	<ul style="list-style-type: none"> <li>• Wipe gently with cloth dampened with cleaning solution. Do not immerse the CO<sub>2</sub> sensor in liquid.</li> </ul>	Mild detergent solution; also consult manufacturer's instructions.	
Masimo SpO <sub>2</sub> cables Nellcor SpO <sub>2</sub> cables	<ul style="list-style-type: none"> <li>• Wipe gently with cloth dampened with isopropyl alcohol.</li> </ul>	Isopropyl alcohol	
Other accessories	<ul style="list-style-type: none"> <li>• Consult manufacturer's instructions.</li> </ul>	Consult manufacturer's instructions.	

a. **Do not use these cleaning solutions (they may damage the monitor):** Butyl alcohol, Denatured ethanol, Freon<sup>™</sup>, Mild chlorine bleach solution, Isopropyl alcohol, Trichloroethane<sup>®</sup>, Trichloroethylene, Acetone, Vesphene II, Enviroquat<sup>®</sup>, Staphene<sup>®</sup>, Misty<sup>®</sup>, Glutaraldehyde, Fantastik<sup>®</sup>, Formula 409<sup>®</sup>, Cidex<sup>®</sup>.

b. The monitor may be disinfected to comply with OSHA requirements for cleaning and decontaminating spills of blood and other body fluids. (Federal OSHA Standard on bloodborne pathogens: 29 CFR 1910.1030, 12/6/91.) **Wex-cide** (Wexford Labs, Inc., Kirkwood, MO) and **T.B.Q.** (Calgon Vestal Lab., Calgon Corp., St. Louis, MO) are disinfectants that meet OSHA requirements, and are EPA approved. Wipe away disinfectants with a water-dampened cloth after the manufacturer's recommended period of time.

c. If liquid gets into the right side panel connectors, it will drain out. If moisture gets into a left side panel connector, dry the connector with warm air, then check the monitoring functions for proper operation.

d. The Mainstream CO<sub>2</sub> sensor may also be disinfected with Wex-cide. Follow the disinfectant manufacturer's instructions. Do not leave Wex-cide on sensor longer than 30 minutes. Thoroughly clean off residue with water-dampened cloth. Prolonged exposure of the sensor to Wex-cide will damage the sensor.



**Caution** Do not autoclave this product or its accessories. Do not immerse the monitor in liquid when cleaning. Do not immerse accessories in liquid when cleaning unless the accessory manufacturer's cleaning instructions instruct you to do so.

## Maintenance

### Service Interval Recommendations

At the intervals recommended below, verify the Propaq Encore for proper operation of all channels and internal circuitry. Such checks and verifications should only be carried out by a qualified biomedical service person.

Other Propaq Encore service information, including calibration procedures, is described in the *Propaq Encore Service Manual* (P/N 810-0696-XX). Refer to it for more information.

Use the following intervals for a guideline. Service may be needed more often in extreme environments (heat, cold, dust, etc.).

Recommended Interval	Service Action
Six months to two years	<ul style="list-style-type: none"> <li>• Complete functional verification; see <i>Propaq Encore Service Manual</i></li> <li>• Inspect the Propaq Encore for mechanical and functional damage</li> <li>• Inspect safety labels for legibility</li> <li>• Inspect the side panel fuse for compliance to specified rating</li> <li>• Verify that visual and acoustic alarms are functioning properly</li> <li>• Test patient leakage current according to IEC 601-1/1988</li> <li>• Test patient leakage current with mains voltage on patient-applied parts according to IEC 601-1/1988: limit 50<math>\mu</math>A<sup>a</sup></li> </ul>
Minimum every three years	<ul style="list-style-type: none"> <li>• Check battery capacity</li> </ul>

a.NOTE: The leakage current should never exceed the 50 $\mu$ A limit. The data should be recorded in an equipment log. If the device is not functioning properly or fails any of the above tests, do not attempt to repair the device. Please return the device to the manufacturer or to your distributor for any required repairs.

### Product Recycling

You can return a Propaq Encore to Welch Allyn for recycling when the monitor reaches the end of its life (7 to 10 years). You can also return an Encore's internal battery to Welch Allyn for recycling when it reaches the end of its life (3-4 years).

## Monitor Care

### Environmental Operating and Storage Limits

Whenever possible, store the Propaq Encore at room temperature in a dry environment. For environmental operating instructions, see page 153.



**Warning** The monitor may not meet its performance specifications if stored or used outside the specified temperature and humidity ranges.

### Extended Storage Precautions

#### Battery Removal



**Caution** Storing the Propaq Encore for extended periods (more than three months) without being connected to the ac power adapter can cause damage to the battery. Even when the Propaq Encore is turned off, a very small amount of current is drawn from the battery. For long-term storage, remove the battery from the Propaq.

See the *Propaq Encore Service Manual* for procedures on removing the battery.



Removing the battery will cause programmed settings to be lost, but they can be reprogrammed when the battery is replaced.

## ***Printer Paper Removal***



**Caution** If a Propaq Encore has a battery installed or ac power connected and is stored for an extended period without use, the printer paper can cause damage to the printhead. Before storing a Propaq Encore for more than two months without use, remove the roll of printer paper.

## ***Printer Maintenance***

### **Loading Paper**

Paper is loaded through the bottom of the printer.

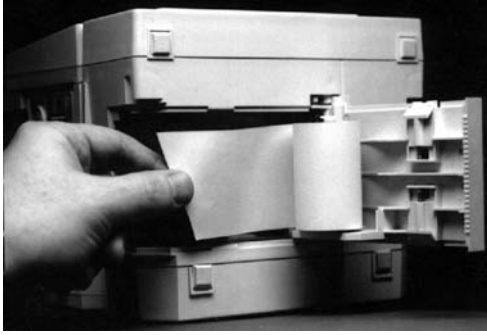


**Caution** Use only low-debris printer paper purchased from or recommended by Welch Allyn. Use of other paper can cause unclear printing of patient data, damage to printing head, and eventual printer failure. Store all paper (including a monitor loaded with paper) in an environment that meets the paper storage specifications listed in Appendix B. Failure to properly store paper can result in paper discoloration and damage to the printer.


1. Lay the monitor on its back to gain access to the bottom of the printer.
2. Squeeze the locks on the paper door toward each other and pull the door toward you to open it.
3. Lift the paper roll from the holder and pull out any paper remaining in the printing mechanism.
4. Place the new paper roll onto the holder, as shown below, and pull out several inches of paper.



## Printer Maintenance




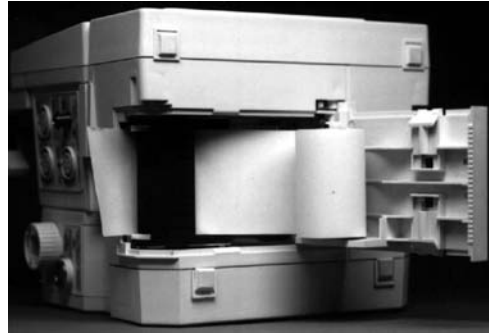Load the new paper roll onto the spindle on the door.



Feed the paper through the printer mechanism.

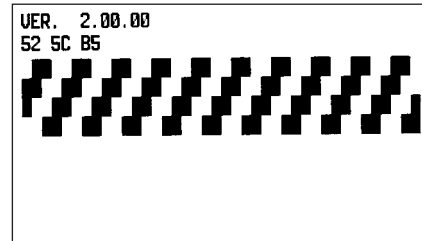5. Slide the end of the paper into the slot of the printing mechanism until it extends out of the paper exit slot.
6. Close the paper door.
7. Place the monitor on its feet.
8. Simultaneously press the **START/STOP** button and the **PRINT TRENDS** button to produce a test print.



## ***Customer Services***

### **Ordering and Customer Service**

For ordering information, for the location of your nearest Welch Allyn sales representative or service center, or for more information on other Welch Allyn products, contact:

Welch Allyn  
8500 SW Creekside Place  
Beaverton, OR 97008-7107 USA  
Worldwide: (503) 530-7500  
In the USA, toll-free: (800) 289-2500  
FAX: (503) 526-4200

### **Technical Service**

If you need technical assistance on troubleshooting, are interested in customer technical training on Welch Allyn products, or help with ordering replacement parts, contact Welch Allyn's Technical Services Department at:

Worldwide: (503) 530-7500  
In the USA, toll-free: (800) 535-6663  
FAX: (503) 526-4970  
Internet: <http://www.monitoring.welchallyn.com/service>  
Internet E-mail: [solutions@monitoring.welchallyn.com](mailto:solutions@monitoring.welchallyn.com)

### **Repacking**

Before returning the monitor for service, call Welch Allyn for return authorization and instructions on shipping. Repack the Propaq Encore in its original shipping container, if possible. The container is designed to protect the monitor from possible damage during shipment.

# Defibrillator Synchronization

*These instructions are not intended to replace existing hospital procedures relative to the provision of cardiac electrical therapy and operation of the specific models of defibrillators. Use all safety standards and clinical protocols as defined by your institution.*

The Defibrillator Synchronization feature is designed to operate only with the Physio-Control LIFEPAK 5® and LIFEPAK 6s® defibrillators.

## ***Defib Sync Connector and Cables***

The DEFIB SYNCHRO connector is on the Propaq's right side panel.



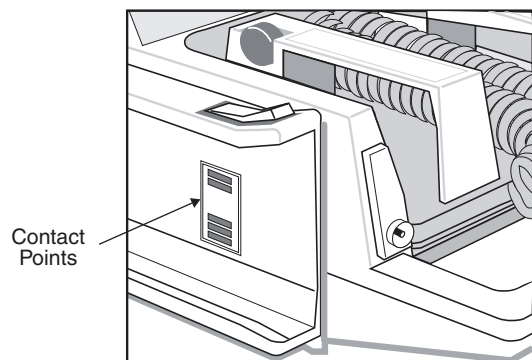
**Warning** Use only the correct Welch Allyn cable with the LIFEPAK 5 or LIFEPAK 6s. (These cables contain circuitry in addition to their wiring.) The use of any other cable will result in incorrect operation. The Welch Allyn *Products and Accessories Booklet* lists correct cables for use with Welch Allyn products.

Defibrillator Application	Welch Allyn Part Number
LIFEPAK 6s	008-0154-00
LIFEPAK 5	008-0136-XX


## ***Installing the Defib Sync Interface to the LIFEPAK 5***

The Propaq Encore Defib Sync interface provides signal transmission between the Propaq Encore and the LIFEPAK 5. The interface connects to the LIFEPAK 5 along the left side of the defibrillator. Use the following instructions to install and remove the interface.


1. Before installing the interface, check that the contacts on the left side of the LIFEPAK 5 are clean in order to ensure signal transmission between the LIFEPAK 5 and the Propaq.



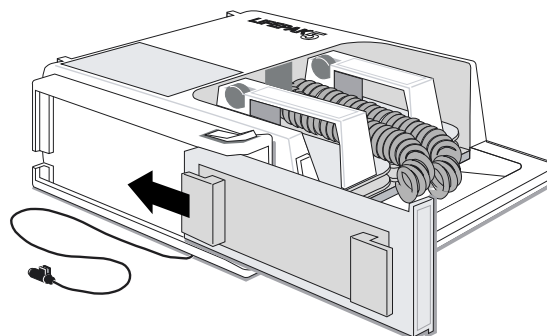2. Slide the Propaq Encore Defib Sync interface onto the left side of the LIFEPAK 5 until it snaps in place.




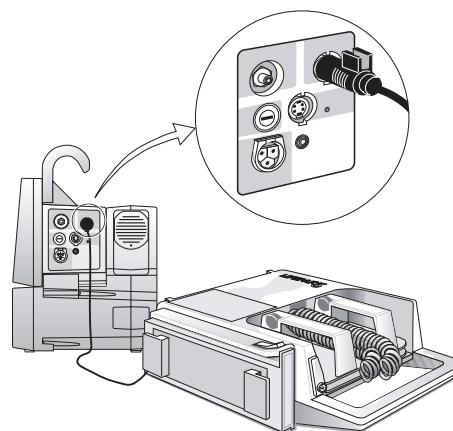3. Connect the cable end to the Propaq's DEFIB SYNCHRO connector on the right side panel.



## ***Removing the Defib Sync Interface from the LIFEPAK 5***

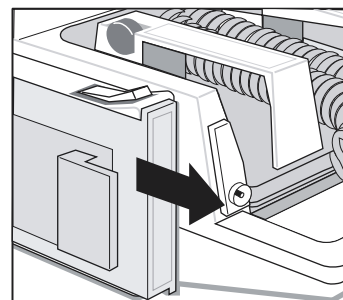1. To remove the interface, disconnect the cable end from the Propaq Encore.
2. Press the lever on the side of the LIFEPAK 5 and slide the interface forward until it is free from the defibrillator.
3. Store the interface in its static-protected plastic bag when not connected to the LIFEPAK 5.



## ***Synchronized Cardioversion Using the LIFEPAK 5***

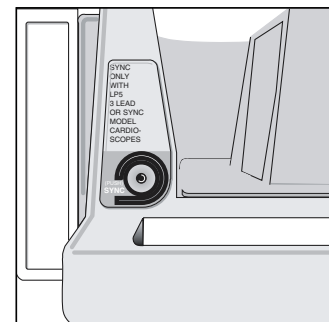
1. Set up the LIFEPAK 5 Defibrillator and any other instrumentation according to institutional procedures and manufacturer's operating instructions. Use the instructions above for installing the Propaq Encore Defib Sync interface.
2. Verify the integrity of the ECG patient electrodes and the fidelity of the ECG waveform on the Propaq Encore for tall, distinct R-waves and minimal artifact.

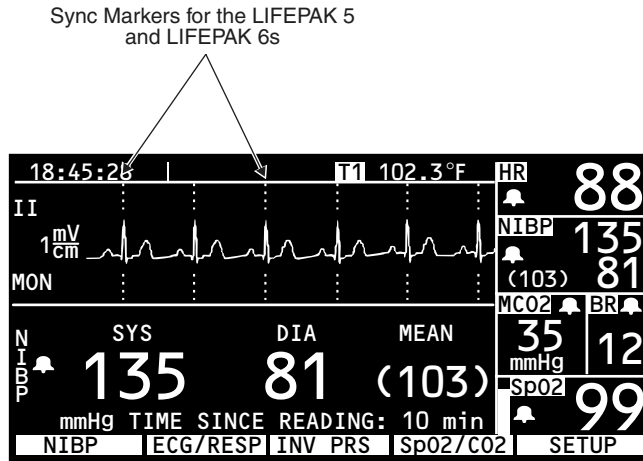


**Warning** The R-wave amplitude must be at least 0.5 mV (5 mm tall when the Propaq Encore ECG sensitivity—**size** button—is set to 1 mV/cm) to guarantee that the defibrillator sync pulse will occur no later than 35 milliseconds after the peak of an R-wave. Reposition the patient electrodes or change the Propaq Encore lead selection as necessary to ensure sufficient ECG waveform amplitude. However, make sure the R-wave amplitude is not so high that it obscures the displayed sync markers.

As a visual gauge for estimating R-wave amplitude, the 'V' of the mV/cm label to the left of the ECG waveform is about 4 mm in height. With the Propaq Encore ECG sensitivity set to 1 mV/cm, compare the letter 'V' with the height of the R-wave, which should be at least 5 mm tall.

3. With the LIFEPAK 5 turned on, press the LIFEPAK's SYNC button on the front-left of the LIFEPAK 5. The button lights when activated. Make sure it lights.
4. After the SYNC button is activated, check the Propaq Encore display for dashed lines above and below each R-wave, occurring in near-synchronization with each R-wave. These are synchronization markers. Check that the LIFEPAK 5's SYNC button also flashes with each R-wave.





**Warning** If the R-wave synchronization markers do not appear to be nearly simultaneous with the R-wave peaks on the Propaq Encore display or are not present, do not proceed with synchronized cardioversion.

You must press the LIFEPAK 5's SYNC button and check for appropriate synchronization markers on the Propaq Encore before each attempt at cardioversion. Welch Allyn cannot guarantee the delay from the sync marker to the defibrillator discharge.



A fault in the cable between the Propaq Encore and the defibrillator, or unplugging the cable, will prevent showing markers on the Propaq and can prevent setting the defibrillator to synchronized mode.

5. Follow hospital procedures and LIFEPAK 5 instructions for cardioversion.
6. If subsequent cardioversion must be performed, repeat steps 3 through 5.

## ***Synchronized Cardioversion Using the LIFEPAK 6s***

1. Set up the LIFEPAK 6s defibrillator and any other instrumentation according to institutional procedures or manufacturer's operating instructions.



The Physio-Control LP6s Defibrillator Sync Connector/Cover (Physio-Control part number 801297-00) must be installed before you can connect it to the Propaq Encore monitor.

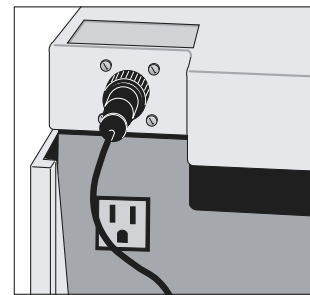
2. Verify the integrity of the ECG patient electrodes and the fidelity of the ECG waveform for tall, distinct R-waves and minimal artifact.



**Warning** The R-wave amplitude must be at least 0.5 mV (5 mm tall when the Propaq Encore ECG sensitivity—**SIZE** button—is set to 1 mV/cm) to guarantee that the defibrillator sync pulse will occur no later than 35 milliseconds after the peak of an R-wave. Reposition the patient electrodes or change the Propaq's lead selection as necessary to ensure sufficient ECG waveform amplitude. However, make sure the R-wave amplitude is not so high that it obscures the displayed sync markers.


As a visual gauge for estimating R-wave amplitude, the 'V' of the mV/cm label to the left of the ECG waveform is about 4 mm in height. With the Propaq Encore ECG sensitivity set to 1 mV/cm, compare the letter 'V' with the height of the R-wave, which should be at least 5 mm tall.

3. Connect the Propaq Encore end of the Welch Allyn Defib Sync cable to the Propaq's right side panel DEFIB SYNCHRO connector.
4. Connect the other end of the cable to the LIFEPAK 6s SYNC connector at the top rear of the defibrillator.
5. With the LIFEPAK 6s turned on, press the SYNC button on the front control panel.



The button lights when activated.

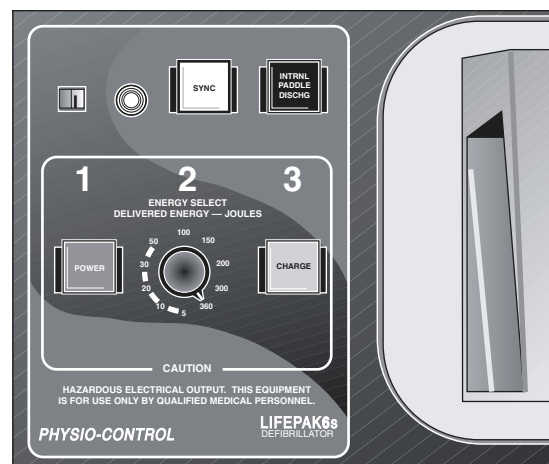6. After the SYNC button has been activated, check that dashed lines appear above and below each QRS on the Propaq Encore display. The LIFEPAK 6s SYNC button will flash with each QRS. (See page 133.)





**Warning** If the R-wave synchronization markers do not appear to be nearly simultaneous with the R-wave peaks on the Propaq Encore display or are not present, do not proceed with synchronized cardioversion.

You must press the LIFEPAK 6s SYNC button and check for appropriate synchronization markers on the Propaq Encore before each cardioversion. Welch Allyn cannot guarantee the delay from the sync marker to the defibrillator discharge.



A fault in the cable between the Propaq Encore and the defibrillator, or unplugging the cable, will prevent showing markers on the Propaq and can prevent setting the defibrillator to synchronized mode.

7. Follow hospital procedures and LIFEPAK 6s instructions for cardioversion.
8. If subsequent cardioversion must be performed, repeat steps 5 through 7.

## ***Defib Sync Message***

The Defib Sync message can appear in an equipment alert window.

**DEFIB FAULT/CHECK INTERFACE CABLE** indicates a probable shorted cable.



# Specifications

## ECG

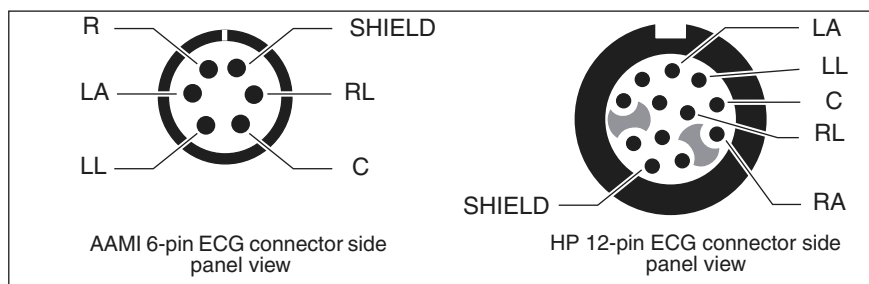
The ECG channel meets all the requirements for Cardiac Monitors Heart Rate Meters and Alarms specified ANSI/AAMI EC13-1992, except for Standardizing Voltage (section 3.2.9.9). The channel also meets the American National Standard, Safe Current Limits for Electromedical Apparatus (ANSI/AAMI ES1-1993).

### ECG Specifications

Characteristic	Specification
Connector	AAMI 6 pin or Hewlett-Packard compatible 12-pin style connector (optional). (See the illustration at the bottom of this table.)
Selectable Leads	I, II, III, aVR, aVL, aVF, V
Lead Fault Indicator	LA, LL, RA, RL, C, multiple
ECG Size (sensitivity) in mV/cm	4, 2, 1, 0.5, 0.2
Display Sweep Speeds	12.5, 25, and 50 mm/sec
QRS Tone Volume	High, Low, Medium, Off
QRS Tone Frequency	900 Hz. for Propaq Encore without Expansion Module, 665 Hertz when equipped with SpO <sub>2</sub> but SpO <sub>2</sub> not being monitored; variable pitch with SpO <sub>2</sub> option and SpO <sub>2</sub> being monitored
Freeze Buffer	3.9 seconds at 25 mm/sec
Bandwidth	0.5 to 40 Hz in Monitor Mode; 0.05 to 40 Hz in Extended Mode (see Real-Time ECG Analog/Defib Sync specification).
Sample Rate	364 Hz
Input Protection	Electrosurgery and defibrillator protected when used with specified ECG cables. All models also include electrosurgery interference suppression.
Lead Fail Sense Current	50 nA dc for active leads 100-200 nA dc for driven lead, depending on number of electrodes attached
Tall T-wave Rejection	Meets and exceeds AAMI (USA) EC13-1992, section 3.1.2.1.c, for 1.2 mV T-wave and 1 mV QRS using AAMI test waveform.
Common Mode Rejection	<1 mV p-p RTI for 10V rms, 50/60 Hz input, 200 pF source impedance, input unbalanced, FILTER function OFF <0.1 mV p-p RTI for 10V rms, 50/60 Hz input, 200 pF source impedance, input unbalanced, FILTER function ON
Input Impedance	>2.5 M $\Omega$ differential @ 60 Hz
Input Range (ac)	10 mV peak to peak
Input Range (dc)	Up to $\pm$ 300 mV
System Noise	$\leq$ 30 $\mu$ V peak-to-peak, R.T.I., with all inputs = 47K in parallel with .047 $\mu$ F.

Characteristic	Specification
QRS Detector	Adult or Pediatric Amplitude Range: 0.22 to 5.0 mV (RTI) Neonatal Amplitude Range: 0.1 to 5.0 mV (RTI) Neonatal and Pediatric Width Range (Duration): 40 to 120 msec Adult Width Range (Duration): 70 to 120 msec
Heart Rate Range	25 to 350 beats per minute (measurement) 25 to 300 beats per minute (display)
Heart Rate Meter Response Time	Responds to change in heart rate within 5 to 9 seconds depending on physiological waveform. (As measured per AAMI standard EC 13-1992 clause 4.1.2.1 (f), including 3.1.2.1 parts f. and g. waveforms.) Includes 1 second readout update interval.
HR Accuracy	$\pm 3$ beats per minute or 3%, whichever is greater NOTE: AAMI Test 4.1.4 part f: Accuracy is affected (i.e., rate drops) when QRS and pacer spikes are nearly simultaneous as occasionally is the case during this AAMI test.
Heart Rate Averaging Method	Heart rate = 60 / latest average interval in seconds For higher heart rates, latest average interval = 7/8 of previous average interval + 1/8 of latest interval. For lower heart rates, latest average interval = 3/4 (previous average interval) + 1/4 latest interval. Transition rates for choice of formula include hysteresis and are 70 and 80 beats per minute.
Drift Tolerance (AAMI Specification EC13-1992, 3.2.6.3)	80 beats per minute indicated for 80 beats per minute ECG plus drift waveform
Pacer Display	Pacer indicator shown on screen if PACER function turned on; pacer spike always shown if of sufficient amplitude.
Pacer Pulse Rejection	Pacer detection range (i.e., will show the dashed vertical marker) for 0.1 ms pulses is $\pm 3$ mV to $\pm 700$ mV, and drops linearly to $\pm 2$ mV to $\pm 700$ mV for 0.2 to 2 ms pulses. Will not count as heartbeats approximately 95% of pacemaker pulses within pacer detection range, with or without AAMI (EC13 1992) tails of 4, 25, 50, 75, or 100 ms decay time constant, whose tail amplitudes are 2.5% or 25%, 2mV maximum, whether ventricular only, or A-V sequential pulses, all per AAMI tests 3.1.4.1 and 3.1.4.2
Response to Irregular Rhythm (AAMI specification EC13-1992, 3.1.2.1. Part e.)	
Ventricular Bigeminy (VB)	78 to 81 bpm (80 bpm expected)
Slow Alternating VB	57 to 65 bpm (60 bpm expected)
Rapid Alternating VB	118 to 123 bpm (120 bpm expected)
Bidirectional Systole	88 to 93 bpm (90 bpm expected)
1mV Ventricular Tachycardia	197 to 198 bpm (206 bpm expected)
2mV Ventricular Tachycardia	193 to 197 bpm (206 bpm expected)

## Real-Time ECG Analog/Defib Sync




## Real-Time ECG Analog/Defib Sync

Special cables are required to interface the defib sync connector to the Physio-Control LIFEPAK 5 or LIFEPAK 6s defibrillator.

### Real-Time ECG Analog/Defib Sync Signals

Signal	Specification
Sync Output	0 to 5V pulse, 100 $\pm$ 5ms wide, starts within 35 ms after peak of R-wave. 15 mA short circuit current.
Real-time ECG Output	Range = $\pm$ 6V minimum, centered about 0V, Gain = 1000X, noninverting for lead II, inverting for all other leads, delay <3 msec, 0.05-100 Hz, going to -5.9V $\pm$ 5% during ECG lead fail. V lead has no Real-Time analog output.
Marker Input (Defib Sync only)	Normally 0V in, a pulse either + or -3 to $\pm$ 15V for 10-70 ms puts a marker in ECG trace. $\sim$ 5 k $\Omega$ input res.
Shield	Common terminal for other signals

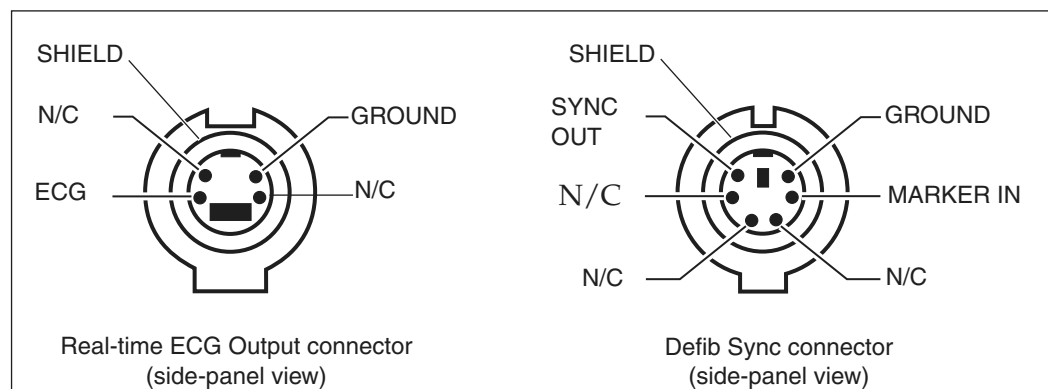The sync and real-time ECG outputs do not operate during In-service mode.



## Impedance Pneumography (RESP)

### Impedance Pneumography (RESP) Specifications

Characteristic	Specification
<b>Display Characteristics</b>	
Sweep speed	3.13, 6.25, 12.5 mm/sec; user-selectable
Amplitude range	1x, 2x, 4x, 8x, 16x
<b>Performance Characteristics</b>	
Excitation signal characteristics	65 $\mu$ A RMS $\pm$ 5% at 63.0 kHz pseudo sine wave
Sensing electrodes	User selectable RA-LA or RA-LL
Base impedance (in addition to 1K resistors in ECG cables)	100 to 1200 ohms is normal monitoring range, approx. 1200-1500 ohms range produces a "NOISY SIGNAL, CHECK ELECTRODES" equipment alert. Above approx. 1500 ohms produces a "RESP FAULT, LEAD FAIL" equipment alert. Thresholds are dependent on ECG cable type.
Impedance dynamic range	20 ohms
Signal bandwidth after detection	0.06 Hz (single pole) to 3.2 Hz (2 pole)
Breath detection threshold	140 milliohms or 2x CVA, whichever is greater
Respiration rate range	Adult/Ped: 0 (apnea), 2 to 150 breaths/min Neonate: 0 (apnea), 3 to 150 breaths/min
Respiration rate accuracy	$\pm$ 2 breaths/min or $\pm$ 2%, whichever is greater
Respiration rate source (RR)	When CO <sub>2</sub> is active, CO <sub>2</sub> is the BR source. Otherwise, RESP from ECG is the RR source.
Apnea alarm delay accuracy	+1 second
Resolution	5 seconds
Apnea alarm delay settings	Central apnea only - alarm delay is set by the user Adult/Ped = 6, 10, 15, 20, 25, 30 Neonate = 6, 10, 15, 20 seconds
Cardiovascular artifact rejection (CVA)	Presence of CVA is detected automatically. Breaths will be picked in the presence of CVA unless the Breath Rate is within 5% of the Heart Rate or a sub-multiple of the heart rate.
Motion artifact rejection	not rejected
Obstructive apnea	not detected
Cardiovascular artifact rejection (CVA)	Presence of CVA is detected automatically. Breaths will be picked in the presence of CVA unless the Breath Rate is within 5% of the Heart Rate or a sub-multiple of the heart rate.
Motion artifact rejection	not rejected
Obstructive apnea	not detected

## Invasive Pressure

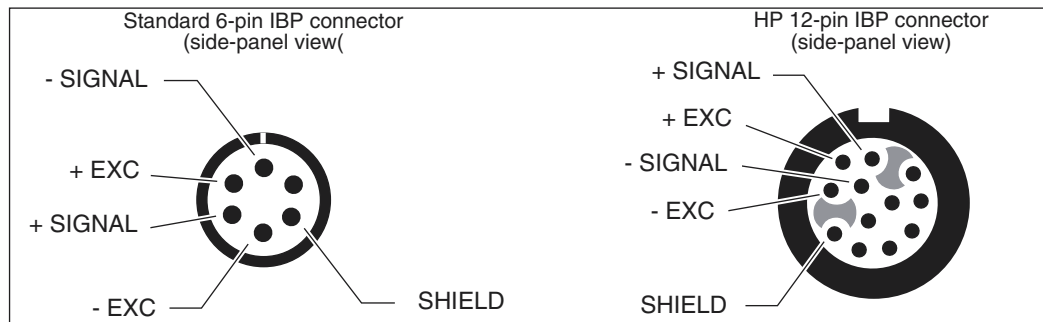
Applies only to models 204 and 206.

### Invasive Pressure Specifications

Characteristic	Specification
Transducer Type	Strain-gauge resistive bridge, or HP quartz (with HP Option). <sup>a</sup>
Transducer Excitation Impedance Range	200 to 2000 $\Omega$
Transducer sensitivity	5 $\mu\text{V/V/mmHg}$
Excitation Voltage	4.85V Pulsed dc @ 181 Hz <sup>b</sup>
Connector	ITT-Cannon plug MS3106F-14S-6P Std. Hewlett-Packard compatible 12-pin connector (optional).
Bandwidth	Digital filtered, dc to 20 Hz
Zero Drift	$\pm 1$ mmHg without transducer drift
Zero Adjustment	$\pm 200$ mmHg including transducer offset
Numeric Accuracy	$\pm 2$ mmHg or 2% of reading, whichever is greater, plus transducer error
Pressure range	-30 to 300 mmHg
Pulse range	25 to 250 beats per minute
Leakage Current	Meets ANSI/AAMI risk (leakage) requirements
Electrosurgery interference suppression	Included in all models

a. Transducers with 40  $\mu\text{V/V/mmHg}$  sensitivity are not compatible.

b. Duty factor depends on transducer impedance. For 200 to  $\sim 900 \Omega$ , duty factor is  $\approx 11\%$ . Above  $\sim 900 \Omega$ , the duty factor increases to  $\approx 91\%$ .



# NIBP

## NIBP Specifications

Characteristic	Specification
Method	Oscillometric
Control	Automatic and manual measurement control
Auto Intervals	1, 2, 3, 5, 10, 15, 30, and 60 minutes
Turbocuf	Maximum measurements allowable in a 5-minute period
Displayed Pressures	Systolic, Diastolic, and Mean plus on-screen manometer
Systolic Range	Adult: 30 to 260 mmHg Ped: 30 to 160 mmHg Neonate: 25 to 120 mmHg
Diastolic Range	Adult: 20 to 235 mmHg Ped: 15 to 130 mmHg Neonate: 10 to 105 mmHg
Mean Range	Adult: 20 to 255 mmHg Ped: 15 to 140 mmHg Neonate: 10 to 110 mmHg
Static Manometer Accuracy	±3 mmHg
Minimum Inflation Pressure	Adult: 100 mmHg Ped: 80 mmHg Neonate: 50 mmHg
Maximum Allowable Pressure	Adult: 270 mmHg Ped: 170 mmHg Neonate: 132 mmHg
Default Inflation Pressure	Adult: 160 mmHg Ped: 120 mmHg Neonate: 90 mmHg
Normal Overpressure Limit (results in up to 2 retries)	Adult: 280 mmHg Ped: 200 mmHg Neonate: 141 mmHg
Single Fault Overpressure Limit	Adult: 308 mmHg Ped: 220 mmHg Neonate: 154 mmHg
Leak Rate	After a 1 minute settling period, leak rate is ≤4mm/Hg over a 3-minute period at 270 mm/Hg.
Pulse Rate Range	30 to 220 beats per minute
Maximum Determination Time (with retries)	Adult: 4.5 minutes Ped: 4 minutes Neonate: 3 minutes
Maximum Determination Time (no retries)	Adult: 3 minutes Ped: 2 minutes Neonate: 1.5 minutes
Typical Determination Time without Artifact	30 to 45 seconds

## Temperature

Characteristic	Specification
Minimum Time between automatic measurements	30 seconds (Auto Mode) 2 seconds (Turbo Mode)
Electrosurgery Interference Suppression	Included in all models: 202EL, 204EL, 206EL
NIBP Performance	Per EN 1060-1, EN 1060-3, and ANSI-AAMI SP10-1992
NIBP Safety	Per IEC 601-2-30

## Temperature

### Temperature Specifications

Characteristic	Specification	
Range	0° to +50°C; 32° to +122°F	
Displays	T1, T2, and $\Delta T$	
Probes	Compatible with YSI Series 400 and 700 probes. HP side panel only compatible with YSI 400 and has HP connector.	
Units	°C and °F selectable	
Channel Accuracy	Temperature Range	Tolerance
	0° to +10°C	±0.2°C
	>10° to +50°C	±0.1°C
	+32° to +50°F	±0.4°F
	>50° to +122°F	±0.2°F
Resolution	0.1°C or °F	
Electrosurgery interference suppression	Included in all models: 202EL, 204EL, 206EL	

## Pulse Oximetry (SpO<sub>2</sub>)

### Masimo SpO<sub>2</sub> Specifications

Characteristic	Specification
Saturation (% SpO <sub>2</sub> ) Range Resolution Alarm Limits	1% to 100% 1% 52% to 100% (upper) 50% to 98% (lower)
Probe Accuracy (25° to 41° C) Adults, Pediatrics: No motion  Neonates: No motion  Adults, Pediatrics, Neonates: During Motion <sup>a,b</sup>	70% to 100% ±2 counts 0% to 69% unspecified  70% to 100% ±3 counts 0% to 69% unspecified  70% to 100% ±3 counts 0% to 69% unspecified
Pulse Rate Range: No motion Range: During motion <sup>a,b</sup> Resolution Alarm Limits	26 to 239 beats per minute, ±3 counts 26 to 239 beats per minute, ±5 counts 1 beat per minute 27 to 250 beats per minute (upper) 25 to 248 beats per minute (lower) Note: Any pulse rate above 239 will activate the pulse rate alarm, even if the upper alarm limit is set above 239. If the lower alarm limit is set to 25, a pulse rate of 25 will activate the pulse rate alarm due to the limitation of the displayable numeric range.
Pulse Rate Accuracy No Motion During Motion <sup>a,b</sup>	±3 beats per minute ±5 beats per minute
Measurement averaging time	8 seconds
Alarm Hold-Off Time Period	10 seconds; resets if the sensor reports levels within limits before 10 seconds elapses
Circuitry	Microprocessor controlled Automatic self-test of oximeter when powered on Automatic setting of default parameters Automatic alarm messages
Electrosurgery interference suppression	Yes
Sensor Compatibility	Compatible only with Masimo sensors listed in the Welch Allyn <i>Products and Accessories</i> booklet.
Sensor LEDs RED Wavelength INFRARED Wavelength	660 nm (nominal) 905 nm (nominal)
Sensor Energies (Radiant Power)	0.13 mW to 0.79 mW at 50 mA pulsed

a. Motion for adults and pediatrics is defined as rubbing and tapping motions at 2 to 4 Hz at an amplitude of 1 to 2 cm and a non-repetitive motion between 1 to 5 Hz at an amplitude of 2 to 3 cm in induced hypoxia studies in the range of 70-100% SpO<sub>2</sub> against a laboratory co-oximeter and ECG monitor. This variation equals ±1 standard deviation which encompasses 68% of the population.

b. Motion for neonates is defined as foot motions at 2 to 4 Hz at an amplitude of 1 to 2 cm against a laboratory co-oximeter and ECG monitor. This variation equals ±1 standard deviation which encompasses 68% of the population.



## Nellcor SpO<sub>2</sub> Specifications

Characteristic	Specification <sup>a</sup>
Saturation (% SpO <sub>2</sub> ) Range Resolution Alarm Limits <sup>b</sup>	0% to 100% 1% 52% to 100% (upper) 50% to 98% (lower)
Probe Accuracy <sup>c</sup> (saturation levels between 70% and 100%, 28° to 42°C) Adult/Pediatric Neonatal	Digit accuracy: ±2 counts Digit accuracy: ±3 counts
Pulse Rate Range <sup>d</sup>  Alarm Limits	Motion tolerant option: 25 to 249 beats per minute Option without motion tolerance: 25 to 250 beats per minute 27 to 250 beats per minute (upper) 25 to 248 beats per minute (lower)
Pulse Rate Accuracy No Motion During Motion	±3 beats per minute ±5 beats per minute
Alarm Hold-Off Time Period	10 seconds; resets if the sensor reports levels within limits before 10 seconds elapses
Circuitry	Microprocessor controlled Automatic self-test of oximeter when powered on Automatic setting of default parameters Automatic alarm messages
Electrosurgery interference suppression	Yes
Sensor Compatibility	Compatible only with Nellcor sensors listed in the Welch Allyn <i>Products and Accessories</i> booklet.
Sensor LEDs RED Wavelength INFRARED (IR) Wavelength	660 nm (nominal) 890 nm (nominal)
Sensor Energies (Radiant Power) Electrical Power Optical Power	52.5 mW max. 15 mW max.
Pulse Rate Accuracy No Motion During Motion	±3 beats per minute ±5 beats per minute

a. When performing SpO<sub>2</sub> monitoring, the monitor SpO<sub>2</sub> channel may not (in rare instances) initially detect a damaged SpO<sub>2</sub> sensor or extension cable (as described in EN865:1997, Section 51.109), but may simply display a flat line and no numerics. If this occurs, try connecting a new SpO<sub>2</sub> sensor (and extension cable, if needed) and restart the monitor to resume normal monitoring.

b. Minimum difference between upper and lower alarm limits is 2%.

c. Refer to the Welch Allyn *Products and Accessories* guide (810-0409-XX) for accuracy specifications for all Nellcor SpO<sub>2</sub> probes recommended for use.

d. When using the pulse oximetry option with motion tolerance, a substantial and rapid (<2 seconds) drop in pulse rate may result in erroneous pulse rate readings and loss of the audible pulse indicator.

# Capnography (CO<sub>2</sub>)

## General CO<sub>2</sub> Specifications

Characteristic	Specification
<b>CO<sub>2</sub> Display</b>	
Screen Display	CO <sub>2</sub> waveform and ETCO <sub>2</sub> and INCO <sub>2</sub> (when in alarm) numerics
Numeric Display Ranges	ETCO <sub>2</sub> 0-99 mmHg, 0-13.2 kPa, 0-23.1% INCO <sub>2</sub> 8 <sup>a</sup> -25 mmHg, 1.1 <sup>a</sup> -5 kPa, 1.1 <sup>a</sup> -5%
Waveform Scale (Maximum)	0-100 mmHg, 0-14 kPa, 0-14%
Units	mmHg, kPa, %; user-selectable
Sweep Speed	3.13, 6.25, 12.5 mm/sec; user-selectable
Response Modes	Fast 15 sec sampling time period Normal 30 sec sampling time period Slow 45 sec sampling time period
Gas Compensation	OFF CO <sub>2</sub> value = calculated CO <sub>2</sub> value; O <sub>2</sub> > 50%, No N <sub>2</sub> O CO <sub>2</sub> value = calculated CO <sub>2</sub> value x 1.03; N <sub>2</sub> O > 50% CO <sub>2</sub> value = calculated CO <sub>2</sub> value x 0.952
Alarm Limit Ranges	ETCO <sub>2</sub> 0-99 mmHg, 0-13.2 kPa, 0-13.2% INCO <sub>2</sub> 2-25 mmHg, 0.2-5 kPa, % (no lower limit)
Resolution	1 mmHg
Accuracy	Mainstream <sup>b</sup> 0-30 mmHg, ±3 mmHg 31-99 mmHg, ± 10% of value Sidestream <sup>c</sup> 0-30 mmHg, ±3 mmHg 31-99 mmHg, ± 10% of value
Altitude Error	±0.4%/1,000 ft (304.8 m)
<b>Breath Rate Display</b>	
Screen Display	Numeric
Breath rate (BR) source	When CO <sub>2</sub> is active, CO <sub>2</sub> is the BR source. Otherwise, RESP from ECG is the RR source.
Units	Breaths/Minute
Range	Adult/Ped 0 (apnea), 2 to 150 breaths/min Neonate 0 (apnea), 3 to 150 breaths/min
Resolution	±1 breaths/min
Accuracy	±1 breaths/min or ±5%, whichever is greater <sup>d</sup>
Alarm Limits Range	Adult/Ped 2 to 150 breaths/min Neonate 3 to 150 breaths/min
<b>CO<sub>2</sub> Performance</b>	
Specification	Per ISO 9918:1993 (E) / EN 864:1996
<b>Apnea Alarms and Tickets</b>	
Apnea Ticket	Set to auto print after apnea event and after 1 minute continued apnea
Apnea Alarm Accuracy	± 2 sec
Apnea delay setting	Adult/Ped = 6, 10, 15, 20, 25, 30 seconds Neonate = 6, 10, 15, 20 seconds

## Capnography (CO<sub>2</sub>)

Characteristic	Specification
<b>Barometric Pressure</b>	
Pressure Compensation	Automatic
Operating Range	-2,000 to 15,000 ft (-610 to 4572 m), 817 to 429 mmHg
Screen Display	Numeric (CO <sub>2</sub> Status Window)
Units	mmHg or kPa
Accuracy	±3 mmHg or 2.5% of difference from calibration pressure, whichever is greater

a. Lower if in alarm.

b. Based on these airway conditions: sensor 42° airway adapter temperature = 33°C, water vapor pressure = 38 mmHg; standard gas mixture = CO<sub>2</sub> in balance air, fully hydrated at 33°C; barometric pressure = 760 mmHg and flow = 60 ml/min.

c. Based on the following additional airway conditions: Sample line = 7 ft, 0.055 in ID (2.13 m, 1.4 mm ID); Sample flow rate = 175 ml/min; Welch Allyn watertrap (new/unused); Respiratory rate ≤ 50 bpm, stable to ±3 breaths/min; Inspired/Expired time ratio = 1:2. Barometric pressure = 760 mmHg.

d. For Sidestream CO<sub>2</sub>, this applies only for BR ≤ 50.

Mainstream CO<sub>2</sub> Specifications

Characteristic	Specification
<b>Mainstream CO<sub>2</sub> Sensor</b>	
Sensor Type	Mainstream
Principle of Operation	Non-dispersive, infrared, single-beam, single path/wavelength, ratiometric
Warm-up time (CO <sub>2</sub> sensor and monitor)	45 sec typical, 3 min maximum
Response Time	30 ms typical, 60 ms maximum
Waveform Rise Time	<120 ms to 90% after step change
Calibration	Verify semi-annually, calibrate only as required
Sensor Housing Temperature	42°C nominal
<b>Mainstream CO<sub>2</sub> Sensor and Cable Dimensions and Weight</b>	
Sensor Height <sup>a</sup>	1.003 in
Sensor Width <sup>a</sup>	1.036 in
Sensor Depth <sup>a</sup>	0.78 in
Sensor Weight <sup>a</sup>	< .39 oz
Cable Length	10 ft (3.05 m) nominal
<b>Mainstream CO<sub>2</sub> Airway Adapter</b>	
Type	Per ISO 3040, single-use
Size	15 mm ID, (meets ISO specifications)
Material	clear polycarbonate, with sapphire windows
Added Deadspace	< 6cc (.37 cubic inches) for adult model, <0.6 cc (.037 cubic inches) for low deadspace model
<b>Mainstream CO<sub>2</sub> Sensor Environmental Specifications</b>	
Operating Ambient Temperature	10° to 40°C
Storage Temperature	-20° to 60°C
Operating Altitude	-2,000 to 15,000 ft (-610 to 4,572 m), 817 to 429 mmHg
Storage Altitude	-2,000 to 40,000 ft (-610 to 12,192 m), 817 to 141 mmHg
Operating and Storage Humidity	0% to 95%, noncondensing
Shock	100 g for 4 msec
Vibration	5-35 Hz, 0.015 in peak-to-peak, 35-100 Hz, 1 g acceleration
Drop	36 inches free fall to floor (tile over concrete, one drop each face, one drop each edge/corner)

<sup>a</sup>.not including cable

## Sidestream CO<sub>2</sub> Specifications

Characteristic	Specification
Sensor Type	Sidestream, internal
Principle of Operation	Non-dispersive, infrared, single-beam, single path/wavelength, ratiometric
Operating Ambient Temperature	5° to 40°C
Startup Time	30 seconds typical, 3 minutes maximum
Rise Time	240 ms (10% to 90%) at 175 ml/min
Delay Time	1.12 seconds maximum <sup>a</sup>
Total System Response Time	1.36 seconds maximum (Rise Time and Delay Time)
Calibration	Verify semi-annually, calibrate only as required
Sampling Chamber	Internal (replaceable by service technician)
Pneumatic and Exhaust System	Integral
Barometric Pressure Compensation	Automatic
BTPS, ATPS, STPD <sup>b</sup>	CO <sub>2</sub> value = calculated CO <sub>2</sub> value x 0.977
Sampling Line	7-foot sampling line, ID 0.055 in (1.4 mm), for use with disposable single-use cannula (CO <sub>2</sub> only or CO <sub>2</sub> sampling/O <sub>2</sub> delivery)
Watertrap	Disposable single-use
Flow Rate	90 or 175 ml/min, user-selectable

a. Based on the following additional airway conditions: Sample line = 7 ft, 0.055 in ID (2.13 m, 1.4 mm ID); Sample flow rate = 175 ml/min; Welch Allyn watertrap (new/unused).

b. BTPS (Body Temperature and Pressure, Saturated), ATPS (Ambient Temperature and Pressure, Saturated), STPD (Standard Temperature and Pressure, Dry).

# Alarms

## Alarms Specifications

Characteristic	Specification
Indicators	Flashing red ALARM light indicates patient alarm; continuously on ALARM light indicates patient alarms are suspended. Continuously on ALARM(S) OFF light indicates one or more alarm limits have been disabled. A flashing ALARM(S) OFF light indicates an equipment alert.
Tone Frequency	900 Hertz Tone is steady for a patient alarm and sounds for 1 second every 4 seconds for an equipment alert.
Selectable Tone Volume	Low, Medium, High
Limits	Settable on all parameters
Control	Automatic preset or manual settings
Alarm on Tachycardias	Most tachycardias will alarm in less than 8 seconds. These include AAMI 3.1.2.1 part f. waveforms. Certain multifocal tachycardias may initially alarm as “low rate.”
Apnea delay setting	Adult/Ped = 6, 10, 15, 20, 25, 30 seconds Neonate = 6, 10, 15, 20 seconds
Alarm Holdoff Time Period <sup>a</sup>	HR/PR = 3 seconds (except NIBP PR) SpO <sub>2</sub> = 10 seconds RR/BR = 5 seconds
Audio Alarm Holdoff with Acuity	When a Propaq Encore in Adult or Pediatric Mode is connected to an Acuity System, the audio alarms at the bedside Propaq can be delayed up to 4 minutes and 15 seconds. The delay time is selected in Acuity software at the time of Acuity installation. Visual alarm indications are not delayed.

a. Alarm holdoff time period is reset if the vital sign returns to acceptable limits before an alarm occurs.

## Trends

### Trends Specifications

Characteristic	Specification <sup>a</sup>
Model 202 Parameters	NIBP, T1, T2, $\Delta T$ , HR (heart rate/pulse rate), SpO <sub>2</sub> , End-tidal CO <sub>2</sub> , Inspired CO <sub>2</sub> , Breath Rate
Model 204 Parameters	NIBP, P1, T1, T2, $\Delta T$ , HR (heart rate/pulse rate), SpO <sub>2</sub> , End-tidal CO <sub>2</sub> , Inspired CO <sub>2</sub> , Breath Rate
Model 206 Parameters	NIBP, P1, P2, T1, T2, $\Delta T$ , HR (heart rate/pulse rate), SpO <sub>2</sub> , End-tidal CO <sub>2</sub> , Inspired CO <sub>2</sub> , Breath Rate
Duration	5 hours for non-NIBP trends (up to 150 readings) A maximum of 128 readings (up to 8 hours) for NIBP trends
Resolution	All channels except NIBP sample data at 2-minute intervals For NIBP trends, a new entry is placed in the table each time an NIBP determination is made.

a. Assumes SpO<sub>2</sub> and CO<sub>2</sub> functions are present.

## Display

### Display Specifications

Characteristic	Specification
Matrix	552 x 256 pixels EL display
Active Viewing Area	145.75 mm x 67.56 mm
Pixel Size	0.203 mm x 0.203 mm
Pixel Pitch	0.264 mm x 0.264 mm
Character Height	Large: 11.03 mm (0.434 in) Medium: 7.34 mm (0.289 in) Small: 3.64 mm (0.143 in)
Viewing Angle	>160° Horizontal and Vertical
Contrast Ratio	>45 ("On" pixel luminance/"Off" pixel luminance)
Display Color	Amber
Display Background Color	Black
"On" Pixel Luminance	>9.0 fL (area of amber pixel; includes protective window)
"Off" Pixel Luminance	<0.2 fL (black pixel)
Refresh Rate	109 Hz

## Monitor (Environmental)

### Monitor Environmental Specifications

Characteristic	Specification
Operating Temperature	0° to 40° C
Shipping and Storage Temperature	-20° to 60° C
Operating Altitude	-2,000 to 15,000 ft (-610 to 4,572 m)
Shipping and Storage Altitude	-2,000 to 40,000 ft (-610 to 12,192 m)
Operating Relative Humidity	15% to 95%, noncondensing per MIL STD 810E, Procedure 1-natural
Shipping and Storage Relative Humidity	15% to 95%, noncondensing per MIL STD 810E, Procedure 1-natural
Shock	50 g
Vibration, Random	0.02g <sup>2</sup> /Hz from 10 to 500 Hz, ramping down to 0.002g <sup>2</sup> /Hz at 2000 Hz. Operating 1 hour per axis, 3 hours per test. Designed to meet RTCA DO-160C, Category C.
Vibration, Sinusoidal	0.10 peak to peak inches 5 to 17 Hz, sloping to 0.01 peak to peak inches at 55 Hz, then sloping to 0.0001 peak to peak inches at 2000 Hz. Operating 1 hour per axis, 3 hours per test. Designed to meet RTCA DO-160C, Category N.
Electromagnetic Compatibility (EMC)	Per IEC 601-1-2, which is a collateral standard of IEC 601-1, for electromagnetic compatibility. Designed to meet RTCA DO-160C, Section 21, Category B. Propaq 200 Series SN EA000225 and higher: MIL-STD 461D:CE102, CS101, CS114, CS115, CS116, RE102, RS103 (United States Air Force).



**Caution** The monitor may not meet performance specifications if it is not used or stored within these environmental specifications.



## Monitor (Physical)

### Monitor Physical Specifications

Characteristic	Specification
<b>Protection Classifications, all Configurations <sup>a</sup></b>	
Type of Protection against Electric Shock:—Power Adapter	Power adapter class 1
Type of Protection against Electric Shock—Monitor (when connected to power adapter or powered by internal battery)	Protective earth not available in monitor. Monitor designed and tested to meet Double Insulation Requirement.
Degree of Protection Against Electric Shock, for Parts Applied to Patients	See monitor labels
Method of Disinfection	Not suitable for autoclaving <sup>b</sup>
Flammable Anesthetics	Not suitable for use with flammable anesthetics
<b>Monitor Only</b>	
Height	6.65 in (16.9 cm)
Width	8.25 in (20.9 cm)
Depth	5.10 in (12.9 cm)
Weight	6.25 lb (2.8 kg)
<b>Monitor with SpO<sub>2</sub> Module</b>	
Height	6.65 in (16.9 cm)
Width	8.25 in (20.9 cm)
Depth	7.50 in (19.10 cm)
Weight	9.12 lb (4.10 kg)
<b>Monitor with Expansion Module (Printer / SpO<sub>2</sub> / CO<sub>2</sub>)</b>	
Height	9.65 in (24.5 cm)
Width	8.25 in (20.9 cm)
Depth	7.56 in (19.2 cm)
Weight with Printer, SpO <sub>2</sub> , and CO <sub>2</sub>	13.5 lb (6.1 kg)

a. Per EN 60601-1 unless otherwise stated.

b. See Chapter 7 for cleaning instructions.

# Printer

## Printer Specifications

Characteristic	Specification
<b>Operation</b>	
Operating Modes	Continuous, Snapshot, Freeze Print, Auto Print, Auto Trend, Tabular Trend, Alarm Print, NIBP Ticket, Apnea Ticket, OxyCRG, OxyCRG on Alarm
Auto Print Intervals	15 min, 30 min, 1 hour, 2 hours, 4 hours
Auto Trend Shifts	Once every 4 hours
Number of Waveforms	Up to three: ECG, P1, P2, SpO <sub>2</sub> , CO <sub>2</sub> , RESP
Grid	5 mm and 1 mm gradations
Annotation	Date, Time, Print Mode, Speed, Heart Rate, Systolic, Diastolic, Mean, SpO <sub>2</sub> , Breath Rate, ETCO <sub>2</sub> , INCO <sub>2</sub> , Temperature, $\Delta T$ , Pacer Status, Company Logo, ECG Bandwidth, Patient Mode, scale factors for all traces and, if Acuity is connected, patient name and identification.
Printing Speeds	6.25, 12.5, 25.0 mm/sec, simulated 6.25 mm/sec for CO <sub>2</sub> and RESP in Snapshot mode
<b>Printer Mechanism</b>	
Printing Method	Thermally sensitive dot method
Dot structure	320 dots per line
Printing width	53 mm
Horizontal Dot Pitch	0.165 mm, 6 dots/mm
Vertical Dot Pitch	0.165 mm
Paper Feed Method	Friction Feed
Paper Feed Precision	$\pm 2\%$ @ 25° C and 60% Relative Humidity
Paper Width	60 mm
Reliability	30 million pulses/dot
<b>Environmental</b>	
<b>Monitor/Expansion Module</b>	
Operating Temperature	+5° to 40° C
Shipping and Storage Temperature	-20° to 60° C
Operating Relative Humidity	35% to 85% noncondensing
Shipping, Storage Relative Humidity	15% to 90% noncondensing
Operating Altitude	-2,000 to 15,000 ft (-610 to 4,572 m)
Shipping and Storage Altitude	-2,000 to 40,000 ft (-610 to 12,192 m)
Shock	30 g
Vibration, Random	0.02g <sup>2</sup> /Hz from 10 to 500 Hz, ramping down to 0.002g <sup>2</sup> /Hz at 2000 Hz. Operating 1 hour per axis, 3 hours per test.

## Power

Characteristic	Specification
Electromagnetic Compatibility (EMC)	Per IEC 601-1-2, which is a collateral standard of IEC 601-1, for electromagnetic compatibility.
<b>Paper</b>	
Short-term Storage Environment (up to 7 days)	-20 to 40°C; 5% to 80% noncondensing
Long-term Storage Environment (up to 5 years)	25°C (optimal), 65% noncondensing

## Power

### Monitor Power Specifications

Characteristic	Specification
Mode of Operation	Continuous
Battery Pack Type	Sealed, gel-type lead acid
Battery Pack Capacity	Monitor only: 8 volts, 3 Ampere-Hours; Monitor with Expansion Modules: 8 volts, 6 Ampere-Hours
Battery Recharger Circuitry	Internal, powered by external power adapter
DC Input Power Required	12 to 28 Volts, 25 Watts
Input Fuse Rating	3A/250V, Type 2AG (0.57x 0.177 in)
Operating Times on Battery	
Monitor only	2.5 hours
Monitor and SpO <sub>2</sub> (Baqpaq)	4.5 hours
Monitor with Expansion Module with printer, SpO <sub>2</sub> and CO <sub>2</sub> Options	3 hours
Battery Recharge Time with instrument on	Range of 8 hours to 12 hours typical, depending upon product configuration
Battery Recharge Time with instrument off	Range of 6 hours to 8 hours depending upon product configuration
Recharge time until monitor is usable, starting with discharged but non-faulty battery	≤ 2 minutes typically (longer time required before NIBP, printer, and CO <sub>2</sub> are available)

## Power Adapters

### Power Adapter General Specifications

Characteristic	Specification
<b>Protection Classifications, all Adapters <sup>a</sup></b>	
Type of Protection Against Electric Shock	Class I, (Protectively Earthed)
Degree of Protection Against Harmful Ingress of Water	For ordinary, indoor locations only.
Method of Disinfection	Not suitable for autoclaving <sup>b</sup>
Flammable Anesthetics	Not suitable for use with flammable anesthetics
<b>Environmental Specifications, All Adapters</b>	
Operating Temperature	0° to 50° C
Shipping and Storage Temperature	-20° to 60° C
Operating Altitude	-2,000 to 15,000 feet (-610 to 4,572 m)
Shipping and Storage Altitude	-2,000 to 40,000 feet (-610 to 12,192 m)
Operating Relative Humidity	15% to 95%, noncondensing
Shipping, Storage Relative Humidity	15% to 95%, noncondensing
Shock	50 g
Vibration	Random Vibration, 0.02g <sup>2</sup> /Hz from 10 to 300 Hz, ramping down to 0.002g <sup>2</sup> /Hz at 500 Hz. Operating 1 hour per axis, 3 hours per test.
Water Resistance	For ordinary, indoor locations only.

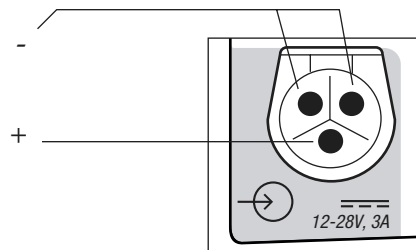
a. Per EN 60601-1 unless otherwise stated.

b. See Chapter 8 for cleaning instructions.

## Power Adapters

### Power Adapter Physical Specifications

Characteristic	Specification
<b>Universal Power Adapter, Part No. 503-0054-00, 503-0093-XX</b>	
Length	5.0 in (12.7 cm)
Width	3.6 in (9.1 cm)
Height	3.1 in (7.9 cm)
Weight	3.1 lb (1.4 kg)
Rated Input	100V-120V ac, 500 mA, 50/60 Hz
Rated Fuses	T800 mA/250V, Time-Delay, 5x20mm
Rated Output (Continuous)	16-24V dc, 25 VA
Additional Features	Detachable power cord, pilot light
<b>Universal Power Adapter, Part No. 503-0054-01, 503-0092-XX</b>	
Length	5.0 in (12.7 cm)
Width	3.6 in (9.1 cm)
Height	3.1 in (7.9 cm)
Weight	3.1 lb (1.4 kg)
Rated Input	200V-240V ac, 250 mA, 50/60 Hz
Rated Fuses	T400 mA/250V, Time-Delay, 5 x 20mm
Rated Output (Continuous)	16-24V dc, 25 VA
Additional Features	Detachable power cord, pilot light



Power Adapter  
connector

## Factory Default Settings

### Factory Default Settings

Characteristic	Specification
Date <sup>a</sup>	MO/DA/YR, DA.MO.YR, or YR/MO/DA
Decimal <sup>a</sup>	. (Period)
HR/PR Sweep	25 mm/sec
RR/BR Sweep	6.25 mm/sec
Alarm Tone	MEDIUM
HR/PR TONE	LOW
HR/PR SOURCE	ECG
RR/BR Source	CO <sub>2</sub> if ON or ECG (if CO <sub>2</sub> OFF)
Patient Mode	Adult
ECG Bandwidth	Monitor
ECG Size	1 mV/cm
ECG Lead	II
ECG Filter <sup>a</sup>	60 Hz
ECG Pacer	ON
RESP size	2X
RESP lead	Ld2
RESP sweep	6.25 mm/sec
RESP ON/OFF	ON
RESP window	ON
IBP Range	0 to 180 mmHg
IBP Rescale	0 to 140 mmHg
IBP Mode	RESCALE
Invasive Pressure Formats	Label dependent
NIBP Mode	MANUAL
NIBP Auto Time	15 min
SpO <sub>2</sub> SIZE	2x
SpO <sub>2</sub> C-LOCK	OFF
SpO <sub>2</sub> Response	NORMAL
TEMP F/C <sup>a</sup>	Celsius
CO <sub>2</sub> Range	0 to 60 mmHg
CO <sub>2</sub> Sweep	6.25 mm/sec
CO <sub>2</sub> Response	NORMAL
CO <sub>2</sub> Units <sup>a</sup>	mmHg
CO <sub>2</sub> Gas Compensation	OFF
Sidestream CO <sub>2</sub> Flow Rate	Adult: 90 ml/minute Ped: 90 ml/minute Neonate: 90 ml/minute (The flow rate cannot be programmed to a different value in a Custom Patient Mode. See Custom Patient Modes.)

## Factory Default Settings

Characteristic	Specification
Display Wave Select	All waves are on except NIBP
Trend Group	NIBP
Alarm Limits	All are ON except P2
HR Limits	Adult: 50, 120 beats per minute Ped: 50, 150 beats per minute Neonate: 100, 200 beats per minute
NIBP Limits - Systolic	Adult: 75, 220 mmHg Ped: 75, 145 mmHg Neonate: 50, 100 mmHg
NIBP Limits - Diastolic	Adult: 35, 110 mmHg Ped: 35, 100 mmHg Neonate: 30, 70 mmHg
NIBP Limits - Mean	Adult: 50, 120 mmHg Ped: 50, 110 mmHg Neonate: 35, 80 mmHg
P1, P2 Limits - Systolic	Adult: 75, 220 mmHg Ped: 75, 145 mmHg Neonate: 50, 100 mmHg
P1, P2 Limits - Diastolic	Adult: 35, 110 mmHg Ped: 35, 100 mmHg Neonate: 30, 70 mmHg
P1, P2 Limits - Mean	Adult: 50, 120 mmHg Ped: 50, 110 mmHg Neonate: 35, 80 mmHg
SpO <sub>2</sub> Limits	Adult: 85%, 100% Ped: 85%, 100% Neonate: 80%, 95%
RR/BR	Adult: 5, 30 BrM Ped: 10, 45 BrM Neonate: 10, 60 BrM
TEMP Limits	35.0°, 37.8° C
ΔT Limits	0.0°, 2.8° C
ETCO <sub>2</sub> Limits	25, 60 mmHg (3, 8 for % and kPa)
INCO <sub>2</sub> Limits	N/A, 5 mmHg (0.7 for % and kPa)
Apnea Delay	Adult/Ped: 20 seconds Neonate: 15 seconds
<b>PRINTER SETTINGS</b>	
Printer Alarm Print	OFF
Printer Auto Print	OFF
Printer NIBP Ticket	OFF
Printer Apnea Ticket	ON
Printer Print Speed	25 mm/sec
Printer Auto Trend	OFF
Printer Trend Selections	NIBP and P1 = ON; all others = OFF
Printer OxyCRG on Alarm	OFF

a. Any time you change the Date, Filter, Temp F/C, Decimal, or CO<sub>2</sub> Units setting, the new setting also becomes the powerup default setting.

## ***In-Service Simulated Values***

### **In-Service Values**

Channel	Display	Initial Value	Alternate Value
ECG	Waveform	Normal sinus rhythm, 1mV, Lead II	Normal sinus rhythm
ECG	Heart Rate	80 beats per minute	125 beats per minute
RESP	Respiration Rate	12 breaths/minute	31 breaths/minute
P1	Waveform	Arterial	Same as Initial Value
P1	Pulse Rate	80 pulses per min	125 pulses per min
P1	Systolic	121 mmHg	120 mmHg
P1	Diastolic	79 mmHg	85 mmHg
P1	Mean	96 mmHg	103 mmHg
P2	Waveform	Pulmonary Artery	Same as Initial Value
P2	Pulse Rate	80 pulses per min	125 pulses per min
P2	Systolic	25 mmHg	25 mmHg
P2	Diastolic	9 mmHg	12 mmHg
P2	Mean	15 mmHg	18 mmHg
NIBP	Mode	Manual (Auto cancels in-service)	Same as Initial Value
NIBP	Numerics	Actual values from patient	Actual values from patient
T1	Numeric	37.0° C	39.1° C
T2	Numeric	36.4° C	37.4° C
$\Delta T$	Numeric	0.6° C	1.7° C
SpO <sub>2</sub>	Waveform	Normal, 2x	Same as Initial Value
SpO <sub>2</sub>	Rate	80 pulses per min	125 pulses per min
SpO <sub>2</sub>	Numeric	97%	88%
CO <sub>2</sub>	Waveform	Normal	Hyperventilating
CO <sub>2</sub>	ETCO <sub>2</sub> Numeric	38 mmHg	60 mmHg
CO <sub>2</sub>	INCO <sub>2</sub> Numeric	0 mmHg	8 mmHg
CO <sub>2</sub>	Breath Rate	12 breaths/minute	31 breaths/minute



# Glossary

---

## **$\Delta T$**

Difference temperature. The difference between T1 and T2.

## ***AAMI***

Association for the Advancement of Medical Instrumentation (United States of America).

## ***AC Power Adapter***

The device that plugs into the 12-28V dc receptacle on the Propaq Encore's side panel to allow operation and battery charging from ac mains.

## ***Acuity***

Welch Allyn's trade name for its central station patient monitoring system.

## ***Altimeter***

A sensor, internal to the Propaq Encore, that measures absolute atmospheric pressure, and is used to correct CO<sub>2</sub> numerics for varying altitudes.

## ***ANSI***

American National Standards Institute

## ***Apnea***

Condition of no respiration occurring during a prescribed time interval.

## ***ART***

Arterial (label for an invasive blood-pressure channel).

## ***Arterial Blood Gas Measurements***

Laboratory value reporting acid-base, oxygenation and ventilation status.

## ***Artifact***

An unwanted disturbance to or by the patient or attached sensors that adds errors (usually erratic) to the measured parameters, e.g., muscle motion or shivering, electrical interference, vibration of the cuff, etc.

## ***Auto Interval***

The interval at which NIBP measurements are initiated when operating in the automatic mode.

## ***Bell***

The symbol that appears in a window to indicate alarm limits status. If alarm limits have been set, a bell appears.

## ***BP***

Blood pressure

## ***bpm***

Beats per minute

## ***Blood Pressure Numerics Windows***

The two larger windows below the heart rate. These windows can display invasive pressures and NIBP pressures.

## ***BR***

Breath rate, expressed in units per minute or 1/min. BR is derived from CO<sub>2</sub>. See also RR.

## ***Buttons***

The five buttons along the bottom-front of the Propaq Encore. A menu appears above each button identifying what each button will do when pressed.

## ***Capnogram***

Hard copy of the ET<sub>CO</sub><sub>2</sub> waveform over time.

## ***Capnometer***

Analyzer used to measure CO<sub>2</sub>, specifically ET<sub>CO</sub><sub>2</sub>.

## ***Channel***

See Patient Channel.

## ***C-Lock***

A processing scheme used in SpO<sub>2</sub> that uses QRS timing to improve the noise tolerance of SpO<sub>2</sub> measurements.

## ***CO<sub>2</sub>***

A patient channel indicating the by-product of respiration, carbon dioxide, which is exhaled by the lungs.

## ***Configuration***

The patient channels included with each Propaq Encore model. A table in Chapter 1 lists the configuration of each Propaq Encore model.

## ***Cursor***

The highlighted block in a status window that indicates the selection you make by pressing the **NEXT** button.

## ***CVA***

Cardiovascular artifact.

## ***CVP***

Central venous pressure (label for an invasive blood pressure channel).

## ***DC Offset***

The DC voltage difference between ECG electrodes. DC offset is caused by using dried out electrodes or electrodes of dissimilar metal types.

## ***Difference Temperature***

The difference between T1 and T2. Also called delta T ( $\Delta T$ ).

## ***Digital Filter***

A computer program in the Propaq Encore that removes unwanted noise that can be induced into the ECG signal from ac mains.

## ***EL (Electroluminescent) Display***

The display screen used in the Propaq Encore.

## ***EMI***

An acronym for Electromagnetic Interference.

## ***Endotracheal Tube***

Plastic breathing tube placed into the patient's windpipe.

## ***Equipment Alert***

Occurs when the Propaq Encore detects an equipment condition requiring operator assistance. A message describing the condition is displayed.

## ***Equipment Alert Window***

The window that appears during an equipment alert.

## ***Error Message***

The message that appears when the monitor detects a malfunction requiring factory service.

## ***Error Message Window***

The window that appears when the monitor detects a malfunction requiring factory service. This window contains error messages and numbers.

## ***Error Number***

The number that identifies a problem encountered during operation.

## ***ESD***

An acronym for Electrostatic Discharge (from static electricity).

## ***ESIS***

An acronym for Electrosurgery Interference Suppression.

## ***ETCO<sub>2</sub>***

An acronym for end-tidal CO<sub>2</sub>. Amount of CO<sub>2</sub> breathed out at the end of an exhalation.

## ***Factory Default Settings***

The current values for all Propaq Encore settable functions when the monitor was shipped from the factory.

## ***Freeze***

The action taken by the **FREEZE** button to stop the display. If three waveforms are displayed, all waveforms are frozen. If less than three waveforms are displayed, the current waveforms are frozen and the top waveform is also shown in real-time. See also Unfreeze.

## ***Gas Compensation***

A correction factor required to obtain accurate CO<sub>2</sub> readings when elevated levels of O<sub>2</sub> or N<sub>2</sub>O are present in respired gases.

## ***Heart Rate Source***

See Heart Rate/Pulse Rate Source.

## ***Heart Rate/Pulse Rate***

The heart rate derived from the heart rate/pulse rate source and expressed in units per minute or 1/min. See also Heart Rate/Pulse Rate Source.

## ***Heart Rate/Pulse Rate Source***

The source from which heart rate/pulse rate is derived. This source can be ECG, any pressure, including NIBP, or SpO<sub>2</sub>. When the monitor is first turned on, the Propaq Encore determines the most likely source for heart rate: ECG (first), P1 (second), SpO<sub>2</sub> (third), P2 (fourth), and NIBP (last).

## ***Highlight***

The method of identifying a selected item on the display. Highlighted selections appear as light characters on a dark background or dark characters on a light background. See also Cursor.

## ***HR***

An acronym for heart rate and expressed in units per minute or 1/min. This is displayed when the heart rate/pulse rate source is ECG.

## ***ICP***

Intracranial pressure (label for an invasive blood pressure channel).

## ***Impedance Pneumography***

A method of detecting respiratory effort by measuring the AC impedance between selected ECG leads.

## ***INCO<sub>2</sub>***

An acronym for inspired CO<sub>2</sub>. The amount of CO<sub>2</sub> measured during inhalation.

## ***In-service Mode***

A user training aid built into all Propaq Encores that provides simulated signals for all patient parameters so that function of the display, alarms, and printer can be explored easily. The in-service mode is activated by the **INSERV** button.

## ***Invasive Pressure Label***

The two or three-character label that appears in the Invasive Pressure Numerics Window identifying the source of blood pressure.

## ***Labels***

The names appearing above the buttons.

## ***Mainstream***

A respiratory CO<sub>2</sub> measurement technique which uses a noninvasive sensor located at the endotracheal tube. This technique avoids signal delays and fluid problems associated with other techniques.

## ***Menu***

A group of labels above the bottom front row of buttons on a Propaq Encore.

## ***NIBP Status Window***

The window that appears when the **NIBP** button is pressed. This window displays NIBP information.

## ***Numerics***

The numbers that appear along the top and right side of the display for heart rate, blood pressure, temperature, etc.

## ***OxyCRG***

An oxygen cardiorespirogram, a graph showing heart rate, **SpO<sub>2</sub>**, and a condensed respiratory waveform.

## ***P1***

A generic label for invasive pressure channel one.

## ***P2***

A generic label for invasive pressure channel two.

## ***PA***

Pulmonary artery (label for an invasive blood pressure channel).

## ***Parameter***

See Vital Sign Parameter.

## ***Patient Alarm***

The condition that exists when a vital sign parameter numeric violates an alarm limit.

## ***Patient Channel***

ECG, P1, P2, T1, T2, SpO<sub>2</sub>, CO<sub>2</sub>, NIBP and RESP.

## ***Patient Mode***

Selects Adult, Pediatric, or Neonatal mode settings for the monitor. These settings determine default alarm limits, maximum cuff inflation pressure, and other internal settings.

## ***Pinout***

The signal descriptions for each pin of a connector.

## ***Polarization***

The activity that occurs when dissimilar metals between ECG electrodes and leads meet. This can cause dc offset and other signal problems.

## ***PR***

Pulse rate, expressed in units per minute or 1/min. This is displayed when the heart rate/pulse rate source is from a pressure channel or SpO<sub>2</sub>.

## ***Pulse Rate***

The heart rate determined from either a pressure channel, SpO<sub>2</sub>, or NIBP; expressed in units per minute or 1/min.

## ***Pushbutton***

See Buttons.

## ***Range Mode***

The method used in invasive pressure display to show two waveforms against the same pressure scale.

## ***Rescale Mode***

The method used in invasive pressure display to show each waveform against its own scale. The scale is automatically selected for best viewing of the entire waveform.

## ***Respiration***

The exchange of oxygen and carbon dioxide in the lungs and with the cells of the body.

## ***RR***

Respiration rate, a measure of the frequency of respiration. See also Impedance Pneumography.

## ***Sensors***

The electrodes, transducers, probes, etc. used to obtain patient information.

## ***Serial Number***

The unique number assigned to the monitor. It is located on the rear panel label.

## ***Sidestream***

A respiratory CO<sub>2</sub> measurement technique which can be used for intubated or non-intubated patients.

## ***Software Version Number***

The unique number assigned to the version of the Propaq Encore's internal programming. This number appears in the Startup window.

## ***SpO<sub>2</sub>***

The standard term assigned to measuring oxygen saturation using a pulse oximeter. The SpO<sub>2</sub> patient channel noninvasively measures oxygen saturation of arteriolar hemoglobin at a peripheral measurement site, such as a finger, toe, or the bridge of the nose.

## ***Startup Window***

The information window that appears while the monitor performs its powerup test just after you turn on the Propaq Encore. This information includes the Propaq Encore model number and software version number.

## ***Status Window***

A window that appears and contains information about the Propaq Encore.



## ***SYNC***

Synchronization. Two uses apply:

- A digital output pulse from the right side panel that starts within 35 msec of the peak of a QRS complex and is used for cardioversion.
- A message in the SpO<sub>2</sub> display indicating successful C-Lock.

## ***Temporary Patient Alarm***

An alarm limit violation that occurred and was corrected without operator intervention.

## ***Trend***

The accumulation of several hours of data at two-minute intervals.

## ***Trend Parameter***

Heart Rate/Pulse Rate, P1, P2, SpO<sub>2</sub>, INCO<sub>2</sub>, ETCO<sub>2</sub>, temperature, NIBP and RR.

## ***Turbocuf Mode***

The mode used to acquire as many NIBP measurements as possible in five minutes.

## ***UA***

Umbilical artery (label for an invasive blood pressure channel).

## ***Unfreeze***

Returns the waveforms to active display. See also Freeze.

## ***UV***

Umbilical vein (label for an invasive blood pressure channel).

## ***Vital Sign Parameter***

The measurements obtained from patient channels (such as, heart rate, systolic, diastolic, mean, pulse rate, SpO<sub>2</sub>, CO<sub>2</sub>, etc.).

## ***Waveform Window***

The area in which waveforms are displayed.

## ***Waveform/Status Window***

See Waveform Window or Status Window.

***Window***

An area on the display screen in which information is displayed.

***YSI***

An acronym for Yellow Springs Instrument Company.

***Zeroing***

The process by which an invasive pressure zero reference is obtained so that pressures can be related to atmospheric pressure. This process also nulls out any residual pressure indicated by a transducer with zero pressure applied.

# Index

## A

- ac adapter connector 18
- Accessories
  - Products & Accessories book 13
- Acuity System 115
  - Acuity menu 117
  - connecting 115
  - Intended use 115
  - Network Fault message 118
  - printing at Acuity 118
- Adult Mode 56
- Alarm Holdoffs 97, 150
- ALARM light 19
- Alarm Limits Menu 27
- ALARM PRINT setting 35, 111
- ALARM TONE setting 29
- ALARM(S) OFF light 19
- Alarms
  - bells 91
  - changing individual limits 99
  - definitions and indications 93
  - flashing numerics 95
  - in-service mode 91
  - life-threatening 94
  - lights 93
  - limits for CO<sub>2</sub> 77
  - Menu 27, 91
  - programming limits 100
  - setting all 99
  - setting quickly (STAT SET) 98
  - specifications 150
  - status window 91
  - temporary 94
  - tone volume 91
  - tones 93
  - visual indications 95
- ALARMS button 29
- Alarms Menu
  - ALL ALRM button 99
  - STAT SET button 98
- ALL ALRM button 99
- Apnea 77
  - alarm limit 83
  - Ticket 110
- APNEA TICKET setting 35
- ART, label for arterial pressure 54
- Artifact affecting ECG 45
- Artifact reduction with C-LOCK (SpO<sub>2</sub>) 72
- Artifact with NIBP 56, 60, 65, 104, 109, 113
- AUTO PRINT setting 35

- AUTO/MAN button 61
- Automatic intervals (NIBP) 62
- Automatic NIBP measurements 61, 62
- Automatic trend prints 114

## B

- Battery
  - care and maintenance 123
  - charging light 19
  - removal, effect on programmable settings 127
  - replacement 126
  - storage 127
  - voltage 20, 124
- Bells
  - alarm limit status 91
- Breath rate measurements 77
- Buttons
  - menu 26
  - monitor 25

## C

- Cable circuitry
  - defibrillator synchronization 131
- CANCEL button 54
- Caution messages 101
- CHANGE button 29, 35
- Changing printer paper 128
- Cleaning instructions 125
- C-LOCK, SpO<sub>2</sub> measurement 72
- CO<sub>2</sub> 16
  - accessories 80, 84
  - alarm limits 77
  - and intubated patients 76, 86
  - and N<sub>2</sub>O 82
  - gas compensation 75
  - intended use 75
  - Mainstream CO<sub>2</sub> 80
  - Mainstream CO<sub>2</sub> messages 86
  - measurements 76
  - Menu 1 26
  - Menu 2 26
  - N<sub>2</sub>O 75
  - printing waveforms 108
  - Sidestream CO<sub>2</sub> 83
  - Sidestream CO<sub>2</sub> connector 85
  - Sidestream CO<sub>2</sub> messages 88
  - specifications 146
  - sweep speed 82

## Index

- CO2 Menu 1
  - mm/sec button 78, 82
  - RANGE button 78, 81
- CO2 Menu 2
  - FLOWRATE button 79
  - GAS COMP button 78, 79, 82
  - RESPONSE button 78, 79, 82
  - SOURCE button 78, 79
- Connectors
  - ac adapter 18
  - defibrillator synchronization 131, 135
  - ECG 41
  - INV PRS 50
  - NIBP 58
  - patient 22
  - power input 18
  - temperature probes 66
- CONTINUOUS (print speed) setting 35
- Cuff (NIBP)
  - placement 59
  - sizes 58
- CURRENT SOURCE setting 29
- Custom patient modes 33
- Customer services 130
- CVP, label for central venous pressure 54

## D

- DATE setting 32
- dc power source 18
- DECIMAL setting 33
- Default alarm limits
  - Patient Mode 20
- Defibrillator synchronization
  - cable circuitry 131
  - cardioversion using LIFEPAK 5 133
  - cardioversion using LIFEPAK 6 134
  - connector 131
  - installing interface 131
  - markers 133, 136
  - message 136
  - right side panel connector 132
  - specifications 139
  - SYNC button 133, 135
- Defibrillators
  - use with ECG 40
- Disinfection instructions 125
- Display
  - Propaq screen 24
  - specifications 151
  - SpO2 69, 71
  - waveform priorities 31
- Display trends 104
- Displaying

- ECG 45
- INV PRS 53
- NIBP 60

## E

- ECG
  - accessories 41
  - artifact 45
  - connector 41
  - displaying 45
  - electrode placement 43
  - filter 48
  - lead fault 48
  - Menu 1 26
  - Menu 2 26
  - messages 48
  - pacemaker patients 47
  - pre-gelled electrodes 42
  - printing waveforms 107
  - specifications 137
  - Status Window 46
  - waveform priorities 31
- ECG Bandwidth 40, 46
- ECG BANDWIDTH setting 46
- ECG LEAD button 45
- ECG Menu 2
  - CHANGE button 46
  - NEXT button 46
- ECG SIZE button 45, 46
- ECG/RESP
  - Menu 1 45
  - Menu 2 45
- ECG/RESP button 44
- ECG/RESP Menu 1
  - ECG SIZE button 45, 46
- ECG/RESP Menu 2
  - ECG LEAD button 45
  - RESP SZE button 45, 46
- Electrodes
  - placement 43
  - pre-gelled 42
- Electrosurgical machines
  - use with ECG 40, 44
- EMP (*see* Expansion Module with Printer)
- Environmental (monitor) specifications 152
- Environmental limits
  - operating 127
  - storage 127
- Equipment Alerts
  - definitions and indications 101
  - Printer Fault 36
  - PROGRAM FAULT 21, 102
  - tone 101

## Index

ESIS 40  
Expansion Module  
    CO2 16  
    Printer 16  
    SpO2 16  
Extended Mode 40, 46  
Extended storage 127

## F

Factory defaults  
    specifications 158  
Factory patient modes 33  
False alarms with SpO2 70, 72  
FILTER setting 32, 48  
Flashing numerics 95  
FLOWRATE button 79  
FORMAT (P1/P2) buttons 54  
FREEZE/UNFREEZE button 25, 36  
Fuses  
    replacing, Universal Power Adapter 120

## G

GAS COMP button 78, 79, 82  
Gas compensation for CO2 75, 82

## H

HP connector-compatible option 17  
HR/PR TONE setting 30, 46

## I

ICP, label for intracranial pressure 54  
Impedance Pneumography (*see* RESP)  
INSERT button 28, 37  
In-service mode 37  
    alarms 91  
    simulated trend data 37  
    using with NIBP 38  
    what you can do 37  
    what you cannot do 38  
Inspect the monitor 125  
INTERVAL button 61  
Introduction 15  
Intubated patients 76  
INV PRS button 51, 52  
INV PRS Menu 1  
    CANCEL button 54  
    RANGE button 54  
    RESCALE button 54  
    ZERO (P1/P2) buttons 54  
INV PRS Menu 2  
    FORMAT (P1/P2) buttons 54

LABEL (P1/P2) buttons 54  
Invasive pressure  
    accessories 50  
    connectors 50  
    displaying 53  
    intended use 40, 50  
    Menu 1 26  
    Menu 2 26  
    messages 51, 54-55  
    numerics 53  
    numerics formats 54  
    preparing for monitoring 51  
    printing waveforms 107  
    range 50  
    rezeroing a transducer 52  
    scales 54  
    specifications 141  
    transducer sensitivity 50  
    transducer zeroing 51  
Invasive pressure labels  
    ART (arterial) 54  
    CVP (central venous pressure) 54  
    ICP (intracranial pressure) 54  
    P1 54  
    P2 54  
    PA (pulmonary artery) 54  
    selecting 52  
    UA (umbilical artery) 54  
    UV (umbilical vein) 54

## L

LABEL (P1/P2) buttons 54  
Lead fault, ECG 48  
Learning the Propaq 37  
Left side panel, patient connectors 22  
Life-threatening alarms 94  
Lights  
    ALARM 19  
    ALARM(S) OFF 19, 93  
    battery charging 19  
Limits  
    programming 100  
    setting individual 99  
Loading printer paper 128  
Loudness  
    heart tone 46  
Loudness (*see* Volume)

## M

Main Menu 26  
    ECG/RESP button 44  
    INV PRS button 51, 52

## Index

- NIBP button 60, 61
- SETUP button 29
- SpO<sub>2</sub>/CO<sub>2</sub> button 26
- MAIN MENU button 25
- Mainstream Capnography (*see* CO<sub>2</sub>)
- Mainstream CO<sub>2</sub>
  - monitoring 80
- Maintenance
  - monitor care 127
  - power adapters 119
  - printer 128
  - replacing monitor fuse 124
  - service intervals 126
- Marriott Configuration (MCL), electrode placement 43
- Menus
  - Alarm Limits 27
  - Alarms 27, 91
  - CO<sub>2</sub> Menu 1 26, 78
  - CO<sub>2</sub> Menu 2 26
  - CO<sub>2</sub> Menu2 78
  - ECG Menu 1 26
  - ECG Menu 2 26
  - ECG/RESP 45
  - INV PRS Menu 1 26, 53
  - INV PRS Menu 2 26, 53
  - Main 26
  - NIBP 26, 60
  - Patient Alarm 95
  - Printer Setup 28
  - Service 28
  - Setup Menu 1 27
  - SpO<sub>2</sub> 69, 71
  - SpO<sub>2</sub> Menu 1 26
  - SpO<sub>2</sub> Menu 2 26
  - Trends 27
  - Wave Select 28, 31
- Messages
  - caution 101
  - defibrillator synchronization 136
  - ECG 48
  - INV PRS 54–55
  - Mainstream CO<sub>2</sub> 86
  - NIBP 62–64
  - Printer Fault 36
  - RESP 49
  - Sidestream CO<sub>2</sub> 88
  - SIMULATING 37
  - TEMP 67
- mm/sec button 78, 82
- Modem-Propaq option 17
- Monitor buttons
  - FREEZE/UNFREEZE button 25, 36
  - MAIN MENU button 25

- START/STOP button 25
- SUSPEND/RESUME button 25
- Monitor care 127
  - fuse replacement 124
- Monitor Mode 40, 46
- Monitor settings
  - ALARM TONE 29
  - CURRENT SOURCE 29
  - DATE 32
  - DECIMAL 33
  - ECG BANDWIDTH 46
  - FILTER 32, 48
  - HR/PR TONE 30, 46
  - PACER DISPLAY 46
  - PATIENT MODE 30
  - RESP (ON/OFF) 46
  - RESP LEAD 46
  - SELECTED SOURCE 29
  - SWEEP (mm/s) 29
  - TEMP F/C 32
  - TIME/DAY 31
- MORE button 29, 35
- Motion artifact and NIBP 56, 60, 65, 104, 109, 113

## N

- Neonatal Mode 56
- Network Fault message 118
- NEXT button 29, 35
- NIBP
  - accessories 58
  - adult mode 56
  - Artifact 56, 60, 65, 104, 109, 113
  - automatic intervals 62
  - automatic measurements 62
  - cardiopulmonary bypass 56
  - connector 58
  - cuff placement 59
  - cuff sizes 58
  - displaying 60
  - Menu 26, 60
  - messages 62–64
  - Motion artifact 56, 60, 65, 104, 109, 113
  - neonatal mode 56
  - pediatric mode 56
  - preparing for measurements 59
  - specifications 142
  - Symbol in display and prints 60, 104, 109, 113
  - Ticket 109
  - waveform window, displayed in 61
- NIBP button 60, 61
- NIBP Menu

AUTO/MAN button 61  
 INTERVAL button 61  
 START button 60  
 STOP button 60  
 TURBOCUF button 61  
 NIBP TICKET setting 35  
 NO INSRV button 28  
 Noninvasive blood pressure (*see* NIBP)  
 Numerics, flashing 95  
 NXT TRND button 103, 105, 113

## O

Options  
   CO2 16  
   HP-compatible side panel 17  
   Printer 17  
   RESP 17  
   SpO2 16  
 Ordering 130  
 OXYCRG button 103  
 OXYCRG ON ALARM setting 35  
 Oxygen saturation display 69, 71

## P

PA, label for pulmonary artery pressure 54  
 Pacemakers  
   detection 48  
   use with patients 47  
 PACER DISPLAY setting 46  
 Paper, printer 128  
 Patient alarms, definitions and  
   indications 93  
 Patient artifact 45  
 Patient connectors 22  
   CO2 22  
   ECG 22  
   HP-option 22  
   INV PRS 22  
   NIBP 22  
   SpO2 22  
   TEMP 22  
 Patient Mode 42, 77  
   Adult 56  
   Custom 33  
   default alarm limits 20  
   factory default 20  
   Neonatal 56  
   Pediatric 56  
 PATIENT MODE setting 30  
 Pediatric Mode 56  
 Physical (monitor) specifications 153  
 Physical inspection 125

Power  
   input connector 18  
   specifications 155  
 Power adapters 119  
   specifications 156  
 Powerup patient modes 33  
 Powerup tones  
   expansion module/SpO2 module 20  
   monitor 20  
 PR TREND button 35, 113  
 PREV MENU button 35  
 PRINT button 103, 113  
 PRINT TRENDS button 36, 113, 129  
 Printer 17  
   loading paper 128  
   maintenance 128  
   print speed (continuous  
     measurements) 35  
   PRINT TRENDS button 36, 113, 129  
   setup 35  
   SIMULATED DATA 38  
   SNAPSHOT button 36, 107  
   specifications 154  
   START/STOP button 36, 107, 129  
   test strip 36  
 PRINTER button 29, 109  
 Printer Fault equipment alert 36  
 Printer Menu  
   CHANGE button 35  
   MORE button 35  
   NEXT button 35  
   PR TREND button 35  
   PREV MENU button 35  
 Printer settings  
   ALARM PRINT 35  
   APNEA TICKET 35  
   AUTO PRINT 35  
   CONTINUOUS (print speed) 35  
   NIBP TICKET 35  
   OXYCRG ON ALARM 35  
 Printer Setup Menus 28  
 Printing 109  
   alarm prints 111  
   Apnea Ticket 110  
   automatic trends 114  
   automatically 107, 114  
   CO2 waveforms 108  
   ECG waveforms 107  
   INV PRS waveforms 107  
   NIBP Ticket 109  
   NIBP Trend 109  
   OxyCRG 112  
   OxyCRG on Alarm 112  
   PR TREND button 113



## Index

PRINT button (*see also* Trends Menu) 113  
Printer Setup window 110  
RESP waveforms 108  
Setup Menu 2 109  
speed 107  
SpO2 waveforms 107  
sweep speeds with SNAPSHOT 108  
trends 113  
trends according to shift 114  
with an expansion module 107  
Probes, temperature 66  
Products & Accessories book 13  
Programmable settings  
effect of battery removal 127  
lost 21, 102  
Propaq  
buttons 25  
documentation set 13  
screen display 24  
Propaq CS  
Physical inspection 125  
Propaq-to-Acuity option overview 17  
Pulse oximetry (*see* SpO2)

## R

RANGE button 54, 78, 81  
Real-time ECG analog  
specifications 139  
RESCALE button 54  
RESP 17  
intended use 39  
messages 49  
printing waveforms 108  
specifications 140  
sweep speed 82  
RESP LEAD setting 46  
RESP setting 46  
RESP SZE button 45, 46  
RESPONSE button 78, 79, 82  
RESPONSE time (SpO2) 71  
Returning monitor for service 130  
Rezeroing a transducer 52  
Right side panel  
fuse replacement 124  
system controls 18  
R-wave amplitude 133, 135

## S

Safety summary 9  
Screen display 24  
SELECTED SOURCE setting 29  
Service intervals 126

Service Menu 28  
Setting  
all alarms 99  
individual alarm limits 99  
SETUP button 29  
Setup Menu 1 27  
ALARMS button 29  
MORE button 29  
TRENDS button 29  
Setup Menu 2  
CHANGE button 29  
MORE button 29  
NEXT button 29  
PRINTER button 29  
WAVE SEL button 29  
Side panels  
left, patient connectors 22  
right, system controls 18  
Sidestream CO2 83  
SIMULATING, on monitor display 37  
Size, waveforms 45  
SNAPSHOT button 36, 107  
SOURCE button 78, 79  
Specifications  
alarms 150  
CO2 146  
defibrillator synchronization 139  
display 151  
ECG 137  
environmental 152  
factory defaults 158  
INV PRS 141  
NIBP 142  
physical 153  
power 155  
power adapters 156  
printer 154  
real-time ECG analog 139  
RESP 140  
SpO2 144  
TEMP 143  
trends 151  
SpO2 16  
Adjust waveform size 70, 71  
C-LOCK 72  
Display 69, 71  
Menu 1 26  
Menu 2 26  
NIBP and SpO2 69, 70  
plethysmograph 107  
printing waveforms 107  
Reducing false alarms 70, 72  
relative gain 107  
Response time selection 71



- Sensor selection 69, 70
- size indicator 107
- specifications 144
- SpO2 Menus 69, 71
- Spot-Check monitoring 73
- Standby Mode 73
- SpO2/CO2 button 26
- Spot-Check SpO2 monitoring 73
- STANDBY and SpO2 73
- Standby Mode, SpO2 73
- START/STOP button (NIBP) 25, 60
- START/STOP button (Printer) 36, 107, 129
- STAT SET 98
- STAT SET button 98
- Status Windows
  - ECG 46
- STBY message, SpO2 STANDBY 73
- Storage
  - battery 127
  - monitor 127
  - printer paper 128
- SUSPEND/RESUME button 25
- SWEEP (mm/s) setting 29
- Sweep speed (CO2, RESP) 82
- Symbol
  - In NIBP displays and prints 60, 104, 109, 113
- Symbols
  - defined 11
- SYNC button (defibrillator) 133, 135
- Synchronized cardioversion
  - LIFEPAK 5 133
  - LIFEPAK 6 134
- Synchronizing R-waves, C-LOCK 72

## **T**

- Technical service and training 130
- TEMP F/C setting 32
- Temperature
  - accessories 66
  - connectors 66
  - intended use 66
  - messages 67
  - preparing for monitoring 66
  - specifications 143
- TIME/DAY setting 31
- Tone
  - alarm volume 91
  - HR/PR (volume) 46
  - powerup 20
- Trends
  - clearing data 104
  - column labels 105

- continuous monitoring 103
  - how to display 104
  - how to select 105
  - maximum number of samples 104
- Menu 27, 103
- NIBP 104
- NIBP and symbol 104
- NXT TRND button 103, 105, 113
- OXYCRG button 103
- P1 103
- P2 103
- PRINT button 103, 113
- printing each shift 114
- printing one 113
- printing several 113
- programmable default 104
- RESP 103
- specifications 151
- status window 103
- Symbol and NIBP 104
- TEMP 103
- TRENDS button 104
- TRENDS button 29, 104
- TURBOCUF button 61

## **U**

- UA, label for umbilical artery 54
- Universal Power Adapter
  - voltage selection 120
- UV, label for umbilical vein 54

## **V**

- Verifying Patient Mode 20
- Verifying powerup tone 20
- Voltage
  - battery 20
  - checking battery 124
- Volume
  - alarm tone 91
  - heart tone 46

## **W**

- WAVE SEL button 29
- Wave Select
  - Menu 28
  - window 31
- Wave Select Menu
  - INSERV button 28, 37
  - NEXT button 31
  - NO INSRV button 28
  - ON/OFF button 31
- Waveform

## **Index**

- SpO2 70, 71
- Waveforms
  - display priorities 31
  - ECG 31
  - scales 54, 81
  - size 45
- Windows
  - Propaq screen 24

## **Z**

- ZERO (P1/P2) buttons 54
- Zeroing
  - messages 51
  - transducers 51